



## Universal flexible Keyboard Cover



### e-lastic® keyboard cover – Effective protection against bacteria and germs

The e-lastic® universal keyboard protection cover is suitable for all common keyboards. The hygienic foil is put over your keyboard with a few simple steps and protects effectively against dirt and the Spreading germs and bacteria.


Even accidental spills of liquids such as coffee or water are caught. The high quality of the material also ensures a long use of the e-lastic® keyboard protection cover, even with frequent and intensive cleaning. A flexible fit with integrated elastic band allows easy handling and individual use on all standard keyboards. The operation of the keyboard is not affected.

Material of the e-lastic® cover: The keyboard protector is made of high quality thermoplastic polyurethane (TPU) based on an aromatic polyether.

#### Benefits

- reliable protection against bacteria and germs
- protects the keyboard against liquids
- protects against daily abrasion
- easy to use, easy to replace
- durable and flexible material
- disinfection with isopropyl alcohol 70%

# e-lastic

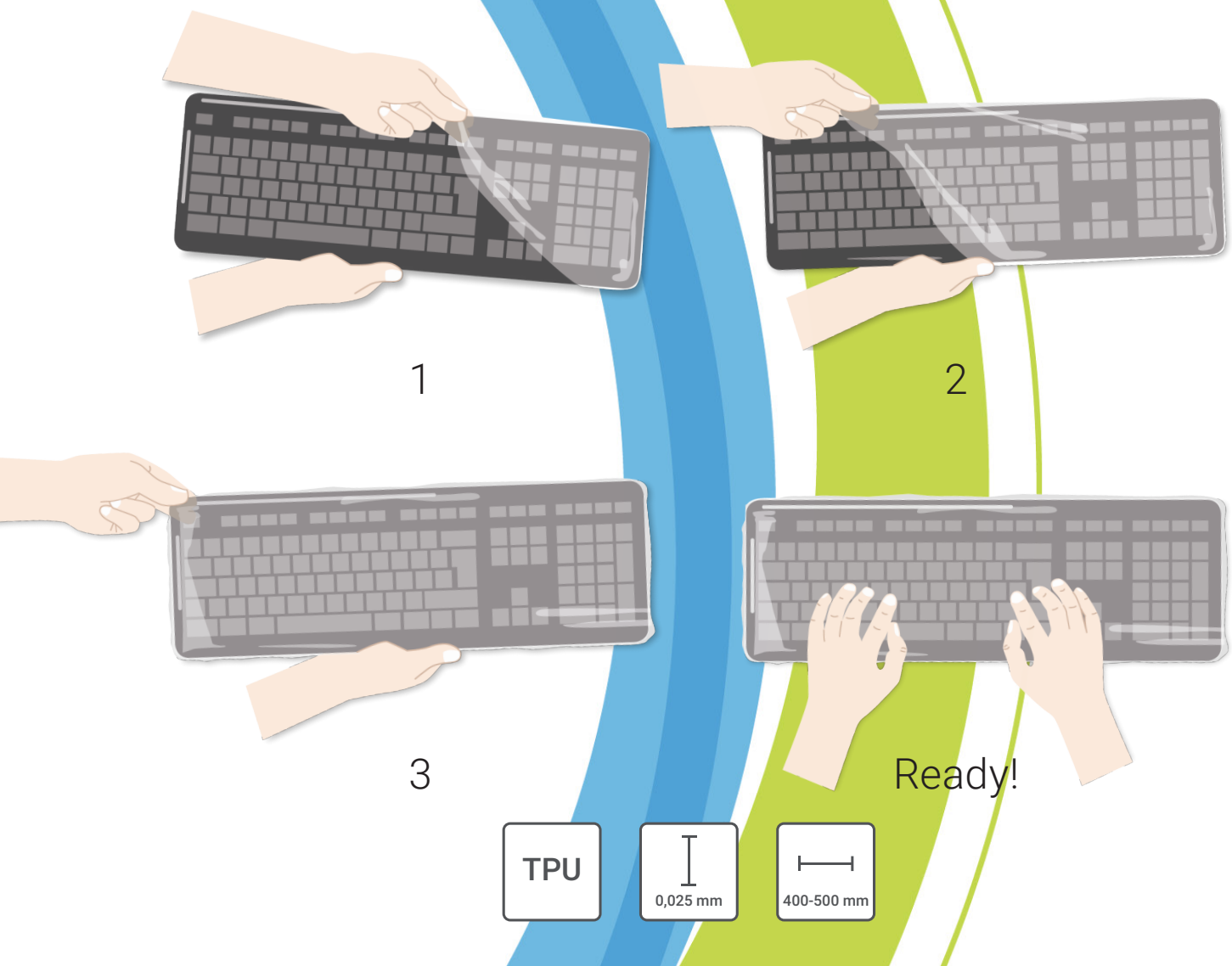

## Usage period | Cleaning | Disinfection

Depending on operational demands, the e-lastic® keyboard cover may be used for up to three months. During its life time, the cover may be used on different keyboards. Cleaning should be done only when the cover is pulled completely over the keyboard. To clean, take a slightly moistened cloth using warm water and soap and wipe the cover thoroughly. For disinfection only use 70% Isopropyl Alcohol. Do not submerge the cover in water or wipe the inside as humidity will occur.

Damaged covers should be replaced immediately. Yellowing may be an indicator of product expiration.

## Technical data

|                    |                                               |
|--------------------|-----------------------------------------------|
| Material           | high-quality thermoplastic polyurethane (TPU) |
| Material thickness | 0.025 mm                                      |
| Usage period       | up to 3 months (depending on use)             |
| Dimensions         | 400 bis 500 mm                                |
| Item no.           | 2005758                                       |