

4K3D Video Recorder

Instructions for Use

Before operating the unit, please read this manual thoroughly and retain it for future reference.

HVO-4000MT



Indications for Use/Intended Use

The Sony HVO-4000MT 4K, 3D and 2D Video Recorder is intended to record 4K video and still images from endoscopic/laparoscopic camera systems, surgical microscope and other compatible medical imaging systems. This medical grade recorder can save video and still images to an internal hard disk drive, a portable USB drive, optical disk or to a compatible network server. The recorded images can be used for the patient record, training or education. The recorder is suitable for use in hospital operating rooms, medical imaging centers, surgical centers, clinics, doctor's offices and similar healthcare environments.

Notes

- Images recorded with this equipment cannot be used for diagnostic use.
- This equipment is for medical professionals.
- This equipment is intended for use in medical environments, such as clinics, laboratory rooms, and operating rooms.

Warning

To reduce the risk of fire or electric shock, do not expose this equipment to rain or moisture.

To avoid electrical shock, do not open the cabinet. Refer servicing to qualified personnel only.

No modification of this equipment is allowed.

Warning

To avoid the risk of electric shock, this equipment must only be connected to a supply mains with protective earth.

Warning

To disconnect the main power, unplug the power plug. When installing the unit, incorporate a readily accessible disconnect device in the fixed wiring, or connect the power plug to an easily accessible socket-outlet near the unit.

Do not position the ME equipment where it is difficult to unplug the power plug.

If a fault should occur during operation of the unit, operate the disconnect device to switch the power supply off, or disconnect the power plug.

Symbols on the product



Consult the instructions for use

Follow the directions in the instructions for use for parts of the unit on which this symbol appears.



This symbol indicates the manufacturer, and appears next to the manufacturer's name and address.



This symbol indicates the Importer, and appears next to the Importer's name and registered office address.



This symbol indicates the European Community representative, and appears next to the European Community representative's name and address.



This symbol indicates the UK Responsible Person, and appears next to the UK Responsible Person's name and address.



This symbol indicates the Swiss authorized representative, and appears next to the Swiss authorized representative's name and address.



This symbol indicates the medical device in the European Community.



This symbol indicates the date of manufacture.



This symbol indicates the serial number.



This symbol indicates the Unique Device Identifier (UDI), and appears next to the bar code representation of the Unique Device Identification.



This symbol indicates the equipotential terminal which brings the various parts of a system to the same potential.



Storage and transport temperature

This symbol indicates the acceptable temperature range for storage and transport environments.



Storage and transport humidity

This symbol indicates the acceptable humidity range for storage and transport environments.



Storage and transport pressure

This symbol indicates the acceptable atmospheric pressure range for storage and transport environments.

For customers in the U.S.A.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense. You are cautioned that any changes or modifications not expressly approved in this manual could void your authority to operate this equipment.

All interface cables used to connect peripherals must be shielded in order to comply with the limits for a digital device pursuant to Subpart B of part 15 of FCC Rules.

For customers in the U.S.A.

Caution

Federal law (United States of America) restricts this device to sale by or on the order of a licensed healthcare practitioner.

Rx
ONLY

For customers in Canada

CAN ICES-3 (A)/NMB-3(A)

For customers in Canada

This unit has been certified according to Standard CAN/CSA-C22.2 No. 60601-1.

Important safeguards and notices for use in the medical environments

1. All devices connected to the unit must be certified or compliant according to IEC 60601-1, IEC 60950-1, and IEC 60065 standards and other IEC/ISO standards applicable to the devices.
2. Furthermore, the system as a whole must comply with IEC 60601-1 standards. All peripheral devices connected to the signal input/output sections of the unit constitute the medical-use system, and therefore, the user is responsible for ensuring that the system as a whole complies with IEC 60601-1 standards. If in doubt, consult qualified Sony service personnel.
3. Connecting the unit to other devices may increase the leakage current.
4. For all peripheral devices connected to the unit that operate on commercial power supplies and do not comply with IEC 60601-1 standards, incorporate an isolation transformer that complies with IEC 60601-1 standards and connect to the commercial power supply via the transformer.
5. The unit generates, uses, and may radiate radio frequency energy. If it is not installed and used in accordance with the instruction manual, it may cause interference on other devices. If the unit causes interference (which can be determined by disconnecting the power cord from the unit), try the following.
 - Relocate the unit with respect to the affected devices.
 - Connect the unit and the affected devices to different branch circuits.For more information, consult qualified Sony service personnel.
(Applicable standard: IEC 60601-1-2)

Important EMC notices for use in medical environments

- The HVO-4000MT needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the instructions for use.
- The HVO-4000MT is intended for use in a professional healthcare facility environment.
- Portable and mobile RF communications equipment, such as cellular phones, can affect the HVO-4000MT.

Warning

- Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the HVO-4000MT. Otherwise, degradation of the performance of this equipment could result.
- If the HVO-4000MT will be used adjacent to or stacked with other equipment, normal operation of the HVO-4000MT under such configurations should be verified via observation.
- The use of accessories and cables other than those specified, with the exception of replacement parts sold by Sony Corporation, may result in increased emissions or decreased immunity of the HVO-4000MT.

Guidance and manufacturer's declaration – electromagnetic emissions		
The HVO-4000MT is intended for use in the electromagnetic environment specified below. The customer or the user of the HVO-4000MT should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The HVO-4000MT uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The HVO-4000MT is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Guidance and manufacturer's declaration – electromagnetic immunity

The HVO-4000MT is intended for use in the electromagnetic environment specified below. The customer or the user of the HVO-4000MT should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, a relative humidity of at least 30% is recommended.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U_T (100% dip in U_T) for 0.5/1 cycles ^a 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25/30 cycles ^a (for 0.5 sec) 0% U_T (100% dip in U_T) for 250/300 cycles ^a (for 5 sec)	0% U_T (100% dip in U_T) for 0.5/1 cycles ^a 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25/30 cycles ^a (for 0.5 sec) 0% U_T (100% dip in U_T) for 250/300 cycles ^a (for 5 sec)	Mains power quality should be that of a typical commercial or hospital environment. If the user of the HVO-4000MT requires continued operation during power mains interruptions, it is recommended that the HVO-4000MT be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: U_T is the a.c. mains voltage prior to application of the test level.

a For example, 10/12 means 10 cycles at 50 Hz or 12 cycles at 60 Hz.

Guidance and manufacturer's declaration – electromagnetic immunity

The HVO-4000MT is intended for use in the electromagnetic environment specified below. The customer or the user of the HVO-4000MT should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz outside ISM bands ^c	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the HVO-4000MT, including cables, than the recommended separation distance calculated from the equation appliance to the frequency of the transmitter. Recommended separation distance $d = 1.2 \sqrt{P}$
	6 Vrms 150 kHz to 80 MHz in ISM bands ^c	6 Vrms	
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.7 GHz	3 V/m	IEC 60601-1-2: 2007 $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz IEC 60601-1-2: 2014 $d = 2.0 \sqrt{P}$ 80 MHz to 2.7 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with following symbol: 

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

- a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the HVO-4000MT is used exceeds the applicable RF compliance level above, the HVO-4000MT should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the HVO-4000MT.
- b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.
- c The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

Recommended separation distances between portable and mobile RF communications equipment and the HVO-4000MT

The HVO-4000MT is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the HVO-4000MT can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the HVO-4000MT as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m				
	IEC 60601-1-2 : 2007			IEC 60601-1-2 : 2014	
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 2.7 GHz $d = 2.0 \sqrt{P}$
0.01	0.12	0.12	0.23	0.12	0.20
0.1	0.38	0.38	0.73	0.38	0.63
1	1.2	1.2	2.3	1.2	2.0
10	3.8	3.8	7.3	3.8	6.3
100	12	12	23	12	20

For transmitters rated a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Guidance and manufacturer's declaration – electromagnetic immunity

The HVO-4000MT is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the HVO-4000MT. Otherwise, degradation of the performance of this equipment could result.

Immunity test	Band ^a	Service ^a	Modulation	IEC 60601 test level	Compliance level
Proximity fields from RF wireless communications equipment IEC 61000-4-3	380 – 390 MHz	TETRA 400	Pulse modulation 18 Hz	27 V/m	27 V/m
	430 – 470 MHz	GMRS 460 FRS 460	FM ±5 kHz deviation 1 kHz sine	28 V/m	28 V/m
	704 – 787 MHz	LTE Band 13, 17	Pulse modulation 217 Hz	9 V/m	9 V/m
	800 – 960 MHz	GSM 800/900 TETRA 800 iDEN 820 CDMA 850 LTE Band 5	Pulse modulation 18 Hz	28 V/m	28 V/m
	1,700 – 1,990 MHz	GSM 1800 CDMA 1900 GSM 1900 DECT LTE Band 1, 3, 4, 25 UMTS	Pulse modulation 217 Hz	28 V/m	28 V/m
	2,400 – 2,570 MHz	Bluetooth WLAN 802. 11 b/g/n RFID 2450 LTE Band 7	Pulse modulation 217 Hz	28 V/m	28 V/m
	5,100 – 5,800 MHz	WLAN 802. 11 a/n	Pulse modulation 217 Hz	9 V/m	9 V/m

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a For some services, only the uplink frequencies are included.

Caution

When you dispose of the unit or accessories, you must obey the laws in the relative area or country and the regulations in the relative hospital regarding environmental pollution.



Warning on power connections

Use a proper power cord for your local power supply.

1. Use the approved Power Cord (3-core mains lead) / Appliance Connector / Plug with earthing-contacts that conforms to the safety regulations of each country if applicable.
2. Use the Power Cord (3-core mains lead) / Appliance Connector / Plug conforming to the proper ratings (Voltage, Ampere).

If you have questions on the use of the above Power Cord / Appliance Connector / Plug, please consult a qualified service personnel.



Warning on power connections for medical use

Customers in the U.S.A. and Canada should use the following type of power cord.

Customers in other countries or regions should use the power cord prescribed by their country or region.

	U.S.A. and Canada
Plug type	HOSPITAL GRADE*
Cord type	Min. Type SJT
	Min. 18 AWG
Minimum rating for plug and appliance couplers	10 A / 125 V
Safety approval	UL Listed and CSA

* Note: Grounding reliability can only be achieved when the equipment is connected to an equivalent receptacle marked "Hospital Only" or "Hospital Grade".

Caution

Danger of explosion if battery is incorrectly replaced. Replace only with the same or equivalent type recommended by the manufacturer.

When you dispose of the battery, you must obey the law in the relative area or country.

Do not dispose of the battery in a fire or a hot oven, or mechanically crush or cut the battery. It may explode or cause a fire.

Do not subject the battery to extremely low air pressure that may result in an explosion or the leakage of flammable liquid or gas.

Do not place the battery in a high temperature place, such as under direct sunlight or near fire. It may ignite, explode, or cause a fire.


Do not immerse or wet the battery in water or seawater. This may cause an electric shock.

Warning

When not using the remote for an extended period of time, remove the batteries to avoid possible damage from battery leakage and corrosion.

Warning

Batteries shall not be exposed to excessive heat such as sunshine, fire or the like.

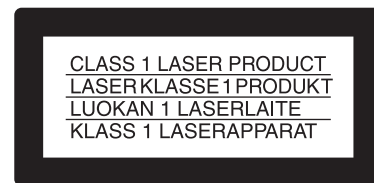

This model (HVO-4000MT) is classified as a CLASS 1 LASER PRODUCT. (IEC 60825-1: 2007 and IEC 60825-1: 2014)

Caution

Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

Caution

The use of optical instruments with this product will increase eye hazard.

Warning

The apparatus shall not be exposed to dripping or splashing. No objects filled with liquids, such as vases, shall be placed on the apparatus.

Caution

When installing, ensure the following space around the periphery of the unit, taking ventilation and servicing into consideration.

- Rear side: 10 cm (4 in.) or more
- Left/Right sides: 10 cm (4 in.) or more
- Top side: 10 cm (4 in.) or more



Caution

Do not come into contact with the unit's internal circuits and patients at the same time.

Doing so may result in a generation of voltage that can be harmful to patients if the unit is malfunctioning.



Caution

Do not use the device in a MR (Magnetic Resonance) environment.

It may cause a malfunction, fire, and unwanted movement.

Caution

This equipment is not suitable for use in locations where children are likely to be present.

Warning

Excessive sound pressure from earphones and headphones can cause hearing loss.

In order to use this product safely, avoid prolonged listening at excessive sound pressure levels.

For the State of California, USA only

Perchlorate Material - special handling may apply, See www.dtsc.ca.gov/hazardouswaste/perchlorate

For the customers in the U.S.A.

SONY LIMITED WARRANTY - Please visit <http://www.sony.com/psa/warranty> for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

For the customers in Canada

SONY LIMITED WARRANTY - Please visit <http://www.sonybiz.ca/pro/lang/en/ca/article/resources-warranty> for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

For the customers in Europe

Sony Professional Solutions Europe - Standard Warranty and Exceptions on Standard Warranty. Please visit <https://pro.sony/support-services/primesupport/support-professional-solutions-europe-standard-product-warranty> for important information and complete terms and conditions.

For the customers in Korea

SONY LIMITED WARRANTY - Please visit <http://bpeng.sony.co.kr/handler/BPAS-Start> for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

Table of Contents

Please Read First	14
Usage Notes.....	16

Chapter 1 Overview

Features	17
System Configuration Example	19
Names and Functions of Parts.....	20
Front	20
Rear.....	23
Infrared Remote Control Unit (RM-M010)	25
Status Displays	26
Setting Screens	28
Using the On-Screen Keyboard (Text Entry)	29
Names and Functions of Parts (On-Screen Keyboard)	29
Handling Discs	30
Notes on Handling	30
Inserting and Removing Discs.....	30

Chapter 2 Preparation

Connections.....	32
Turning the Unit On and Off	32
Configuring System Settings	33
Configuring Recording Settings (User Settings).....	34
Displaying the [User Settings] Screen.....	34
Configuring Input Signal Settings	34
Configuring Image Quality Settings.....	34
Configuring Recording Settings	35
Configuring Save Settings	36
Configuring Print Settings	37
Configuring Other Settings.....	38
Function Keys.....	38
Using Function Keys	38
Factory Assigned Functions	39
Other Assignable Functions	39

4K Convert Modes	40
------------------------	----

Chapter 3 Basic Recorder Operations

Operation Flow	42
Step 1: Record.....	43
Step 2: Capture Still Images	44
Step 3: Quick Playback	44
Step 4: End Recording	45

Chapter 4 Recording and Playback

Recording Preparation	46
Patient Information and Data Storage Settings	46
Using a Hard Keyboard for Patient Information and Data Storage Settings	47
Registering Multiple Patients in Advance.....	47
Registering Patients from MWL	48
Test Recording	49
Recording	50
Manual Recording	50
Capture Still Images	51
Manual Still Image Capture	51
Using the Ch1/Ch2 Simultaneous Recording Function	52
Configuring Settings for Ch1/Ch2 Simultaneous Recording	52
Performing Ch1/Ch2 Simultaneous Recording	52
Using the Ch1 2K Recording Function	52
Configuring Settings for Ch1 2K Recording	52
Performing Ch1 2K Recording	52
Limitations during Playback and Output	52
Playback	53
Playing Back the Most Recent Data (Quick Playback).....	53

Image Search	54
Specifying Search Conditions	54
Viewing Thumbnails of Recorded Data	54
Playback from the Search Results List	55
Sorting the Recorded Data List	55
Processing Recorded Data.....	56
Printing	56
Copying to External Media	57
Editing Patient Information	58
Protecting Recorded Data.....	59
Deleting Recorded Data	59
Viewing Information on Recorded Data	59

Chapter 5 System Administrator Settings

Displaying the [System Admin Settings] Screen	61
Language and Time Settings	62
List of time zones	62
Function Settings	64
[Patient Info.] Tab	64
[Input Detection] Tab	64
[Input Signal] Tab.....	64
[Auto Live] Tab.....	65
[General 1] Tab.....	65
[General 2] Tab.....	66
[Facility] Tab	66
[CMS] Tab.....	67
Device Settings.....	67
[Device 1] Tab.....	67
[Device 2] Tab	69
[Control] Tab	69
[Contact Switch] Tab.....	70
Password Settings	71
Network Settings	72
[Network] Tab	72
[IP Address] Tab.....	72
[DNS Server] Tab.....	72
[File Server] Tab.....	72
[Shared] Tab	73
[NTP] Tab.....	73

Doctor List Registration	73
Registering Doctors.....	73
[General] Tab	73
[File Server] Tab	74
[Shared] Tab.....	74
[Print] Tab	74
[External Media] Tab	74
Editing the Doctor List	75
Sorting the Doctor List.....	75
Editing Doctor Settings	75
Deleting Doctors	75
Editing Cases	76
Registering Cases.....	76
Deleting Cases.....	76
Sorting the Case List.....	76
Editing Case Names	76
Registering Categories	76
Selecting Categories.....	77
Editing Categories	77
Sorting the Category List	77
Editing Category Names	77
Auto Delete Settings.....	78
DICOM Settings.....	78
[MWM Server 1] Tab.....	78
[MWM Server 2] Tab.....	79
[C-Store] Tab	79
[Local Station 1] Tab	80
[Local Station 2] Tab	80

Chapter 6 Touch Panel and Mouse

Overview	81
Configuring Touch Panel / Mouse Settings.....	81
Screen Displays in Touch Panel / Mouse Mode	82
[MENU] Screen	82
[Status] Screen	83
Video Playback Screen	84
Still Image Playback Screen.....	84

Chapter 7 Miscellaneous

Error Messages	85
Troubleshooting	87
Licenses	88
DCMTK.....	88
zlib	93
PostgreSQL.....	93
expat	93
freetype	94
memcached.....	95
openldap	96
OpenSSL.....	96
proj.....	98
zlib	98
gdal, libxm2, libxslt.....	98
libedit, libevent, libmemcached, tiff.....	99
wxwidget	100
ncurses	100
krb5.....	101
curl.....	101
cyrus-sasl	101
libjpeg	102
LibTIFF	102
GNU GPL/LGPL Licensed Software	102

GNU GENERAL PUBLIC LICENSE	103
Preamble.....	103
TERMS AND CONDITIONS FOR COPYING, DISTRIBUTION AND MODIFICATION	103
END OF TERMS AND CONDITIONS	106
How to Apply These Terms to Your New Programs	106
GNU LESSER GENERAL PUBLIC LICENSE	107
Preamble.....	107
TERMS AND CONDITIONS FOR COPYING, DISTRIBUTION AND MODIFICATION	108
END OF TERMS AND CONDITIONS	112
How to Apply These Terms to Your New Libraries	112
Specifications.....	114
Index.....	116

Trademarks

- Blu-ray Disc™, Blu-ray™, and their logos are trademarks of Blu-ray Disc Association.
 - The terms HDMI and HDMI High-Definition Multimedia Interface, and the HDMI Logo are trademarks or registered trademarks of HDMI Licensing Administrator, Inc. in the United States and other countries.
 - The products or system names appearing in this document are trademarks or registered trademarks of their respective owners.
- Further, the © or ™ symbols are not used in the text.

- Reproduction or duplication, in whole or part, of the software or operation manual supplied with the recorder, as well as renting or leasing of the software without the authorization of the right holder is prohibited under copyright law.
 - Sony assumes no responsibility for damages, loss of income, or any claims from a third party arising out of use of the recorder or supplied software.
 - For complete terms and conditions of the warranty for the recorder, refer to the warranty card included in the package.
 - The software supplied with the recorder cannot be used with any other recorders.
 - It is not possible to install any software into the equipment other than the software supplied by Sony specifically for use with the equipment.
 - Note that the specifications of the recorder and supplied software are subject to change for improvement without prior notice.

Please Read First

Copyright

Using this unit for video and/or audio recording, or distribution over the network or otherwise may in some cases require the permission of the copyright holder of the video or audio. To protect copyright, observe the following points carefully when using this unit.

- When connecting a recording device to this and recording video or audio, carefully observe laws relating to copyright.
- Without the permission of the copyright holder, the showing or distribution of video or audio material of which the copyright is held by a third party, or the act of recording on the hard disk of this unit, sharing folders, and permitting of access to a private group or to the public is prohibited by law.
- With a software upgrade or functional extension, with the object of protecting copyright, the specifications for the video and audio signals that can be input may be changed without notice.
- Under copyright law, you may not use recorded video or audio other than for your personal enjoyment without the permission of the copyright holder. Note that at live performances, shows and exhibitions, even for your personal entertainment shooting may be restricted.

Disclaimer of Liability for Recorded Data

Sony Corporation does not accept any liability whatsoever for any problems arising from a failure to record, or from damage or erasure of recorded content on this equipment, for any reason. This includes claims for compensation of recorded content, and for any concomitant and consequential damages. Sony Corporation will not repair, restore, or duplicate any recorded content. Your use of this product is subject to these conditions.

Problems may occur if you perform the following operations:

- When you use a disc recorded using this unit with another BD (Blu-ray Disc)/DVD recorder or the BD/DVD drive of a personal computer. Discs recorded using this unit will not play back on other BD/DVD recorders.
- When you attempt to re-use a disc with this unit, after using it in the operations described above.
- When you use a disc recorded using another BD/DVD recorder or BD/DVD drive of a personal computer with this unit. Discs recorded using other BD/DVD recorders or computers will not play back on this unit.

Deletion of data when disposing of the unit

Before disposing of the unit, you first must delete patient data and other confidential information that is stored therein.

Failure to delete stored data before disposing of the unit could expose confidential information to third parties. If you have any questions in this regard, contact your Sony dealer.

On security

SONY WILL NOT BE LIABLE FOR DAMAGES OF ANY KIND RESULTING FROM A FAILURE TO IMPLEMENT PROPER SECURITY MEASURES ON TRANSMISSION DEVICES, UNAVOIDABLE DATA LEAKS RESULTING FROM TRANSMISSION SPECIFICATIONS, OR SECURITY PROBLEMS OF ANY KIND.

Depending on the operating environment, unauthorized third parties on the network may be able to access the unit. When connecting the unit to the network, be sure to confirm that the network is protected securely.

On condensation

If the unit is suddenly taken from a cold to a warm location, or if ambient temperature suddenly rises, moisture may form on the outer surface of the unit and/or inside of the unit. This is known as condensation. If condensation occurs, turn off the unit and wait until the condensation clears before operating the unit. Operating the unit while condensation is present may damage the unit.

LCD panel

The LCD panel fitted to this unit is manufactured with high precision technology, giving a functioning pixel ratio of at least 99.99%. Thus a very small proportion of pixels may be “stuck”, either always off (black), always on (red, green, or blue), or flashing. In addition, over a long period of use, because of the physical characteristics of the liquid crystal display, such “stuck” pixels may appear spontaneously. These problems are not a malfunction. Note that any such problems have no effect on recorded data.

LCD image display

Due the physical characteristics of LCD panels, there may be a decrease in brightness or change in color temperature over a long period of use. These problems are not a malfunction.

In addition, these occurrences will not affect recorded data.

On consumable parts

- The HDD, fan, battery, and BD/DVD drive are consumable parts that will need periodic replacement. When operating at room temperature, a normal replacement cycle will be about 2 to 3 years. However, this replacement cycle represents only a general guideline and does not imply that the life expectancy of these parts is guaranteed. For details on parts replacement, contact your dealer.
- The life expectancy of the electrolytic capacitor is about 5 years under normal operating temperatures and normal usage (8 hours per day; 25 days per month). If usage exceeds the above normal usage frequency, the life expectancy may be reduced correspondingly.

Precautions for products with built-in HDD

This unit has a built-in hard disk drive (HDD). The HDD is a precision device. If subject to shock, vibration, static electricity, high temperature or humidity, data loss can occur. When installing and using the unit, closely observe the following precautions.

Protect from shocks and vibrations

When subject to shocks or vibrations, the HDD can be damaged and loss of data on the HDD can occur.

- When transporting the unit, use the specified packing material. When transporting on a dolly or similar, use a type which does not transmit excessive vibrations. Excessive shocks and vibrations can damage the HDD.
- Never move the unit while it is powered.
- Do not remove panels or outer parts of the unit.
- When placing the unit on a floor or other surface, make sure that the unit is equipped with the specified rubber feet, and put the unit down carefully. If there are no feet, mount the rubber feet first.
- Do not place the unit near other devices that may become a source of vibrations.

Wait for 30 seconds after turning power off

For a brief interval after the power is turned off, the platters inside the HDD will still keep spinning and the heads will be in an insecure position. During this interval, the unit is more susceptible to shocks and vibrations than during normal operation. For a period of at least 30 seconds after turning power off, avoid subjecting the unit even to very light shocks. After this period, the hard disk will be fully stopped and the unit can be manipulated.

Temperature and humidity related precautions

Use and store the unit only in locations where the specified temperature and humidity ranges are not exceeded. (Be sure to use the unit that conforms fully to the specifications of this unit.)

When the HDD seems to be faulty

Even if the HDD is showing signs of malfunction, be sure to observe all the above precautions. This will prevent further damage from occurring until the problem can be diagnosed and corrected.

Replacement of the HDD and other consumable parts

The HDD and battery are consumable parts that will need periodic replacement. When operating at room temperature, a normal replacement cycle will be about two to three years. However, this replacement cycle represents only a general guideline and does not imply that the life expectancy of these parts is guaranteed. For details on parts replacement, contact your dealer.

Notes on media

Operation with all optical media or USB media is not guaranteed. Contact your dealer regarding the type of media you intend to use.

Use with electrosurgical knives and similar devices

If this unit is used together with an electrosurgical knife, etc., the picture may be disturbed, warped or otherwise abnormal as a result of strong radio waves or voltages from the device. This is not a malfunction.

When you use this unit simultaneously with a device from which strong radio waves or voltages are emitted, confirm the effect of this before using such devices, and install this unit in a way that minimizes the effect of radio wave interference.

Usage Notes

Do not block the vents

When installing, ensure the following space around the periphery of the unit, taking ventilation and servicing into consideration.

- Rear side: 10 cm (4 in.) or more
- Left/Right sides: 10 cm (4 in.) or more
- Top side: 10 cm (4 in.) or more

Do not subject the unit to severe shocks

The internal mechanism may be damaged or the body warped.

Do not cover the unit while operating

Doing so will cause temperatures to rise inside the unit, possibly resulting in failure.

After use

Press the  (on/standby) switch to enter standby mode. If you do not intend to use the unit for an extended period, set the main power switch to the  (off) position after entering standby mode, and then disconnect the power cord.

Shipping

- Remove the BD/DVD disc before transporting the unit.
- If sending the unit by truck, ship, air or other transportation service, pack it in the shipping carton of the unit.

Care of the unit

If the casing, or panel becomes dirty, wipe them gently with a soft, dry cloth.

For stubborn dirt, use a cloth wet with mild liquid detergent to wipe the unit, and then wipe it with a dry cloth. Using alcohol, thinner, benzine, insecticides, or other volatile substances may damage the surface of the unit or remove the finish.

Use isopropyl alcohol with a concentration of 50% to 70% (v/v) or ethanol with a concentration of 76.9% to 81.4% (v/v) to clean the surface of the unit.

In the event of operating problems

If you should experience problems with the unit, contact your Sony dealer.

Features

The Sony HVO-4000MT has video input and output terminals, built-in HDD. It is possible to record input video signal and playback the recording data to outputs terminal.

This system has Blu-ray, DVD drive and user can move the recording data to Blu-ray and DVD-R disc.

Improved efficiency during and after medical procedures

Video and still image recording

Record videos and still images from endoscopic devices and cameras used in the operating room. 4K compatibility allows recording and playback in high definition.

Ch1/Ch2 simultaneous recording

The unit is equipped with a Ch1/Ch2 simultaneous recording function that automatically performs Ch2 recording (2K) during Ch1 recording (4K).

Resolution switching (4K/2K) for Ch1 recording

The recording resolution for Ch1 can be switched between 4K and 2K without restarting the unit. When [Ch1 2K Rec] is configured, reception of 2K (1920 × 1080) input signals is enabled. In such cases, Ch2 recording will be disabled.

Simultaneous saving to external storage media

Data can be saved to external media while it is being recorded to the internal hard disk.

Obtaining patient information from MWL

Patient IDs, patient names, and other patient information can be obtained from Modality Worklist (MWL).

Printer connectivity

Connect the unit to a printer, and print images via simple controls from the [Image List] screen. You can also configure settings to print images as they are captured via a foot switch, for example. The printed images can be used in simple medical reports.

Support for a variety of external storage media

Record simultaneously to USB-compatible external hard drives and USB memory devices via simple controls, or copy the data at a later time. This is useful when you want to edit recorded images on a computer.



Warning

Using this unit for medical purposes

The connectors on this unit are not isolated.

Always use USB-compatible external hard drives and USB memory devices that operate solely on the power supplied by the unit.

Connecting a device that operates on an AC power supply may result in an influx of leakage currents from the connected device, which may in turn result in electric shocks to the patient and operator.

If use of such devices cannot be avoided, be sure to connect an isolation transformer to the power supply of the connected device or connect an isolator between the connection cables to isolate the power supply, and verify that the risk of electric shock from the system has been sufficiently reduced to levels that conform to IEC 60601-1 standards.

Editing recorded data

Record images to external media, and use commercially available software to play back or edit the images on a computer as necessary.

For details on playback and editing software, contact your local Sony representative.

Network file transfer

Transfer recorded data from the operating room directly to computers within the hospital network.



Warning

Using this unit for medical purposes


The connectors on this unit are not isolated.

Do not connect any device other than one which conforms to IEC 60601-1 standards.

When an information technology device or AV device that uses an alternating current is connected, current leakage may result in an electric shock to the patient or operator.

If use of such a device is unavoidable, isolate its power supply by connecting an isolation transformer, or by connecting an isolator between the connecting cables. After implementing these measures, confirm that the reduced risk now conforms to IEC 60601-1 standards.

System Configuration Example

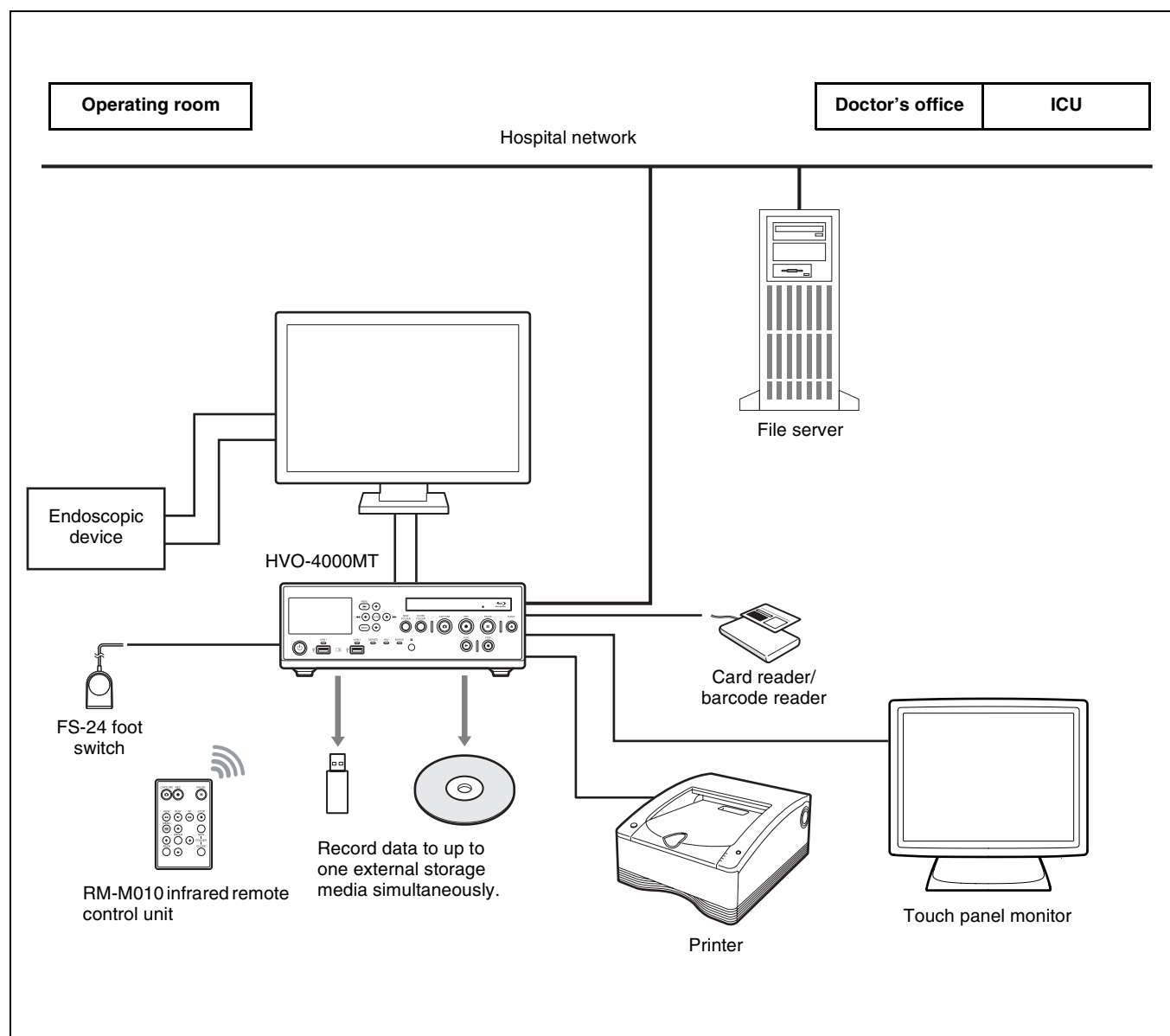


Caution

The FS-24 has an Ingress Protection rating of IPX3. Therefore, do not operate it environments exposed to splashing liquids (e.g., surgical operating rooms).

For safety, use a device with a rating of IPX6 or higher when operating in such environments.

Names and Functions of Parts

Front

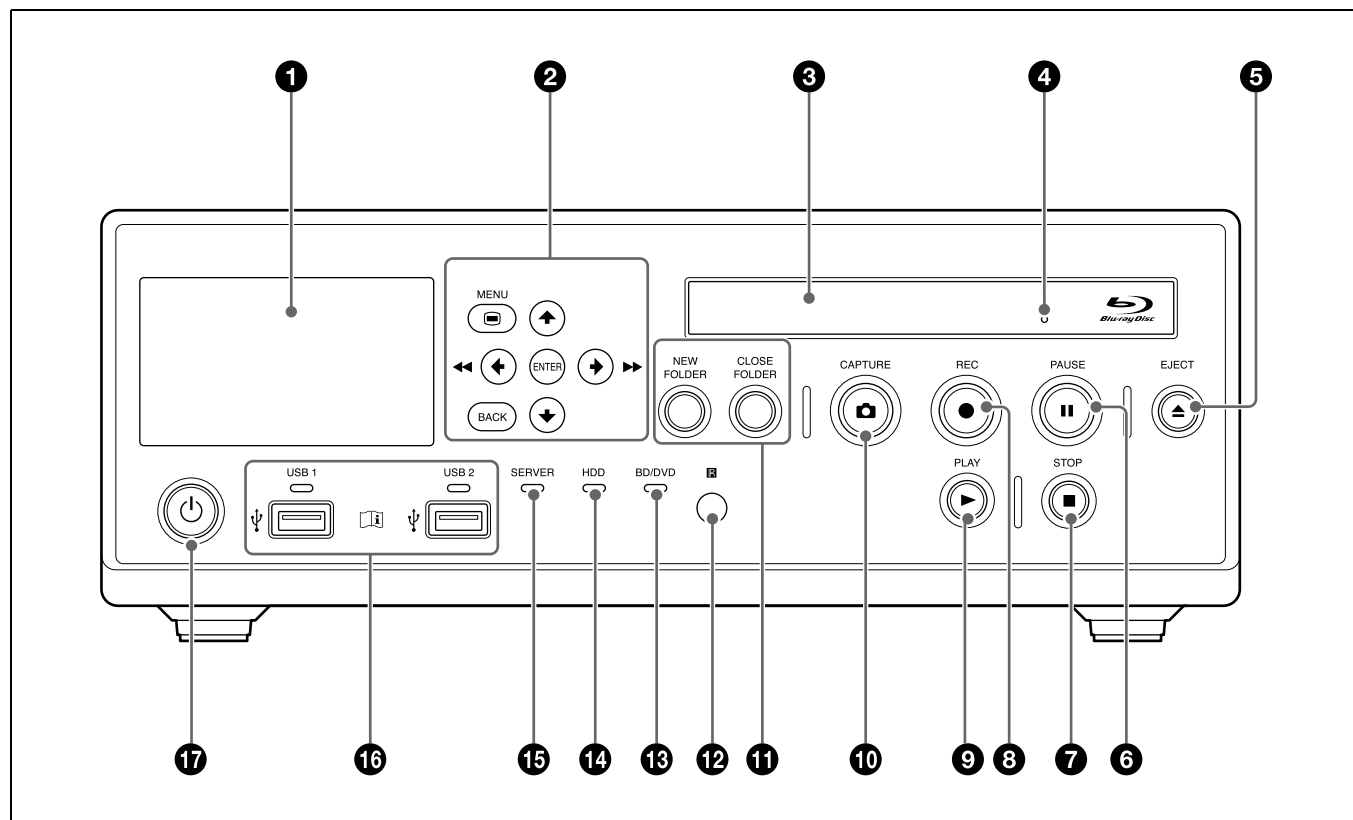

Note

Do not touch the laser pickup inside the BD drive

The discharge of static electricity that may result from touching the pickup may cause the drive to hang up, preventing the ability to write data. If a hang up occurs, try turning the recorder off and turning it on again.

1 Front panel display

Displays Ch1 input/output videos and still images, the status of the unit, error messages, etc.

2 Menu navigation buttons

Use these to navigate menus.

MENU button

Display or hide the menu.

button

Use this to move the cursor up or select an item that is above the current item.

button

Use this to move the cursor left or select an item that is to the left of the current item.

You can also use this button to rewind during playback of recorded data.

button

Use this to move the cursor down or select an item that is below the current item.

button

Use this to move the cursor right or select an item that is to the right of the current item.

You can also use this button to fast forward during playback of recorded data.

ENTER button

Use this to confirm a currently selected menu or item, and execute operations.

BACK button

Use this to return to a previous screen during menu or operation screen displays.

When you press this button while nothing is displayed on the screen, the text that appears on the front panel display will be hidden.

3 Disc tray

Insert discs here.

4 Emergency eject hole

If the disc tray does not open when you press the  EJECT button, insert a long, thin object such as a straightened paper clip into this hole to open the disc tray.

5 EJECT button

Use this to open the disc tray.

6 PAUSE button

Use this while recording is in progress to pause recording. Press this button again to resume recording. A chapter separator is also created at the position in which recording was paused. You can also use this button to pause playback. To resume playback, press this button again or press the  PLAY button. “PAUSE” appears on the front panel display during pause.

7 STOP button

Use this to stop recording or playback. A title separator is also created at the position in which recording was stopped.

8 REC button

Use this to start video/audio recording. During recording, “REC” appears on the front panel display in addition to the recording time.

For details on use, see “Manual Recording” (page 50).

9 PLAY button

- Use this to play back the most recent recorded data stored on the internal hard disk, or to play back the recorded data selected in the [Recording List] or [Image List] screen.
- Press this button while playback is paused to resume playback.

10 CAPTURE button

Use this to capture still images. A chapter separator is also created at the position in which the still image was captured.

For details on use, see “Capture Still Images” (page 51).

11 NEW FOLDER button

Use this to display the [Recording Preparation] screen to start recording a new surgical procedure or examination.

For details on use, see “Recording Preparation” (page 46).

CLOSE FOLDER button

Use this to display the [End Recording] screen which is used to end surgical procedure or examination recordings.

For details on use, see “Step 4: End Recording” (page 45).

12 Infrared receiver

Point the infrared remote control unit toward this receiver.

13 BD/DVD indicator

Displays the status of the BD/DVD drive.

Indicator	Status
Blinking green	Data is being written to the BD/DVD disc.
Lit green	The recording area of the BD/DVD disc is full. Insert a different disc.
Lit orange	An error has occurred.

Note

You cannot remove the disc when this indicator is blinking green.

14 HDD indicator

Displays the status of the internal hard disk.

Indicator	Status
Blinking green	Data is being recorded to the internal hard disk.
Lit orange	An error has occurred. The recording area of the internal hard disk is full. Delete any unnecessary recorded data.

15 SERVER indicator

Displays the server access status.

Indicator	Status
Blinking green	The server is being accessed.
Lit orange	An error has occurred. The recording area of the server is full.

16 USB ports/indicators 1 and 2

Connect USB memory devices, external hard disks, and other USB media here.

The indicators display the status of USB media access.

Indicator	Status
Blinking green	Data is being written to the USB media.

Indicator	Status
Lit green	The recording area of the USB media is full. Delete any unnecessary recorded data, or insert different media.
Lit orange	An error has occurred.

Note

Never remove USB media when these indicators are blinking green.

Tip

Supported USB devices

- Use Sony USB media.
- The connectors do not support all USB devices.
- USB hubs and devices with built-in hubs are not supported.



Warning

Using this unit for medical purposes

The connectors on this unit are not isolated. Always use USB-compatible external hard drives and USB memory devices that operate solely on the power supplied by the unit.

Connecting a device that operates on an AC power supply may result in an influx of leakage currents from the connected device, which may in turn result in electric shocks to the patient and operator.

If use of such devices cannot be avoided, be sure to connect an isolation transformer to the power supply of the connected device or connect an isolator between the connection cables to isolate the power supply, and verify that the risk of electric shock from the system has been sufficiently reduced to levels that conform to IEC 60601-1 standards.

17 (on/standby) switch/indicator

Use this switch to turn on the unit and the front panel display. Press this switch again to turn off the unit and enter standby mode.


The indicator indicates the following.

Indicator	Status
Off	The main power switch on the rear panel is set to the ○ (off) position.
Lit orange	Standby mode. (The main power switch on the rear panel is set to the I (on) position.)
Lit green	The unit is turned on. (The  (on/standby) switch was pressed to turn on the unit.)

Tip

The indicator will turn off about 15 seconds after the main power switch on the rear panel is turned off from the standby state.

Rear

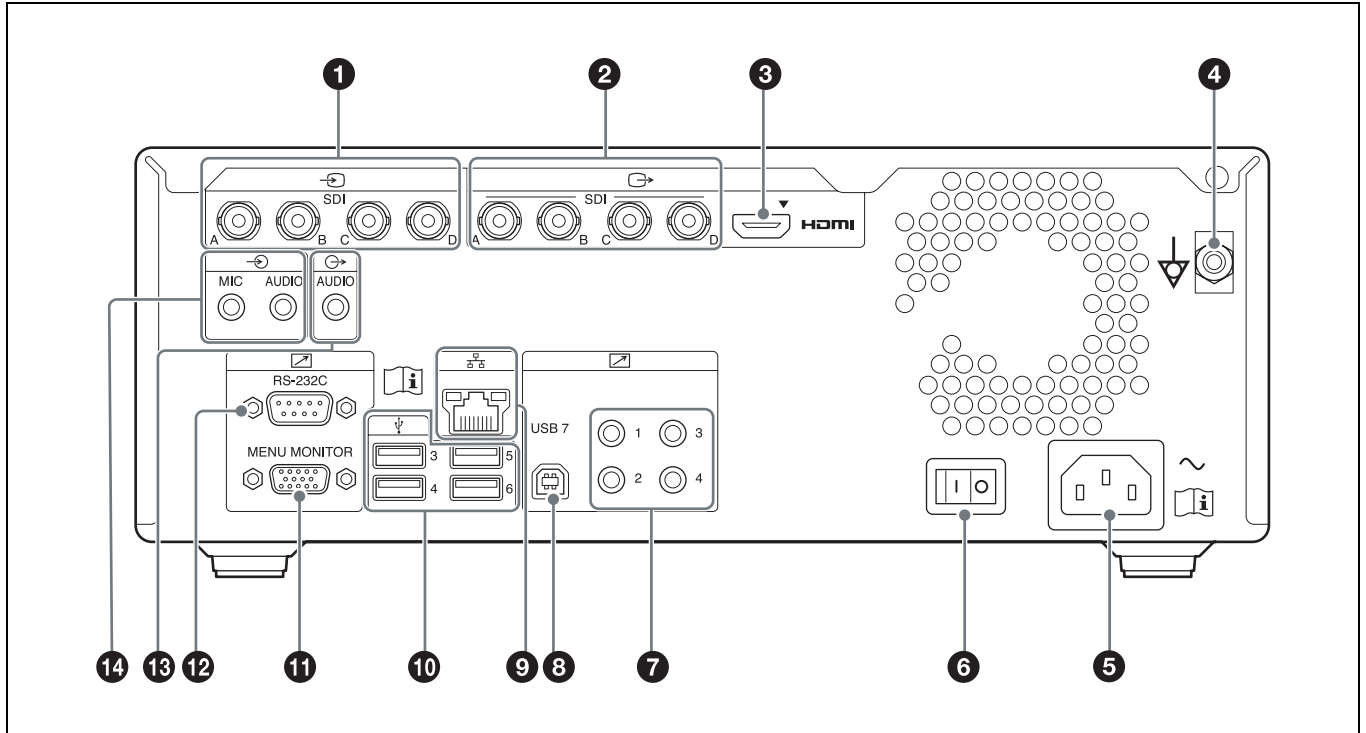

Warning

Using this unit for medical purposes

The connectors on this unit are not isolated.
Do not connect any device other than one which conforms to IEC 60601-1 standards.

When an information technology device or AV device that uses an alternating current is connected, current leakage may result in an electric shock to the patient or operator.

If use of such a device is unavoidable, isolate its power supply by connecting an isolation transformer, or by connecting an isolator between the connecting cables. After implementing these measures, confirm that the reduced risk now conforms to IEC 60601-1 standards.



Caution

Do not come into contact with the terminals of the rear panel connectors and patients at the same time. Doing so may result in a generation of voltage that can be harmful to patients if the unit is malfunctioning. Always disconnect the power cord before connecting and disconnecting connectors.

1 SDI input connectors A to D (BNC type)

Inputs SDI signals (3G SDI).

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [User Settings] screen - [Recording] tab, SDI signals (3G SDI) are only input from SDI input connector C.

Tips

- If SDI input images are not displayed, turn on the unit first and make sure that startup is complete before turning on the video input device.
- When receiving 4K signals, the four SDI inputs must be synchronized.
- Embedded audio is not supported.

2 SDI output connectors A to D (BNC type)

Outputs SDI signals (3G SDI).

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [User Settings] screen - [Recording] tab, SDI signals (3G SDI) are only output from SDI output connector C.

3 HDMI output connector (Type A)

Outputs HDMI signals.

Note

We recommend using a Sony HDMI cable (not supplied).

- 4  **Equipotential grounding terminal**
Use this for equipotential grounding connections.

- 5  **AC IN connector**
Connect a power cord (not supplied) here.

- 6 **I (on)/O (off) main power switch**
Set this to the **I (on)** position to turn on the power. To turn off the power, set the switch to the **O (off)** position.
During regular use of the unit, leave the main power switch in the **I (on)** position, and use the 17  (on/standby) switch on the front panel of the unit to switch between operational mode and standby mode.

Note

When you press the 17  (on/standby) switch on the front panel of the unit while the unit is in operational mode, data is saved before the unit enters standby mode. Always make sure that the unit is in standby mode before setting the main power switch to the **O (off)** position to turn off the power.

- 7  **REMOTE contact switch connectors 1 to 4 (stereo mini jack)**

Use these to control the unit via contact switches.

For details on contact switches, see “[Contact Switch] Tab” (page 70).

- 8  **USB port 7 (Type B / USB 2.0)**
Use this to perform system controls from an external device.
For details on system controls, see “[Control] Tab” (page 69).

- 9  **Network connector (RJ-45)**
Connect a 1000 Base-T or 100 Base-TX network cable here.

Caution

For safety, do not connect the connector to peripheral device wiring that might have excessive voltage. Follow the instructions for use for this port.

- 10  **USB ports 3 to 6 (Type A / USB 2.0)**
Connect USB devices here.

Tip

Supported USB devices

- Use Sony USB media.
- The connectors do not support all USB devices.
- USB hubs and devices with built-in hubs are not supported.



Warning

Using this unit for medical purposes

The connectors on this unit are not isolated. Always use USB-compatible external hard drives and USB memory devices that operate solely on the power supplied by the unit.

Connecting a device that operates on an AC power supply may result in an influx of leakage currents from the connected device, which may in turn result in electric shocks to the patient and operator.

If use of such devices cannot be avoided, be sure to connect an isolation transformer to the power supply of the connected device or connect an isolator between the connection cables to isolate the power supply, and verify that the risk of electric shock from the system has been sufficiently reduced to levels that conform to IEC 60601-1 standards.

- 11  **MENU MONITOR connector (15-pin mini D-sub)**

Connect a touch panel monitor or similar device here.

- 12  **RS-232C connector (9-pin D-sub)**

Use this to perform system controls from an external device.

For details on system controls, see “[Control] Tab” (page 69).

- 13  **AUDIO (analog audio signal) output connector (stereo mini jack)**

Outputs analog audio signals.


- 14  **MIC input connector (stereo mini jack)**

Inputs analog audio signals from the microphone. You can configure settings to supply power from the unit.

-  **AUDIO (analog audio signal) input connector (LINE) (stereo mini jack)**

Inputs analog audio signals.

Infrared Remote Control Unit (RM-M010)

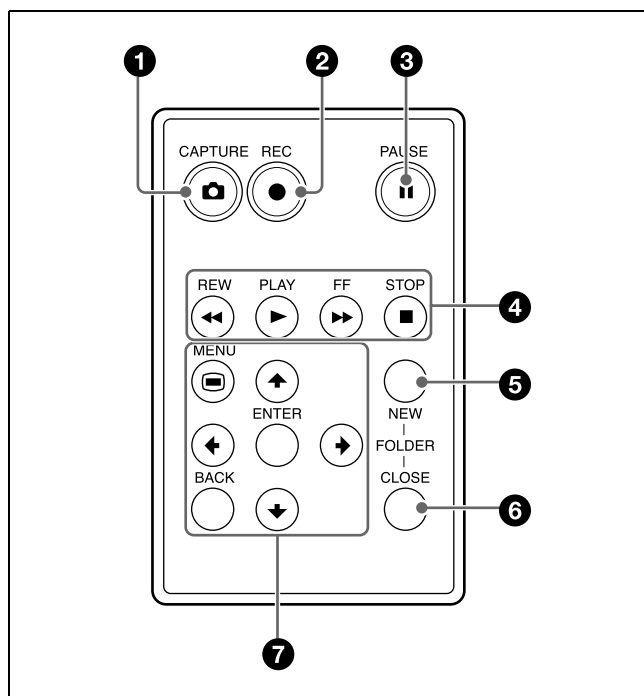

1 **CAPTURE button**

Use this to capture still images.

A chapter separator is also created at the position in which the still image was captured.

For details on use, see “Capture Still Images” (page 51).

2 **REC button**

Use this to start video/audio recording.

During recording, “REC” appears on the front panel display in addition to the recording time.

For details on use, see “Manual Recording” (page 50).

3 **PAUSE button**

Use this while recording is in progress to pause recording. Press this button again to resume recording.

A chapter separator is also created at the position in which recording was paused.

You can also use this button to pause playback.

To resume playback, press this button again or press the ► **PLAY** button.

“PAUSE” appears on the front panel display during pause.

4 **REW button**

Use this to rewind during playback of recorded data.

► **PLAY button**

- Use this to play back the most recent recorded data stored on the internal hard disk, or to play back the recorded data selected in the [Recording List] or [Image List] screen.
- Press this button while playback is paused to resume playback.

►► **FF button**

Use this to fast forward during playback of recorded data.

■ **STOP button**

Use this to stop recording or playback.

A title separator is also created at the position in which recording was stopped.

5 **NEW FOLDER button**

Use this to display the [Recording Preparation] screen to start recording a new surgical procedure or examination.

For details on use, see “Recording Preparation” (page 46).

6 **CLOSE FOLDER button**

Use this to display the [End Recording] screen which is used to end surgical procedure or examination recordings.

For details on use, see “Step 4: End Recording” (page 45).

7 **Menu navigation buttons**

Use these to navigate menus.

☰ **MENU button**

Display or hide the menu.

↓ **button**

Use this to move the cursor down or select an item that is below the current item.

→ **button**

Use this to move the cursor right or select an item that is to the right of the current item.

ENTER button

Use this to confirm a currently selected menu or item, and execute operations.

← **button**

Use this to move the cursor left or select an item that is to the left of the current item.

↑ **button**


Use this to move the cursor up or select an item that is above the current item.

BACK button

Use this to return to a previous screen during menu or operation screen displays.
When you press this button while nothing is displayed on the screen, the text that appears on the front panel display will be hidden.

Before using the remote control unit


Remove the insulation sheet.

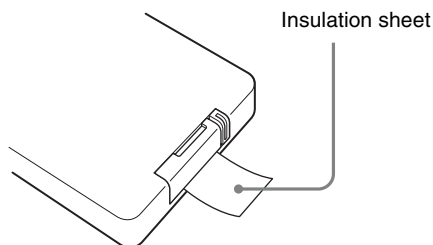


To replace the lithium battery

Use a Sony or Murata Manufacturing CR2025 lithium battery. Do not use any other type of battery with the remote control unit.

- 1 While holding in the lock lever (①), pull out the battery holder (②).

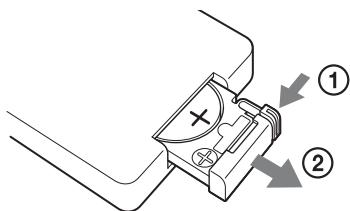

- 2 Insert the replacement battery with the positive (+) side facing up (①), and reinsert the battery holder until it clicks into place (②).

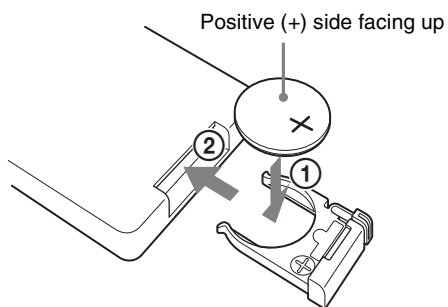

Warning

Batteries shall not be exposed to excessive heat such as sunshine, fire or the like.

Caution

Danger of explosion if battery is incorrectly replaced.
Replace only with the same or equivalent type recommended by the manufacturer.
When you dispose of the battery, you must obey the law in the relative area or country.

Battery lifetime

When the lithium battery output falls, even button presses may not operate. The average lithium battery lifetime is about one year, but this depends on the pattern of use. If pressing the remote control buttons produces absolutely no effect on this unit, replace the battery, then check the operation again.


Status Displays

The following information is displayed when the status display settings are enabled.


The status display settings are configured in the [System Admin Settings] screen – [Function Settings] screen – [General 1] tab. For details on this setting, see “[General 1] Tab” (page 65).

Rec status display / idle status display

Rec status displays

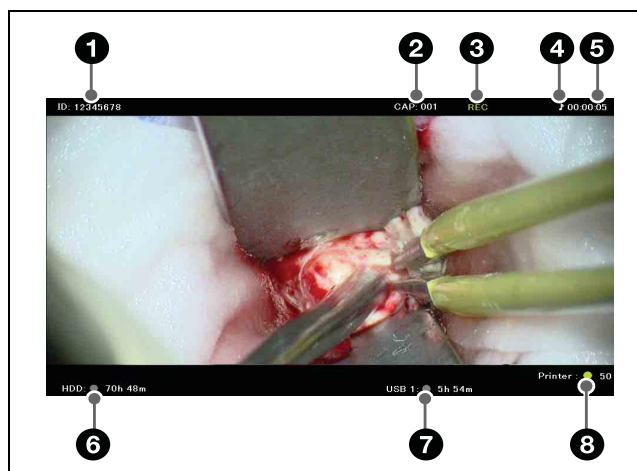

Idle status displays

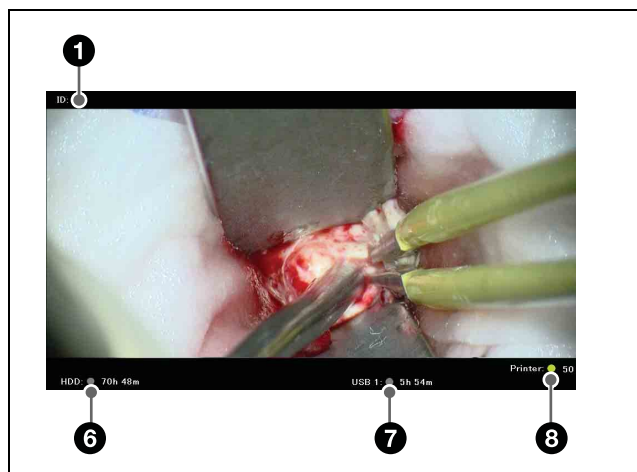

- 1 **Patient ID**
Displays the patient ID.

② Capture counter

Displays the number of still images captured.

③ Recording status

Displays the recording status.

REC: Recording is in progress.

REC PAUSE: Recording has been paused.

REC STOP: Recording preparation is complete or recording has been stopped.

④ Audio recording status

Displays the audio recording status.

♪: Audio recording is in progress.

no display: Audio is not being recorded.

⑤ Recording time

The recording time is counted in seconds and displayed in “HH:MM:SS” a format.

⑥ HDD status

Displays the internal HDD status and the remaining time.

The status color indicates the following.

Green: Recording is in progress or still image recording is in progress.

Orange: The maximum number of recorded data entries has been exceeded or the internal HDD is full.

Gray: Status other than the above.

⑦ External media status

Displays the status of the external media and the remaining capacity time.

The status color indicates the following.

- When transferring to BD/DVD

Green: Transfer or formatting is in progress.

Orange: Additional recording to the BD/DVD is not possible, an unsupported BD/DVD is inserted, the maximum number of BD/DVD exchanges has been exceeded, or a finalized DVD is inserted.

Gray: Status other than the above.

- When transferring to USB

Green: Transfer is in progress.

Orange: An unsupported USB memory device is inserted, the maximum number of USB memory device exchanges has been exceeded, or a transfer error has occurred.

Gray: Status other than the above.

- When transferring to a server

Green: Transfer is in progress.

Orange: Insufficient server capacity or a transfer error has occurred.

Gray: Status other than the above.

⑧ Print status

Displays the printer status and the remaining number of print sheets.


The status color indicates the following.

Green: Still image printing is in progress.

Orange: A print error has occurred.

Gray: Status other than the above.

Playback status displays

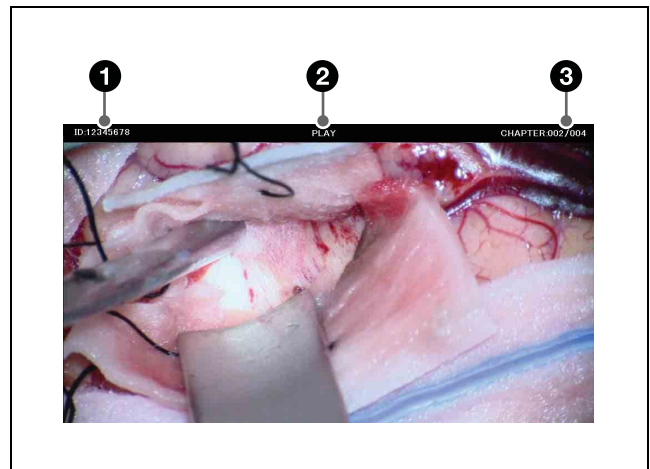

① Patient ID

② Playback status

Displays the playback status.

PLAY: Video playback is in progress.

FF x2 to FF x60: Fast forward at x2 to x60 speed is in progress.

REW x2 to REW x60: Rewind at x2 to x60 speed is in progress.

PAUSE: Playback has been paused.


③ Video chapter number / still image playback status

CHAPTER:xxx/xxx: Displays the number of the chapter being played back.

STILL:xxx/xxx: Displays the number of the still image being played back.

Setting Screens

Example: [User Settings] screen

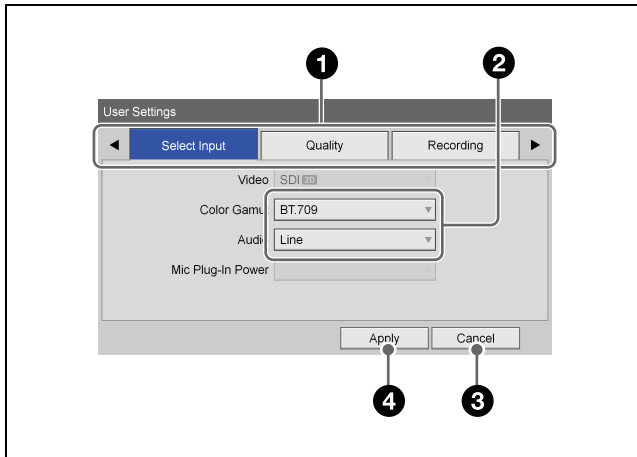

1 Tabs

The setting screens consist of tabs.

To switch between tabs, use the ◀ and ▶ buttons on the front panel or remote control to select the desired tab.

If additional tabs exist before and after the displayed tab, ◀ (previous) and ▶ (next) are also displayed.

The tabs will be displayed in sequence when you select ◀ and ▶.


2 Setting item

Use the ↑, ↓, ←, and → buttons to select an item, and press the ENTER button to display the setting screen for that item. You can select or enter setting values in the screens that appear.

Example: When selecting a setting value

Select the item, and press the ENTER button.


The previous screen will appear again, and the setting value will be entered.

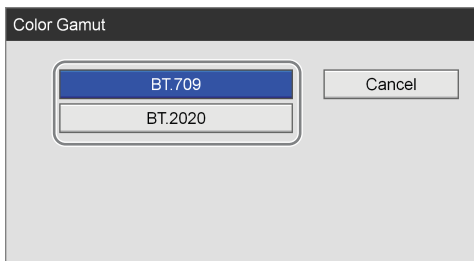

Example: When selecting check boxes

Select the item and press the ENTER button to select the check box. The check box will be selected/cleared with each press of the ENTER button. (Multiple selections can be made.)

When you select [OK] and press the ENTER button, the previous screen appears again with the setting value entered.

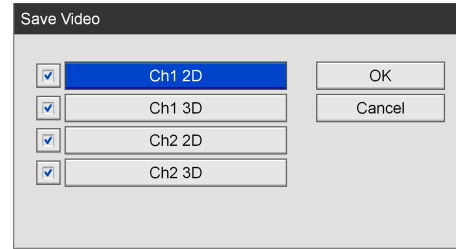

3 [Cancel]

Select this button and press the ENTER button to cancel the settings and close the setting screen.

4 [Apply]

Select this button and press the ENTER button to save the settings and close the setting screen.


Using the On-Screen Keyboard (Text Entry)

When text entry is necessary, an on-screen keyboard appears.

Names and Functions of Parts (On-Screen Keyboard)

You can display an alphanumeric, symbol, or umlaut character set on the on-screen keyboard. The name and function of each part on the keyboard is as follows.

Example: On-screen keyboard during alphanumeric character display

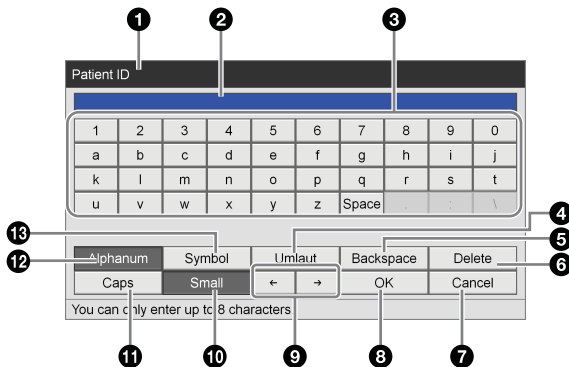

- 1 Title bar**
The name of the setting item for which the on-screen keyboard was opened is displayed here.
- 2 Entry box**
Characters selected on the keyboard are entered here.

Tip

Characters input on a hardware keyboard are also entered here.

- 3 Keyboard**
When you select a character you want to input using the **↑**, **↓**, **←**, and **→** buttons and then press the ENTER button, the character is entered in the **2** entry box.
- 4 [Umlaut]**
Select this and press the ENTER button to switch the keyboard to umlaut character display.

- 5 [Backspace]**
Select this and press the ENTER button to discard the character immediately preceding the cursor.
- 6 [Delete]**
Select this and press the ENTER button to discard the character immediately following the cursor.
- 7 [Cancel]**
Select this and press the ENTER button to cancel text entry and close the on-screen keyboard.
- 8 [OK]**
When you select this and press the ENTER button, the text that appears in the entry box is applied and input in the settings screen.
- 9 **←/→****
Select this and press the ENTER button to move the cursor one space to the left or right.
- 10 [Small]**
Select this and press the ENTER button to switch the keyboard to lower case display.
- 11 [Caps]**
Select this and press the ENTER button to switch the keyboard to upper case display.
- 12 [Alphanumeric]**
Select this and press the ENTER button to switch the keyboard to alphanumeric character display.
- 13 [Symbol]**
Select this and press the ENTER button to switch the keyboard to symbol display.

Handling Discs

Compatible discs

This unit supports the following BD and DVD discs.

- DVD-R
- BD-R / BD-R DL
- BD-RE / BD-RE DL

Tip

Use Sony discs.


Disc formatting

This unit automatically formats new and previously unused discs.

Notes on Handling

Handling

- Do not touch the recording or playback surface of the disc. Handle the disc by its edge.

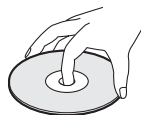

- Do not use the following types of disc. Doing so may cause the unit to malfunction.
 - A disc on which a paper label or a seal is attached
 - A disc on which the glue of cellophane tape or a label, or on which a portion of a label still remains
 - A cracked disc
 - A cracked disc which has been mended using glue (or a similar substance)




Storage

- Do not store discs where they may be subjected to direct sunlight, or in other places where the temperature or humidity is high.
- Store discs in their cases.
- Finger prints or dust accumulated on the disc cause deterioration of picture quality. Keep the disc clean.

Cleaning

- Clean the disc with a soft cloth, wiping it from the center out. If the dust is heavy, wipe it with a soft cloth moistened with water, then wipe off the water with a dry cloth.

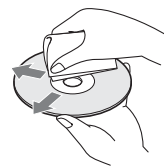

- Do not use solvents such as benzene, thinner, commercially available record cleaners, or anti-static spray. Using these products may damage discs.


Inserting and Removing Discs

To insert a disc

- 1 Turn the unit on.


For details, see “To turn on the unit” (page 32).

- 2 Press the  EJECT button.

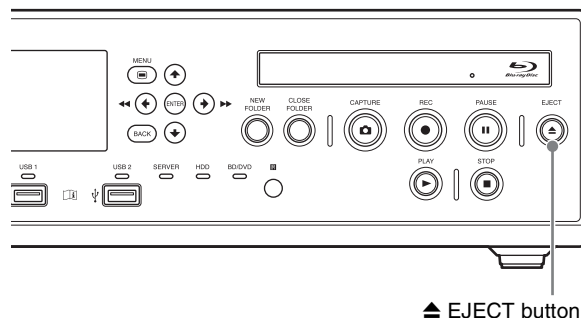

The disc tray opens.

- 3 Insert the disc.

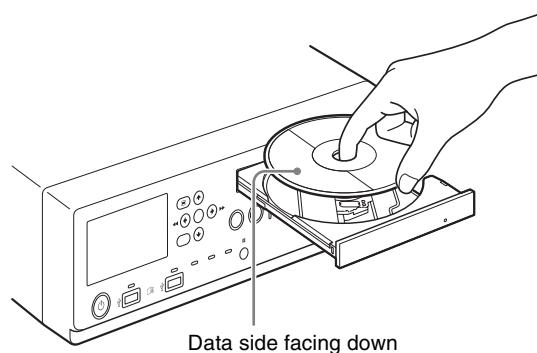

Tip

Press the center of the disc until it clicks into place to ensure that it is placed securely. If the disc is not placed securely, the disc may fall out or be damaged or other malfunctions may occur.

- 4 Gently push in the disc tray.

The disc tray closes.

When an incompatible disc is inserted

The BD/DVD indicator lights orange. A message notifying you of the incompatible media also appears on the menu screen. Remove the disc, and insert a compatible disc.

To remove a disc

- 1 Press the  EJECT button while the unit is turned on.

The disc tray opens.

- 2 Remove the disc.

- 3 Gently push in the disc tray to close it.

Tip

When simultaneously saving to a BD/DVD, the disc tray cannot be opened while recording is in progress or paused.


Connections

- 1 Connect a power cord (not supplied) to the ~ AC IN connector on the rear of the unit.
- 2 Connect the output signal of the medical equipment to SDI-IN connectors A to D on the rear of the unit.
- 3 Connect the plug on the power cord to an AC power supply.
- 4 Turn on the medical equipment that will act as the signal source.

Turning the Unit On and Off

To turn on the unit

- 1 Set the main power switch on the rear panel to the I (on) position.
- 2 Press the ⏻ (on/standby) switch on the front panel of the unit.

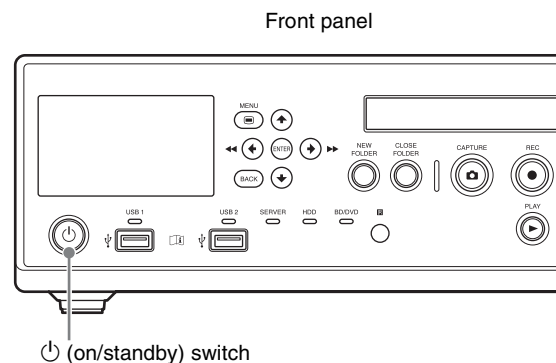

A splash screen appears on the front panel display when the unit turns on.

After “BOOTING...” and then the patient ID, current date, and other information appear on the front panel display, operations can be performed.

To turn on power or enter standby during normal operation

Use the ⏻ (on/standby) switch on the front panel to switch between turning the unit on and entering the standby mode.

“FINISHING” appears on the front panel display while the unit is shutting down. “GOOD BYE” appears and disappears after shutdown is complete, and the unit enters standby mode.

During long periods of disuse

Turn off the main power switch on the rear panel of the unit when not using the unit for long periods of time.

- 1 Press the  (on/standby) switch on the front panel of the unit to enter standby mode.

The  (on/standby) indicator lights orange.

- 2 Set the main power switch on the rear panel of the unit to the  (off) position.

The power is turned off.

Tip

To completely shut off the power supply, set the main power switch to the  (off) position, and then disconnect the power cord.

Configuring System Settings

Be sure to configure the system settings before you begin operating the unit. The system settings should be configured by the system administrator.

For details on system settings, see “System Administrator Settings” (page 61).

Tip


If the system administrator settings are not configured, operations may not be performed as intended in some cases.

Configuring Recording Settings (User Settings)


Configure input signal and image quality settings in the [User Settings] screen.
You can navigate menus using the navigation buttons on the front panel of the unit.

Displaying the [User Settings] Screen

- 1 Press the MENU button.
The [MENU] screen appears.
- 2 Select [Settings], and press the ENTER button.

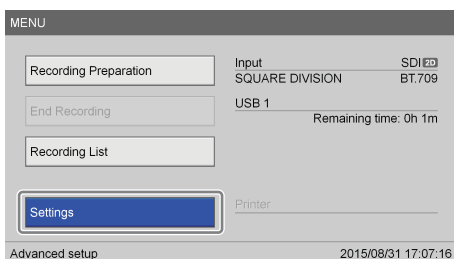

- The [Settings] screen appears.
- 3 Select [User Settings], and press the ENTER button.
The [User Settings] screen appears
 - 4 Select a tab, and configure the settings as necessary.

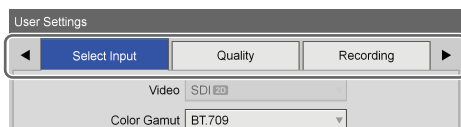

Configuring Input Signal Settings

Select the video and audio input signals.

- 1 Display the [Select Input] tab, and configure each setting.

[Video]

Select the video input signal.
You can select from [SDI 2D] (2D), [SDI 3D] (3D line-by-line), and [3D 3D] (3D top-and-bottom).
This item is available when [2D/3D] is set to [Manual] in the [System Admin Settings] screen –

[Function Settings] screen – [Input Detection] tab.
When [Auto] is selected, the automatically detected input signal will be displayed.
The default setting is [SDI 2D] (2D).

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this is fixed at [SDI 2D].

[Color Gamut]

Select the color gamut for the input signal.
You can select from [BT.709] and [BT.2020].
This item is available when [Color Gamut] is set to [Manual] in the [System Admin Settings] screen – [Function Settings] screen – [Input Detection] tab.
When [Auto] is selected, the automatically detected color gamut will be displayed.
The default setting is [BT.2020].

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this setting is disabled.

[Audio]

Select the audio input format.
You can select from [LINE] and [MIC].
The default setting is [LINE].

[Mic Plug-In Power]

Select whether to supply power to the microphone from the unit.
This item is available when [Audio] is set to [MIC].
The default setting is [Not Used].
You can configure the voltage that is supplied in the [System Admin Settings].

- 2 To configure image quality settings, proceed to “Configuring Image Quality Settings” (page 34).

If you are finished configuring settings, select [Apply] and press the ENTER button.

Configuring Image Quality Settings

Configure image quality settings for videos and still images.

- 1 Display the [Quality] tab, and configure video settings in the [Video] area.

Tip

The recording capacities displayed in [Rec Time] and [Rec Quantity] are estimations for single-layer BD discs (25 GB) based on the specified [Ch1/Ch2 Rec Mode] setting in the [Recording] tab.

[Quality]

Select the video recording quality.

You can select from [Best], [High], and [Standard].
The default setting is [High].

Tip

The setting specified here will become the default setting for [Video Quality] in the [Recording Preparation] screen.

2 Configure still image settings in the [Still] area.

[Format]

Select the still image format (file format).

You can select from [Bitmap], [TIFF], and [JPEG].
The default setting is [Bitmap].

[Quality]

Select the still image (JPEG) recording quality.

You can select from [Best], [High], and [Standard].
This item is available when [Format] is set to [JPEG].
The default setting is [High].

[Capture]

Displays the still image format used during image capture.

This is fixed at [Frame] on the unit.

3 To configure recording settings, proceed to “Configuring Recording Settings (User Settings)” (page 34).

If you are finished configuring settings, select [Apply] and press the ENTER button.

Configuring Recording Settings

1 Display the [Recording] tab, and configure recording settings.

[Ch1/Ch2 Rec Mode]

Select the recording mode.

You can select from [Ch1 Rec], [Ch1/Ch2 Simul Rec], and [Ch1 2K Rec].

When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, only [Ch1/Ch2 Simul Rec] and [Ch1 2K Rec] can be selected.

The default setting is [Ch1/Ch2 Simul Rec].

For details on Ch1/Ch2 simultaneous recording, see “Using the Ch1/Ch2 Simultaneous Recording Function” (page 52).

For details on Ch1 2K recording, see “Using the Ch1 2K Recording Function” (page 52).

Tips

- When [Ch1 2K Rec] is selected, the following [Ch2 Rec Format], [Convert Color Gamut], and [4K Convert Mode] settings will be grayed out.
- When [CMS Mode] is changed from [Not Used] to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, the setting will be configured as follows.
 - If [Ch1 Rec] or [Ch1/Ch2 Simul Rec] was selected before the change, the setting will be configured as [Ch1/Ch2 Simul Rec].
 - If [Ch1 2K Rec] was selected before the change, the setting will be configured as [Ch1 2K Rec].

[Ch2 Rec Format]

Select 3D or 2D for the Ch2 recording format.

Tips

- When [Ch1/Ch2 Rec Mode] is set to [Ch1/Ch2 Simul Rec], the [Link On] (Ch1/Ch2 link) setting for the touch panel monitor will be turned on (green).
- 3D recording will be performed according to the input signal format. However, images will be recorded in top-and-bottom format when the input is a line-by-line format.
- When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab and [Ch1/Ch2 Rec Mode] is set to [Ch1/Ch2 Simul Rec], the [Ch2 Rec Format] setting will be fixed at the current selection.

[Convert Color Gamut]

Select the conversion format for the recorded color depth.

You can select from [Do Not Convert], [Ch1/Ch2 to BT.709] and [Ch2 to BT.709].

The default setting is [Ch2 to BT.709].

[4K Convert Mode]

Select the method used to convert to 3840×2160 resolution when 4096×2160 resolution signals are input.

You can select from [Cut] and [Downconvert].

The default setting is [Cut].

For details on the actual images that are converted when [Cut] or [Downconvert] is configured, see “4K Convert Modes” (page 40).

Note

When using 3D input signals, [Cut] mode will be used regardless of the [4K Convert Mode] setting.

- 2 To configure save settings, proceed to “*Configuring Save Settings*” (page 36).

If you are finished configuring settings, select [Apply] and press the ENTER button.

Configuring Save Settings

Configure settings related to storage destinations for recorded data and simultaneous transfer, and specify whether to output metadata.

- 1 Display the [External Storage 1] tab, and configure settings related to storage destinations for recorded data.

[Ch1 External Storage]

Select the external storage media that will be used when performing simultaneous recording.

You can select from [Not Specified], [BD/DVD], [USB 1], [USB 2], [SERVER], and [Auto-Detect]. The default setting is [Not Specified].

- When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, the setting will be fixed at the current selection.
- When [Prohibit USB/BD/DVD] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the options for selection will be [Not Specified] and [SERVER].
- When an external storage destination other than a server is configured and [Prohibit USB/BD/DVD] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the setting will return to the default setting.

[Ch2 External Storage]

The external storage media selected in [Ch1 External Storage] appears here.

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this item will be grayed out.

[Record Save to Media]

Select the external storage media that will be used when copying recording lists and image lists from the [Recording List] and [Image List] screens.

You can select from [BD/DVD], [USB 1], [USB 2], and [SERVER].

The default setting is [USB 1].

- When [Prohibit USB/BD/DVD] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the setting will be fixed at [SERVER].
- When an external storage media other than a server is configured and [Prohibit USB/BD/DVD] is set to

[Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the setting will be fixed at [SERVER].

- When [Prohibit USB/BD/DVD] is changed from [Use] to [Not Used] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the setting will be set to [USB 1].

[Folder Structure]

Select the directory structure that will be used when writing to external destinations.

You can select from [Date + Patient ID] and [Patient ID].

The default setting is [Date + Patient ID].

When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, the setting will be fixed at the current selection.

Tip

When recorded data with the same patient ID is created, a new recording date folder is created under the patient ID folder, and the recorded data is saved to the new folder.

- 2 Display the [External Storage 2] tab, and configure settings related to simultaneous transfer.

[Save Video]

Select the check boxes of the video formats that will be simultaneously transferred. (Multiple selections can be made.)

You can select from [Ch1 2D], [Ch1 3D], [Ch2 2D], and [Ch2 3D].

All of the check boxes are selected under default settings.

[Save Still]

Select the check boxes of the storage formats for 3D still images. (Multiple selections can be made.)

You can select from [3D], [2D], [Convert to 2D DICOM], [2D Left], [Convert to 2D DICOM Left], [2D Right], and [Convert to 2D DICOM Right].

All of the check boxes are selected under default settings.

Tips

- When [Save Still] is set to [3D], the data will be saved in MPF format.
- When none of the check boxes are selected, [Save Video] and [Save Still] will be blank. When multiple formats are selected, they will be displayed with slashes (/) used as separators. If the entire text cannot be displayed, “...” will be displayed.

- 3** Display the [Metadata] tab, and configure metadata output settings.

[Metadata Output]

Select the check boxes of the metadata items you want to transfer to the external media. (Multiple selections can be made.)

You can select from [Recording Date], [Patient Info.], [Doctor], [Case], and [Facility Name].

None of the check boxes are selected under default settings.

Tips

- Metadata is saved to the uppermost directory of the recorded data folder.
- When [Patient Info.] is selected, the patient's ID, name, gender, and date of birth will be output as metadata.

[Metadata Format]

Displays the format of the metadata outputs.

[Encryption]

Select whether to encrypt the metadata that is output. The default setting is [Not Used].

[Key]

This item is available when [Encryption] is set to [Use].

Enter the encryption key using 32 hexadecimal digits.

Tip

When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, the setting will be fixed at the current selection and cannot be changed.

- 4** To configure print settings, proceed to “*Configuring Print Settings*” (page 37).

If you are finished configuring settings, select [Apply] and press the ENTER button.

Configuring Print Settings

Configure print settings for still images.

The unit supports connection and use of the Sony UP-DR80MD printer.

Tip

You can select the printer for use from the [System Admin Settings] screen – [Device Settings] screen – [Device 1] tab.

For details on this setting, see “[Device 1] Tab” (page 67).

- 1** Display the [Print 1] tab, and configure settings related to print methods and items to be printed.

The name of the currently connected printer appears in the [Printer] box.

[Layout]

Select the number of still images to print on each sheet.

For the portrait orientation

You can select from 1 image (1 × 1), 2 images (1 × 2), 4 images (2 × 2), 6 images (2 × 3), 8 images (2 × 4), 9 images (3 × 3), 12 images (3 × 4), 15 images (3 × 5), or 18 images (3 × 6).

The default setting is [6 images (2 × 3)].

For the landscape orientation

You can select from 1 image (1 × 1), 2 images (2 × 1), 4 images (2 × 2), 6 images (3 × 2), 8 images (4 × 2), 9 images (3 × 3), 12 images (4 × 3), 15 images (5 × 3), or 18 images (6 × 3).

The default setting is [6 images (3 × 2)].

[Auto Print]

Select whether to print automatically at the moment a still image is captured.

The default setting is [Not Used].

[Add. Info.]

Select the check boxes of the additional information items you want to print. (Multiple selections can be made.)

You can select from [Recording Date], [Patient Info.], [Doctor], [Case], [Recording Timestamp], [Logo], [Facility Name], and [Watermark].

None of the check boxes are selected under default settings.

[Watermark]

Select whether to print watermarks.

Tip

To print watermarks or logos, you must import watermark or logo files that were created in the proper format beforehand. You can import watermark or logo files from the [System Admin Settings] screen – [Device Settings] screen – [Device 1] tab.

For details on importing, see “[Import Watermark]” (page 68) and “[Import Logo]” (page 68).

- 2** Display the [Print 2] tab, and configure settings related to paper size and print quantity.

[Paper]

Select the printer's paper size.
The default setting is [Letter].

[Orientation]

Select the orientation for printing.
The default setting is [Portrait].

[Copies]

Specify the number of copies to print (1 to 10).
The default setting is [1].

- 3 To configure screen display settings, proceed to "Configuring Other Settings" (page 38).

If you are finished configuring settings, select [Apply] and press the ENTER button.

Configuring Other Settings

Configure settings related to the screen displays.

- 1 Display the [Other] tab, and configure the following settings.

[Menu Display]

Select 2D or 3D for the menu screen display format.
The default setting is [2D].

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this item will be grayed out.

[Select Right/Left of Image]

Select which image (left or right) to use for 2D recorded data when recording still images or recording simultaneously.
The default setting is [Left].

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this item will be grayed out.

- 2 When you finish configuration, select [Apply] and press the ENTER button.

The [Settings] screen appears again.

Function Keys

When using the unit, you can assign certain functions (recording, playback, etc.) to the function keys (F1 to F12) at the top of a keyboard and perform these functions with a single keystroke.

We recommend assigning frequently used functions based on your working environment.

Tip


To perform controls for only Ch2 without the use of a touch panel monitor, the Ch2 controls must be assigned to function keys beforehand.

Using Function Keys

Tip


The function keys are enabled only when the hard keyboard mode is enabled.

- 1 Press the MENU button.
The [MENU] screen appears.
- 2 Select [Settings], and press the ENTER button.
The [Settings] screen appears.
- 3 Select [System Admin Settings], and press the ENTER button.
The [System Admin Settings] screen appears.
- 4 Select [Device Settings], and press the ENTER button.
The [Device Settings] screen appears.
- 5 Display the [Device 2] tab, and set [Keyboard Mode] to [Hard Keyboard].
- 6 Select [Advanced] for [Keyboard Mode], and press the ENTER button.
The [Function Key Setting] screen appears.
- 7 Select the check boxes of the necessary function keys, select [Edit], and then press the ENTER button.

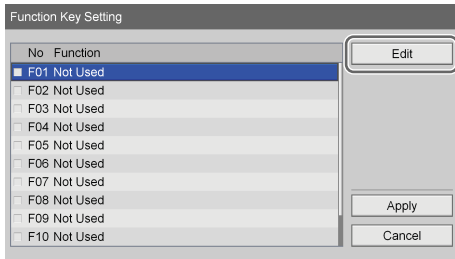

The [Edit] screen appears.

- 8 Select the functions to assign to the function keys under [Function], select [Apply], and then press the ENTER button.

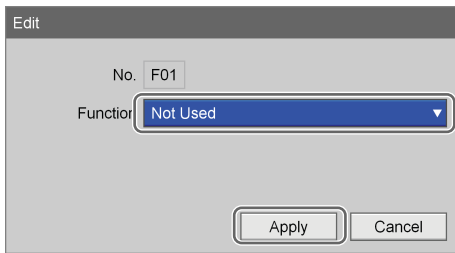

The [Function Key Setting] screen appears again.

- 9 Select [Apply], and press the ENTER button.

The [Function Key Setting] screen closes, and the [Device 2] tab appears again.

- [Record]: Record video and audio.
- [Set Chapter]: Insert a chapter separator without saving a still image of the recording.
- [Close Folder]: End the recording session.

Other Assignable Functions

Aside from the factory assigned functions, the following functions can also be assigned to the function keys.

- [Audio Rec On]: Record audio.
- [Audio Rec Off]: Disable audio recording.
- [Menu Display]: Display the [MENU] screen.
- [Chapter Next]: Jump forward to the next chapter.
- [Chapter Back]: Jump back to the previous chapter.
- [1 GOP Next]: Jump forward to the next GOP.
- [1 GOP Back]: Jump back to the previous GOP.
- [Eject]: Eject the inserted disc.
- [Rec List Display]: Display the [Recording List] screen.
- [Status Display On]: Display the unit's operation status (i.e., recording, play, pause).
- [Status Display Off]: Hide the unit's operation status (i.e., recording, play, pause).
- [Test Record]: Perform a test recording.
- [Ch2 Rec]: Start Ch2 recording.
- [Ch2 Stop]: Stop Ch2 recording or playback.
- [Ch2 Play]: Play back Ch2 recorded data.
- [Ch2 Pause]: Pause Ch2 recording or playback.

Factory Assigned Functions

The following functions are assigned to the function keys under factory default conditions.

Function Keys	Function
F1	Rec Prep / New Folder
F2	Playback
F3	Fast forward
F4	Rewind
F5	Stop
F6	Pause
F7	Record
F8	Capture
F9	Call list
F10	Add to list
F11	Set Chapter
F12	Close Folder

Tip

The following functions operate as follows.

- [Rec Prep / New Folder]: Display the [Recording Preparation] screen. A new folder is created.


4K Convert Modes

The following image conversions are performed when [Cut] and [Downconvert] are selected for [4K Convert Mode]. You can configure the [4K Convert Mode] setting in the [User Settings] screen - [Recording] tab. For details on this setting, see “*Configuring Recording Settings*” (page 35).

[Cut] mode


4K → 4K conversion

- The input image is recorded with 128 pixels cut on both the left and right sides (shaded area in illustration).
- The image is output with black bars (128 pixel width) added to the left and right.

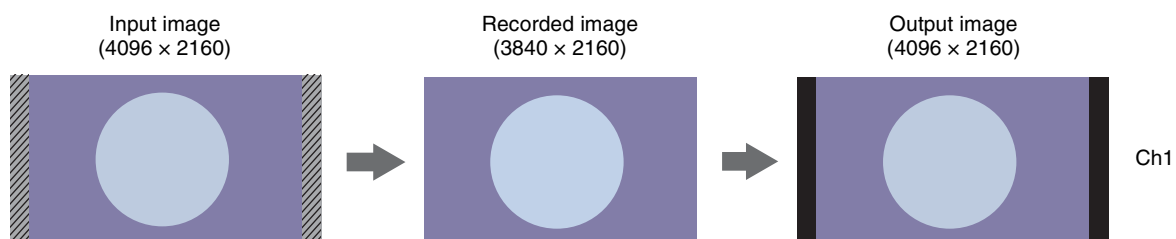

4K → 2K conversion


- The input image is recorded with 128 pixels cut on each of the left and right sides with its height and width reduced to 1/2 size.
- The recorded image is output as it is recorded.

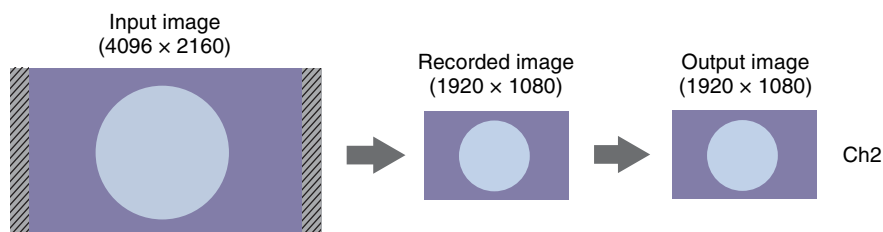

[Downconvert] mode


4K → 4K conversion

- The input image is recorded in reduced size with its original aspect ratio and black bars added to the top and bottom (67 and 68 pixel height respectively).
- For SDI outputs, the black bars are removed and the image is expanded to its original resolution.
- For HDMI outputs, black bars (128 pixel width) are added to the left and right of the recorded image.

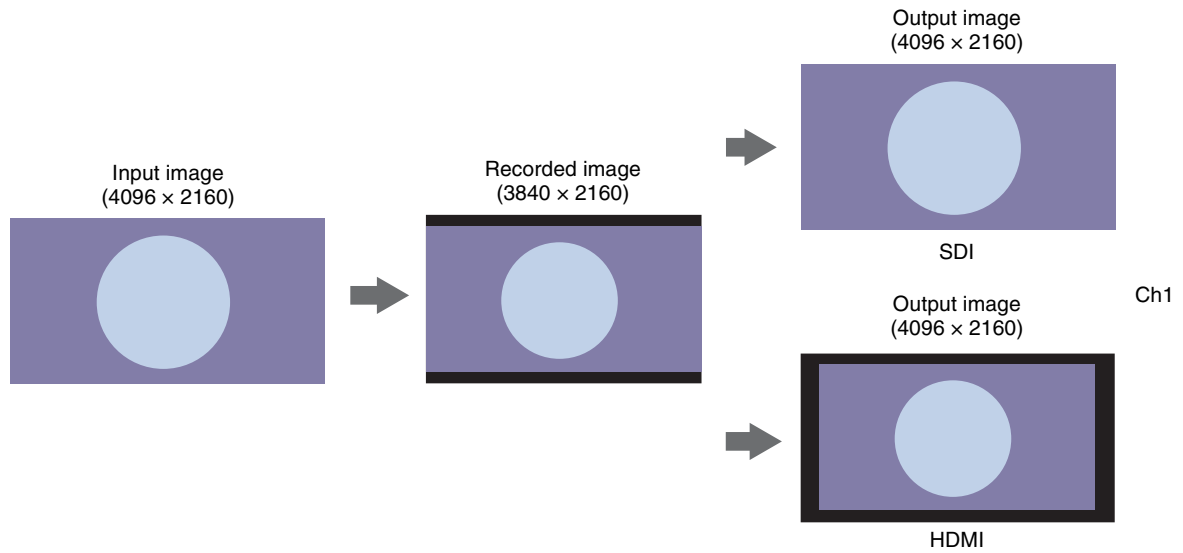

4K → 2K conversion

- The input image is recorded with its height and width reduced to 1/2.1 size and black bars added to the top and bottom (33.5 and 34 pixel height respectively).
- The recorded image is output as it is recorded.

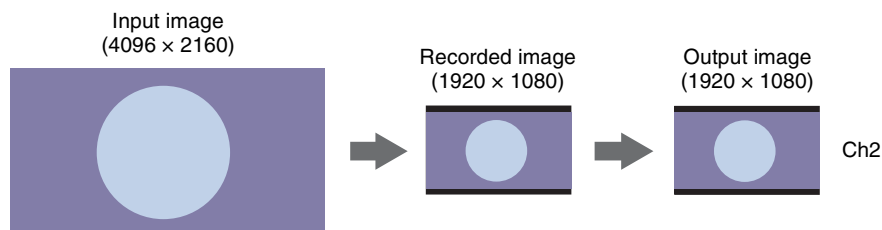

Note

When using 3D input signals, [Cut] mode will be used regardless of the [4K Convert Mode] setting.


Basic Recorder Operations

Chapter

3

Operation Flow

This chapter describes how to record and play back videos and capture still images while viewing the display on the front panel of the unit and using only the buttons on the front panel.

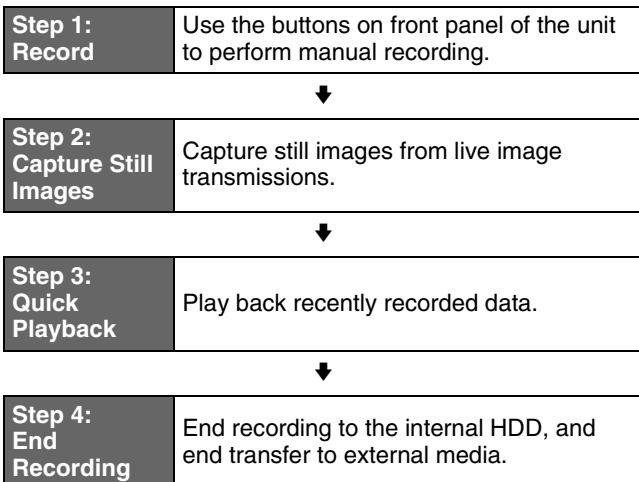

For details on recording preparation, recording, and playback operations, see “Chapter 4 Recording and Playback” (page 46).

Overview of recorded data

“Recorded data” on this unit refers to the data recorded from the time a “new folder” is created (i.e., when recording is started) to the time “close folder” is performed. Recorded data entries include titles. A single recorded data entry can include multiple titles. A title is created whenever a stop operation is performed during recording.

Titles can consist of video data or captured still images. Up to 15,000 recorded data entries can be registered on the unit.


A single recorded data entry can include up to 255 titles.

Tip


If an input signal is interrupted, a chapter separator will be inserted when it is restored. When saving to external media, files will be separated by chapters.

Step 1: Record

- 1 Press the ● REC button.

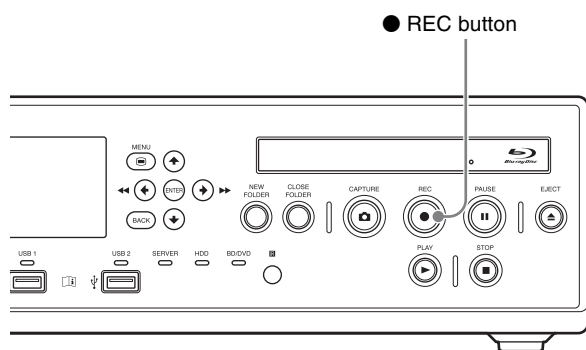


Recording starts.
“REC” appears while recording is in progress.

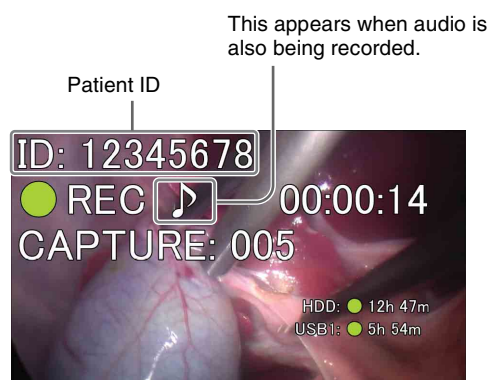

Tip

If images are recorded and patient information has not been entered in the [Recording Preparation] screen, patient IDs will be automatically assigned. The letter “U” appears at the beginning of automatically assigned patient IDs.

- 2 Press the ■ STOP button to stop recording.

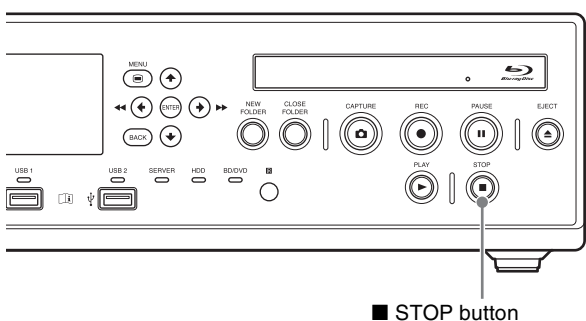

Recording stops, and “STOP” appears.

If the surgical procedure or examination is finished and you want to end the recording session, proceed to “Step 4: End Recording” (page 45) and perform the recording end procedure.


Tip

If you press the ● REC button again, recording will start again as a new recorded data entry.


Step 2: Capture Still Images

Capture still images from live image transmissions. Up to 500 still images can be captured for a single recorded data entry.

Press the  CAPTURE button.

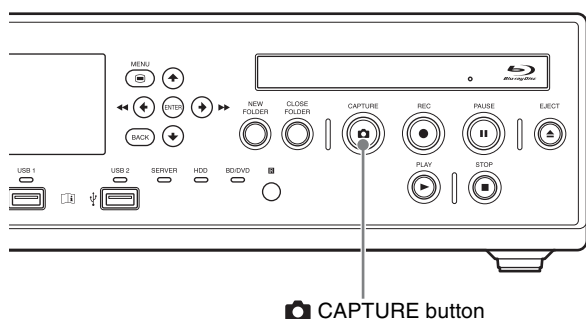

When you capture an image, the capture number appears.

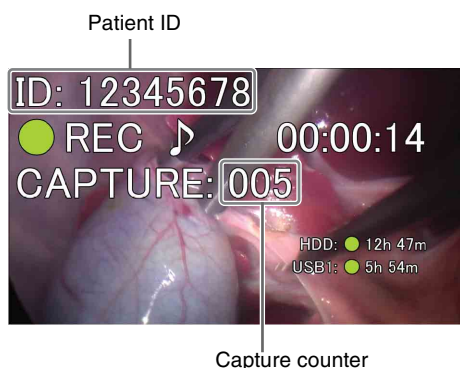

The still image data is saved.

For details on capturing still images, see “Capture Still Images” (page 51).


Recorded data entries are created when still images are captured.

If the surgical procedure or examination is finished and you want to end the recording session, proceed to “Step 4: End Recording” (page 45) and perform the recording end operation.


Step 3: Quick Playback

Play back data recorded in Step 1.

Press the  PLAY button.

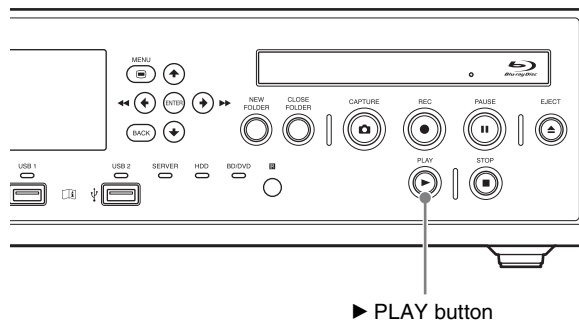

The most recent recorded data plays back. “PLAY” and the number of the chapter being played back appear while playback is in progress.

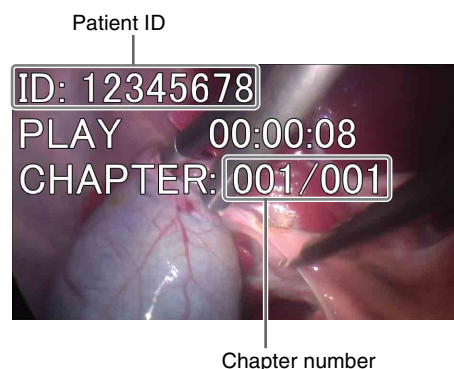

You can also perform the following operations during playback.

Fast forward playback

Press the  FF button.

Rewind playback

Press the  REW button.

Pause playback

Press the  PAUSE button.

To resume playback, press the  PAUSE button again or press the  PLAY button.


Stop playback

Press the  STOP button.


Step 4: End Recording

End the recording session to stop recording to the internal HDD and stop transfer to external media. This operation is referred to as “closing a folder.”

- 1 Press the ■ STOP button to stop recording.
- 2 Press the CLOSE FOLDER button.

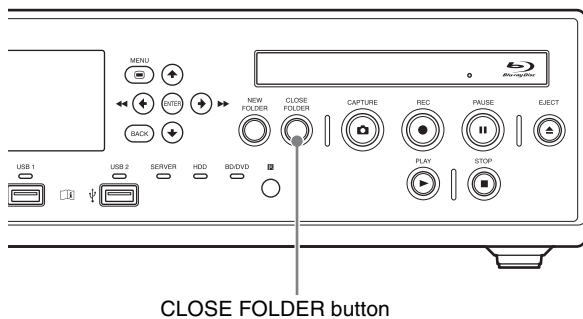

Select [OK] when the confirmation message appears, and press the ENTER button to display the [End Recording] screen.

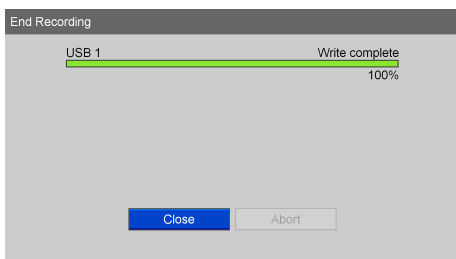

- 3 When the writing process is complete, select [Close] and press the ENTER button.
The [MENU] screen appears again.

Recording Preparation

Before starting to record, enter information regarding the patient of the surgical procedure or examination, configure the storage destination and quality settings for recorded data, and then perform a test recording. You can also register information for multiple patients in advance and recall the appropriate patient information on the day of the procedure or examination, or obtain patient information from MWL. You can prepare for recording and test record in the [Recording Preparation] screen.


Patient Information and Data Storage Settings

Enter information regarding the patient of the surgical procedure or examination, and configure the storage destination and quality settings for recorded data. Use the on-screen keyboard to enter text.

For details on use, see “Using the On-Screen Keyboard (Text Entry)” (page 29).

1 Press the NEW FOLDER button.

Alternatively, select [Recording Preparation] in the [MENU] screen, and press the ENTER button.

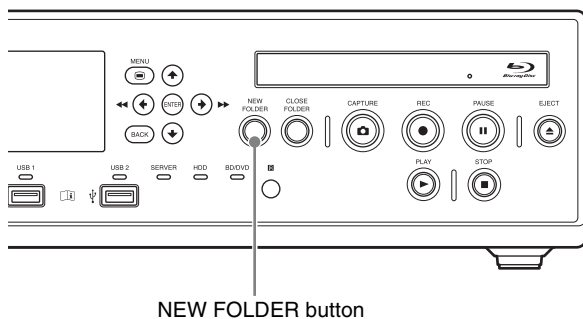

The [Recording Preparation] screen appears.

2 Enter patient information.

Patient information entry can be skipped.

[Patient ID]

Enter an identification number for the patient.

Tip

Certain symbols (¥ / : ? * " < > | . \) cannot be used for the patient ID.

[Name]

Enter up to 64 characters in each field of the patient's name.

For details on the input fields for the patient name, see “[Patient Name Field]” (page 64).

[Gender]

Select the gender of the patient.

You can select from [Male], [Female], or [Unknown], or leave the field empty.

[Birth Date]

Select the year, month, and day.

3 Configure storage settings for recorded data.

[Doctor]

Select the name of the doctor performing the surgical procedure or examination. This item can be skipped.

Tip

You can register doctors in the [System Admin Settings] screen – [Edit Doctor List] screen.

For details on this setting, see “Doctor List Registration” (page 73).

[Case]

Select the case.

Tip

The cases must be registered beforehand.

For details on this setting, see “Registering Cases” (page 76).

[Video Quality]

Select the video image quality.

You can select from [Best], [High], and [Standard]. This setting applies to both Ch1 and Ch2.

Tip

You can configure a default setting for this item in the [User Settings] screen – [Quality] tab.

For details on this setting, see “Configuring Image Quality Settings” (page 34).

[External Media]

When you want to simultaneously save to external media in addition to saving to the unit’s internal hard disk, configure these settings to specify the storage destination for recorded data.

If you do not want to simultaneously record to external media, select [Not Used].

If you want the unit to detect the inserted media automatically, select [Auto].

Tips

- If the external media is currently inserted and recognized, the remaining recording times are also displayed.
- Servers will not be automatically recognized.
- When auto recognition is enabled and multiple media are inserted, data will be recorded to the media in the following order of priority: USB1 → BD/DVD → USB2.
- You cannot change media types in the middle of recording when auto recognition is enabled.
- Recorded data for Ch2 will be recorded to the same external media specified for Ch1.
- When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab, the setting will be fixed at its current setting.
- When [Prohibit USB/BD/DVD] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the only options for the external storage destination setting will be [Not Used] and [SERVER].


- When a storage destination other than a server is configured and [Prohibit USB/BD/DVD] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [General 2] tab, the external storage destination will return to the default setting.

- 4 To perform a test recording, proceed to “Test Recording” (page 49).

If you are finished configuring settings, select [Done] and press the ENTER button.

Using a Hard Keyboard for Patient Information and Data Storage Settings

When [Keyboard Mode] is set to [Hard Keyboard] in the [System Admin Settings] screen – [Device Settings] screen, you can configure patient information and data storage settings using only the connected keyboard. When hard keyboard mode is enabled, a guide appears on the right side of the screen and the following keys function as follows.

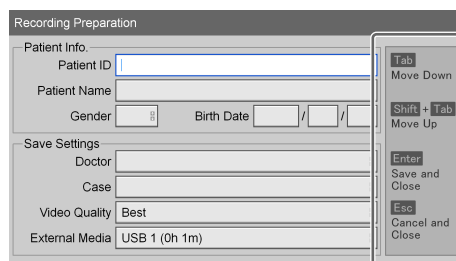

Key	Function
Tab	Move to the next input field.
Shift + Tab	Move to the previous input field.
↑ ↓ (arrows)	Select items.
← → (arrows)	Move within entered text.
Enter	Apply the setting, and close the [Recording Preparation] screen.
Esc	Cancel the setting, and close the [Recording Preparation] screen.

Registering Multiple Patients in Advance

You can register information for multiple patients of procedures and examinations in advance. By doing so, you can simply select the patient from the registered patient list on the day of the procedure or examination and save time on preparation.

For details on registering patients from MWL, see “Registering Patients from MWL” (page 48).

To register patients

Register patients for procedures and examinations to the patient list.

- 1 Configure patient information and storage settings for recorded data in the [Recording Preparation] screen, select [Add to List], and then press the ENTER button.

For details on patient information and save settings for recorded data, see “Patient Information and Data Storage Settings” (page 46).

The patient is added to the patient list.

- 2 Repeat to add any additional patients as necessary.

Tip

You can register up to 100 patients to the patient list. When you exceed 100, existing patients will be overwritten automatically starting with the oldest.

To recall patients

Recall patients that have been registered in advance.

- 1 Select [Call List] in the [Recording Preparation] screen, and press the ENTER button.

Tip

[Call List] appears when [MWM] is set to [Not Used] in the [System Admin Settings] screen – [DICOM Settings] screen – [MWM Server 1] tab.

The [Patient List] screen appears.

- 2 Select the target patient, and press the ENTER button to select its check box.
- 3 Select [Execute], and press the ENTER button.

The [Recording Preparation] screen appears again, and the patient information will be applied in the [Recording Preparation] screen.

To edit registered patient settings

Edit patient information and save settings for recorded data that were registered in advance.

- 1 Select [Call List] in the [Recording Preparation] screen, and press the ENTER button.

The [Patient List] screen appears.

- 2 Select the patient for which you want to edit settings, and press the ENTER button to select its check box.

- 3 Select [Change], and press the ENTER button.

The [Change Patient List] screen appears.

- 4 When you finish editing settings, select [Apply] and press the ENTER button.

The [Patient List] screen appears again.

- 5 Select [Execute], and press the ENTER button.

The [Recording Preparation] screen appears again.

To delete registered patients from the list

Delete patients that were registered in advance from the list.

- 1 Select [Call List] in the [Recording Preparation] screen, and press the ENTER button.

The [Patient List] screen appears.

- 2 Select the patient to delete, and press the ENTER button to select its check box.

- 3 Select [Delete], and press the ENTER button.

- 4 Select [OK] when the confirmation message appears, and press the ENTER button.

The patient entry is deleted.

- 5 Select [Back], and press the ENTER button.

The [Recording Preparation] screen appears again.

Registering Patients from MWL

Patient information can be obtained from MWL and entered in the [Recording Preparation] screen.

Tips

- Patient information obtained from MWL cannot be edited.
- When patient information is obtained from MWL, the [Doctor] field will be grayed out.

Preparation

Set [MWM] to [Use] in the [System Admin Settings] screen – [DICOM Settings] screen – [MWM Server 1] tab.

For details, see “[MWM Server 1] Tab” (page 78).

To access MWL

- 1 Select [Access MWL] in the [Recording Preparation] screen, and press the ENTER button.

Tip

[Access MWL] appears when [MWM] is set to [Use] in the [System Admin Settings] screen – [DICOM Settings] screen – [MWM Server 1] tab.

The [MWL] screen appears.

- 2 Select the target patient, and press the ENTER button.

Tips

- Any invalid characters or symbols in the patient ID and patient name will be replaced and displayed as underscores (_).
- If the number of digits is less than the value configured for the [Patient ID Length] setting in the [System Admin Settings] screen – [Function Settings] screen – [Patient Info.] tab, zeros (0) will be added according to the [Add 0 to Patient ID] setting.
- You can also specify search conditions to search for specific patients.
- If patients are not registered to MWL, an empty list will be displayed.

The [Recording Preparation] screen appears again, and the patient information that was obtained from MWL will be applied.

To search for patients

Tip

Modalities are always used in searches. When search conditions are specified, the search will be performed using the specified conditions and modalities.

- 1 Select [Search] in the [MWL] screen, and press the ENTER button.

The [Search] screen appears.

- 2 Specify the search conditions, select [Start], and press the ENTER button.

[Appt. Date]

Select the year (last two digits), month, and day. You can also specify a period between dates.

[Patient ID]

Enter an identification number for the patient.

[Patient Name]

Enter the patient's name.

Last name only, first name only, and last name and first name (separated by a space) are supported for searches. (Middle names are not supported.)

[AE Title]

Enter the AE title.

[Accession No.]

Enter the accession number.

Tip

Select [Clear] and press the ENTER button to clear all specified search conditions.

A search for patient information matching the specified conditions is performed, and the results are displayed in the [MWL] screen.

Test Recording

Confirm whether images record properly.

When you perform a test recording, images record for 5 seconds and then play back automatically.

- 1 Select [Test Record] in the [Recording Preparation] screen, and press the ENTER button.
- 2 Select [OK] when the confirmation message appears, and press the ENTER button.

Recording starts, and a message indicating that test recording is in progress appears throughout the process.

After 5 seconds of recording, playback starts automatically.

When test recording is finished, the [Recording Preparation] screen appears again.

- 3 Select [Done], and press the ENTER button.

The [MENU] screen appears again.

Recording

Images can be recorded manually or via contact switch. This section describes how to record manually.

For details on configuring patient information and storage destinations for recorded data beforehand, see “Recording Preparation” (page 46).

Note

Signals will only be input from SDI input connector C during [Ch1 2K Rec] mode.

Tip

For details on how to record via contact switch, contact your local Sony representative.

Note on the auto delete function for recorded data

Whenever the remaining capacity on the unit’s internal hard disk becomes low, an auto delete function automatically deletes recorded data, starting with the oldest start dates for recording. This is done to maintain free space on the unit’s internal hard disk at all times. To disable the auto delete function, set [Auto Delete] to [Not Used] in the [System Admin Settings] screen.

Manual Recording

- 1 Press the ● REC button.

Recording starts.
“REC” appears on the front panel display while recording is in progress.

- 2 Press the ■ STOP button to stop recording.

- 3 To end the surgical procedure or examination recording session, press the CLOSE FOLDER button.

A confirmation message for ending recording appears.

- 4 Select [OK], and press the ENTER button.

The [End Recording] screen appears, and the status of data transfer to external media, for example, appears. Wait a moment for the writing process to complete.

- 5 When the writing process is complete, select [Close] and press the ENTER button.

To cancel transfer to external media, select [Abort] and press the ENTER button.

The [MENU] screen appears again.

Tip

You can record up to 24 continuous hours from the moment you start recording or capture a still image. After 24 hours, recording will stop automatically, and the “close folder” is performed.

Recorded data directory

Video data is stored in one of the following directories, depending on the system administrator settings.

Video data directory

<uppermost folder of the external media>/<recording start date (year, month, day, hour, minute, second)_patient ID>/MOVIE/
or

<uppermost folder of the external media>/<patient ID>/<recording start date (year, month, day, hour, minute, second)>/MOVIE/

Tips

- The order of the year, month, and day of the recording date used in the storage directory is based on the date format of the system settings.
- If the power supply is disconnected before the writing of the data is complete, the validity of the written data cannot be guaranteed.
- When recording to external media, the media can be switched out up to 15 times for a single recorded data entry. When this limit is exceeded, the “Exch. media limit exceeded” message appears on the [End Recording] screen.
- If you press the ● REC button again, recording will start again as a new title for the same patient.

Capture Still Images

You can capture a particular shot during a live transmission, and store it as a still image file. Still images can be captured manually or via contact switch. This section describes manual operations.

Tip

- The size of the still images that can be captured is the same as the resolution of the input signal.
 - When [Ch1 Rec] or [Ch1/Ch2 Simul Rec] is selected, still images will be recorded in 4096 × 2160 or 3840 × 2160 size.
 - When [Ch1 2K Rec] is selected, still images will be recorded in 1920 × 1080 size.
- For details on contact switch controls, contact your local Sony representative.

Manual Still Image Capture

To capture a still image, press the  CAPTURE button at the moment you want to capture the still image. A chapter separator is also created at the position in which the still image was captured. Recorded data entries are created when still images are captured. If the surgical procedure or examination is finished and you want to end the recording session, perform the recording end procedure.

For details on use, see “Step 4: End Recording” (page 45).

Tip

You can configure format setting for still images files in the [User Settings] screen – [Quality] tab.

For details on this setting, see “Configuring Image Quality Settings” (page 34).

Still image data directory

Still image data is stored in one of the following directories, depending on the system administrator settings.

Still image directory

<uppermost folder of the external media>/<recording start date (year, month, day, hour, minute, second)_patient ID>/STILL/
or

<uppermost folder of the external media>/<patient ID_recording start date (year, month, day, hour, minute, second)>/STILL/

Tip

The order of the year, month, and day of the recording date used in the storage directory is based on the date format of the system settings.

Using the Ch1/Ch2 Simultaneous Recording Function

This function allows automatic Ch2 recording (2K) during Ch1 recording (4K).

Tip

In [Ch1/Ch2 Simul Rec] mode, the 4K signals input to Ch2 are downconverted to 2K 1920×1080 format and recorded while Ch1 recording (4K) is in progress.

Configuring Settings for Ch1/Ch2 Simultaneous Recording

Configure the following settings in the [User Settings] screen - [Recording] tab.

[Ch1/Ch2 Rec Mode]

Select [Ch1/Ch2 Simul Rec].

[Ch2 Rec Format]

Select 2D or 3D.

For details on the [Recording] tab, see “Configuring Recording Settings” (page 35).

Tip

When using 2D input signals, 2D mode will be used regardless of the [Ch2 Rec Format] setting.

Performing Ch1/Ch2 Simultaneous Recording

Operations for Ch1/Ch2 simultaneous recording are identical to those for Ch1 recording.

When you press the ● REC button, recording for Ch1 starts and triggers recording for Ch2 to start.

Using the Ch1 2K Recording Function

This function records 1920 × 1080, 59.94p/50p signals input to SDI input connector C in 1920 × 1080 (1080p) resolution.

Still images will also be recorded in 1920 × 1080 size.

Notes

- Signals to connectors other than SDI input connector C will be ignored.
- Ch2 is disabled during [Ch1 2K Rec] mode.
- This function does not support color gamut configuration and the color gamut conversion function.

Configuring Settings for Ch1 2K Recording

Configure the following settings in the [User Settings] screen - [Recording] tab.

[Ch1/Ch2 Rec Mode]

Select [Ch1 2K Rec].

For details on the [Recording] tab, see “Configuring Recording Settings” (page 35).

Tip

The recording mode will be fixed at 2D.

Performing Ch1 2K Recording

Press the ● REC button for Ch1 to start recording.

Limitations during Playback and Output

- Videos and still images recorded in [Ch1 2K Rec] mode can only be played back in [Ch1 2K Rec] mode. In addition, videos and still images recorded in modes other than [Ch1 2K Rec] mode cannot be played back in [Ch1 2K Rec] mode.
- The display mode in the [MENU] screen will be fixed at 2D.

Playback

You can play back the most recent recorded data stored on the unit's internal hard disk via simple controls.

Notes

- Videos and still images recorded in [Ch1 2K Rec] mode can only be played back in [Ch1 2K Rec] mode. Similarly, videos and still images recorded in modes other than [Ch1 2K Rec] mode can be played back in modes other than [Ch1 2K Rec] mode.
- Signals will only be output to SDI output connector C during [Ch1 2K Rec] mode. The playback resolution will be 1920 × 1080 for both videos and still images.

You can also play back recorded data from search results. For details, see "Image Search" (page 54).

Tips

- If you press the MENU button during playback, playback will stop.
- When NTSC is selected as the video format under the system administrator settings, PAL-format recorded data cannot be played back properly. Conversely, when PAL is selected, NTSC-format recorded data cannot be played back properly.

Playing Back the Most Recent Data (Quick Playback)

Press the ► PLAY button to play back the most recently recorded data on the internal hard disk.

Tip

Ch1 recorded video data is played back with quick playback. Recorded data that does not include Ch1 video will not be played back.

Playback operations

You can use the buttons on the front panel of the unit to perform the following playback operations.

Fast forward playback

Press the ►► FF button.

Tip

Depending on the structure of the data, fast forward playback may not be possible for the entire recording.

Rewind playback

Press the ◀◀ REW button.

Pause playback

Press the ■■ PAUSE button.

To resume playback, press the ■■ PAUSE button again or press the ► PLAY button.

Stop playback

Press the ■ STOP button.

Image Search

You can search for recorded data stored on the unit's internal hard disk.

Use the on-screen keyboard to enter text.

For details on use, see "Using the On-Screen Keyboard (Text Entry)" (page 29).

Specifying Search Conditions

Search for recorded data by specifying conditions such as patient ID and recording date.

- 1 Press the MENU button.
The [MENU] screen appears.
- 2 Select [Recording List], and press the ENTER button.
The [Recording List] screen appears.
- 3 Select [Search], and press the ENTER button.
The [Search] screen appears.
- 4 Specify the search conditions, select [Search], and press the ENTER button.

[Patient ID]

Enter an identification number for the patient.

[Patient Name]

Enter the patient's name.

[Doctor]

Select the name of the doctor performing the surgical procedure or examination.

[Recording Date]

Select the year (last two digits), month, and day.

[Status]

Select whether the recorded data has been saved to external media.

You can select from [Saved to Media], [Not Saved to Media], or [DICOM Unsent], or leave the field empty.

Tip

The [DICOM Unsent] setting appears only if [DICOM] is set to [Use] in the [System Admin Settings] screen – [DICOM Settings] screen – [Server] tab.

[2D/3D]

Select whether the recorded data is 2D data or 3D data.

Tip

Select [Clear] and press the ENTER button to clear all specified search conditions.

- 5 Select [Search], and press the ENTER button.

A search for recorded data matching the specified conditions is performed, and the results are displayed.

Status of recorded data

The following icons appear based on the status of the recorded data.

Status	Description	Icon
Not saved to external media	The recorded data has not been saved to external media. The icon disappears when saving is completed.	
Protected	The recorded data is protected and cannot be deleted. The icon disappears when the protection is released.	
DICOM transfer failed	DICOM transfer failed. The icon disappears when transfer is successful.	
3D data	The data was recorded in 3D mode.	
2D data	The data was recorded in 2D mode.	
4K data	The recorded data is in 4K format.	
2K data	The recorded data is in 2K format.	

Viewing Thumbnails of Recorded Data

You can view thumbnail images of recorded data in an image list.

- 1 Select [Image List] for the operation box in the [Recording List] screen.
- 2 Select the recorded data in the [Recording List] screen, and press the ENTER button.

The [Image List] screen appears, and thumbnail images of the recorded data are displayed.



- 3 To change the type of images for which to display thumbnails, select the [Display] box and press the ENTER button.

The [Display] screen appears.

- 4 Select [Video], [Still], or [All], and press the ENTER button.

The display changes.

Tips

- If [Video] is selected, the thumbnails are displayed by title or by chapter, depending on the [Image List Display] setting in the [System Admin Settings] screen – [Function Settings] screen – [General 1] tab.
- You can also select recorded data for printing or copying to external media in the thumbnail display screen.

For details, see “Processing Recorded Data” (page 56).

Playback from the Search Results List

Notes

- Videos and still images recorded in [Ch1 2K Rec] mode can only be played back in [Ch1 2K Rec] mode. Similarly, videos and still images recorded in modes other than [Ch1 2K Rec] mode can be played back in modes other than [Ch1 2K Rec] mode.
- Signals will only be output to SDI output connector C during [Ch1 2K Rec] mode. The playback resolution will be 1920 × 1080 for both videos and still images.

- 1 Select [Image List] for the operation box in the [Recording List] screen.
- 2 To play back a title, select the recorded data you want to play back in the [Recording List] screen (search results list), and press the ► PLAY button.

- 3 To play back by chapter, select the recorded data you want to play back, and press the ENTER button.

The [Image List] screen appears, and thumbnail images of the recorded data are displayed.

Tip

To display recorded data by chapter, set [Image List Display] to [Chapter] in the [System Admin Settings] screen – [Function Settings] screen.

- 4 Select the recorded data you want to play back, and press the ► PLAY button.

The recorded data that is highlighted plays back.

Tip

Only the highlighted recorded data is played back.

Sorting the Recorded Data List

You can sort the recorded data list by recording date, patient ID, patient name, size, or save status.

Tip

The sizes displayed in the [Recording List] screen indicate total sizes for Ch1 and Ch2.

- 1 Select [Sort] in the [Recording List] screen, and press the ENTER button.
The [Sort] screen appears.
- 2 Select the [Item] box, and press the ENTER button.
The [Item] screen appears.
- 3 Select the category by which to sort, and press the ENTER button.
The [Sort] screen appears again.
- 4 Select the [Order] box, and press the ENTER button.
The [Order] screen appears.
- 5 Select [Descending] or [Ascending], and press the ENTER button.

Depending on the selected category, the descending and ascending orders are as follows.

Category	Descending	Ascending
[Recording Date]	Start with newest	Start with oldest
[Patient ID]	Alphabetical order → numerical order	Reverse numerical order → reverse alphabetical order
[Patient Name]	Alphabetical order → numerical order	Reverse numerical order → reverse alphabetical order
[Status]	[Not Saved to Media] → [Saved to Media]	[Saved to Media] → [Not Saved to Media]
[Size]	Start with largest	Start with smallest

The [Sort] screen appears again.

- 6** Select [Execute], and press the ENTER button.

The recorded data is sorted according to the specified order.

Processing Recorded Data

You can process recorded data stored on the unit's internal hard disk in the following ways.

- Printing (page 56)
- Copying to External Media (page 57)
- Editing Patient Information (page 58)
- Protecting Recorded Data (page 59)
- Deleting Recorded Data (page 59)
- Viewing Information on Recorded Data (page 59)

Printing

You can select still images for printing from the [Recording List] screen.

Tip

For details on configuring the number of images to print on each sheet and other print settings, see “Configuring Print Settings” (page 37).

- 1 Select [Image List] for the operation box in the [Recording List] screen.
- 2 Select the recorded data in the [Recording List] screen, and press the ENTER button.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

The [Image List] screen appears.

- 3 If thumbnail images of video data are displayed, switch to still image display.

For details on switching the display, see “Viewing Thumbnails of Recorded Data” (page 54).

- 4 Select the check boxes for the recorded data you want to print.



- 5 Select the [Operation] box, and press the ENTER button.

The [Operation] screen appears.

- 6 Select [Print], and press the ENTER button.

The [Image List] screen appears again.

- 7 Select [Execute], and press the ENTER button.

The [Print] screen appears.

- 8 Change the settings as necessary.

If necessary, you can enter a comment of up to 60 characters to be printed with the images.

For details on other settings, see “Configuring Print Settings” (page 37).

Printing starts.

The remaining amount of paper is displayed during printing.

Tip

To cancel printing, select [Abort] and press the ENTER button.

- 9 When printing is complete, press the ENTER button while [Close] is selected.

The [Image List] screen appears again.

Copying to External Media

You can copy recorded data stored on the unit’s internal hard disk onto a BD/DVD disc, USB memory device, or server.

Recorded data can be copied from the [Recording List] screen or [Image List] (thumbnail display) screen.

To select data for copying from the [Recording List] screen

This section describes how to select recorded data for copying onto external media in the [Recording List] screen. You can select whether to copy videos, still images, or both with this procedure.

Tips

- If the power supply is disconnected before the writing of the data is complete, the validity of the written data cannot be guaranteed.
- When recording to external media, the media can be switched out up to 15 times for a single recorded data

entry. When this limit is exceeded, the “Exch. media limit exceeded” message appears on the [Save to] screen.

- 1 If you are copying to a BD/DVD disc or USB memory device, insert the media into the unit.

- 2 Display the [Recording List] screen, select the operation box, and press the ENTER button.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

The [Operation] screen appears.

- 3 Select [Save to], and press the ENTER button.

The [Recording List] screen appears again.

- 4 Select the check boxes for the recorded data you want to copy.

- 5 Select [Execute], and press the ENTER button.

The [Save to] screen appears.

- 6 Configure save settings for the recorded data, select [Execute], and press the ENTER button.

[External Media]

Select the storage destination.

You can select from [BD/DVD], [USB 1], [USB 2], and [SERVER].

[File Server]

If [SERVER] is selected as the external storage destination, select a file server (the server configured for each doctor).

[Save Video]

This is enabled when recorded video data is selected. Select the format of the recorded data you want to save.

You can select from [Ch1 2D], [Ch1 3D], [Ch2 2D], and [Ch2 3D]. (Multiple selections can be made.)

All of the check boxes are selected under default settings.

[Save Still]


Select the check boxes of the storage formats for 3D still images.

You can select from [3D], [2D], [Convert to 2D DICOM], [2D Left], [Convert to 2D DICOM Left], [2D Right], and [Convert to 2D DICOM Right]. (Multiple selections can be made.)

When “convert to DICOM” is selected, the specified still image format is converted to DICOM before saving.

All of the check boxes are selected under default settings.

Copying of recorded data starts.

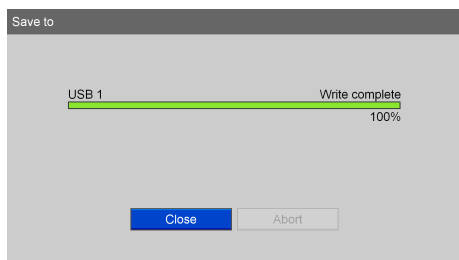

Tip

To cancel copying, select [Abort] and press the ENTER button.

- 7 When copying is complete, press the ENTER button while [Close] is selected.

The [Recording List] screen appears again.

To select data for copying from the [Image List] (thumbnail display) screen

This section describes how to select thumbnail images of recorded data in the [Image List] screen for copying onto external media. The copying procedure is identical for both videos and still images.

- 1 If you are copying to a BD/DVD disc or USB memory device, insert the media into the unit.
- 2 Select [Image List] for the operation box in the [Recording List] screen.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

- 3 Select the recorded data in the [Recording List] screen, and press the ENTER button.

The [Image List] screen appears.

- 4 Switch to the thumbnail display for videos when you want to copy video data, and switch to the thumbnail display for still images when you want to copy still image data.

For details on switching the display, see “Viewing Thumbnails of Recorded Data” (page 54).

- 5 Select the check boxes for the recorded data you want to copy.
- 6 Select the operation box, and press the ENTER button.

The [Operation] screen appears.

- 7 Select [Save to], and press the ENTER button.

The [Image List] screen appears again.

- 8 Select [Execute], and press the ENTER button.

The [Save to] screen appears.

- 9 Configure save settings for the recorded data, select [Execute], and press the ENTER button.

For details on each setting, see step 6 of “To select data for copying from the [Recording List] screen” (page 57).

Tip

To cancel copying, select [Abort] and press the ENTER button.

- 10 When copying is complete, press the ENTER button while [Close] is selected.

The [Image List] screen appears again.

Editing Patient Information

You can edit the ID, name, and other patient information that is attached to recorded data and not yet configured.

Notes

- When [Change Patient Info.] is set [Not Used] in the [System Admin Settings] screen – [Function Settings] screen – [Patient Info.] tab, you can only edit patient information that has not yet been entered. Patient information entered by a user cannot be edited.
- Patient information obtained from MWL cannot be edited.

- 1 Display the [Recording List] screen, select the operation box, and press the ENTER button.

The [Operation] screen appears.

- 2 Select [Change], and press the ENTER button.

The [Recording List] screen appears again.

- 3 Select the check box of the recorded data for which you want to edit patient information.

Tip

You can only edit patient information for one recorded data entry at a time. If multiple recorded data entries are selected, patient information editing will be disabled.

- 4 Select [Execute], and press the ENTER button.

The [Change] screen appears.

- 5 Edit each item as necessary, select [Change], and press the ENTER button.

For details on each setting, see “Patient Information and Data Storage Settings” (page 46).

The patient information is changed, and the [Recording List] screen appears again.

Protecting Recorded Data

You can protect data to prevent accidental deletion. Recorded data can be protected or protection can be removed from the [Recording List] screen.

To select data for protection from the [Recording List] screen

This section describes how to select recorded data for protection in the [Recording List] screen. Both videos and still images will be protected with this procedure.

- 1 Display the [Recording List] screen, select the operation box, and press the ENTER button.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

The [Operation] screen appears.


- 2 Select [Protect], and press the ENTER button.

The [Recording List] screen appears again.

- 3 Select the check boxes for the recorded data you want to protect.

- 4 Select [Execute], and press the ENTER button.

The recorded data is protected, and the  icon appears.

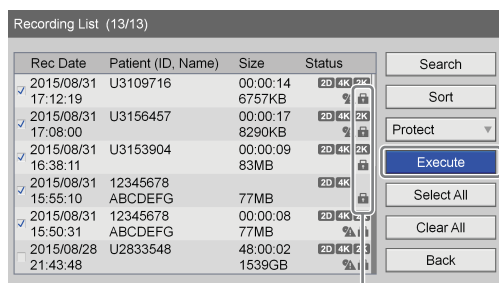

Rec Date	Patient (ID, Name)	Size	Status
✓ 2015/08/31 17:12:19	U3109716	00:00:14 6757KB	2D 4K 5P 2
✓ 2015/08/31 17:08:00	U3156457	00:00:17 8290KB	2D 4K 2K
✓ 2015/08/31 16:38:11	U3153904	00:00:09 83MB	2D 4K 2K
✓ 2015/08/31 15:55:10	12345678 ABCDEFGH	77MB	2D 4K
✓ 2015/08/31 15:50:31	12345678 ABCDEFGH	00:00:08 77MB	2D 4K 5P
2015/08/28 21:43:48	U2833548	48:00:02 1539GB	2D 4K 5P 3

The  icon appears for recorded data that is protected.

To remove protection

Select [UnProtect] in step 2 of the previous procedure, and press the ENTER button. Perform the rest of the procedure as shown.

Deleting Recorded Data

You can delete recorded data stored on the unit's internal hard disk for particular surgical procedures or examinations.

Recorded data can be deleted from the [Recording List] screen.

Note

Recorded data that is deleted cannot be restored.

- 1 Display the [Recording List] screen, select the operation box, and press the ENTER button.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

The [Operation] screen appears.

- 2 Select [Delete], and press the ENTER button.

The [Recording List] screen appears again.

- 3 Select the check boxes for the recorded data you want to delete.

- 4 Select [Execute], and press the ENTER button.

- 5 Select [OK] when the confirmation message appears, and press the ENTER button.

When deletion of the recorded data is complete, a notification message will appear.

To return to the [Recording List] screen, press the ENTER button.

Viewing Information on Recorded Data

You can view information on recorded data stored on the unit's internal hard disk.

- 1 Select [Image List] for the operation box in the [Recording List] screen.
- 2 Select the check box of the recorded data entry you want to view in the [Recording List] screen, and press the ENTER button.

For details on displaying the [Recording List] screen, see “Image Search” (page 54).

The [Image List] screen appears.

- 3** Select [Rec Info], and press the ENTER button.

The information on the recorded data is displayed.

Tip

The sizes displayed in the [Recording List] screen indicate total sizes for Ch1 and Ch2.

Displaying the [System Admin Settings] Screen

You can set the initial configurations for various settings in the [System Admin Settings] screen.

Configurations in the [System Admin Settings] screen must be performed by the system administrator.

- **Language and Time Settings (page 62)**
Configure the display language and the date and time settings.
- **Function Settings (page 64)**
Configure initial settings for recording and general functions.
- **Device Settings (page 67)**
Configure initial settings for connected external devices.
- **Password Settings (page 71)**
Register or edit the password for the system administrator settings.
- **Network Settings (page 72)**
Configure initial settings for the network and server.
- **Doctor List Registration (page 73)**
Register and edit information regarding the doctors performing the surgical procedures and examinations.
- **Editing Cases (page 76)**
Register or edit entries for surgical procedures and examinations.
- **Auto Delete Settings (page 78)**
Configure settings for the auto delete function.
- **DICOM Settings (page 78)**
Configure DICOM settings.

Use the on-screen keyboard to enter text.

For details on use, see “Using the On-Screen Keyboard (Text Entry)” (page 29).

- 1 Press the MENU button.

The [MENU] screen appears.

- 2 Select [Settings], and press the ENTER button.

The [Settings] screen appears.

- 3 Select [System Admin Settings], and press the ENTER button.


Tip

If password protection is enabled, a password entry screen will appear.

- 4 Press the ENTER button when the system administrator settings notification screen appears.

The [System Admin Settings] screen appears

- 5 Select the settings you want to configure, and press the ENTER button.

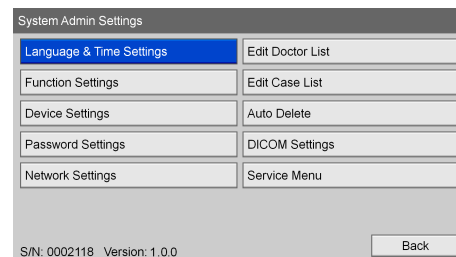

The setting screen for the selected item appears.
The serial number of the unit and the current software version appear at the bottom left.

About the [Service Menu]

The [Service Menu] is a function reserved for service personnel.

It cannot be used by general users.

Language and Time Settings

Configure the display language and the display format for the current date and time.

- 1 Select [Language & Time Settings] in the [System Admin Settings] screen, and press the ENTER button.

The [Language & Time Settings] screen appears.

- 2 Configure each setting in the [Language] tab.

[Language]

Select the language for screen displays.

You can select from [English], [日本語] (Japanese), [Deutsch] (German), [Français] (French), [Italiano] (Italian), and [Español] (Spanish).

[Time Zone]

Select the time zone.

Only the difference from GMT (Greenwich Mean Time) is displayed for each time zone setting.

For details on time zones, see “List of time zones” (page 62).

[Use Daylight Saving Time]

Select whether to enable DST (daylight saving time).

The default setting is [Not Used].

- 3 Configure each setting in the [Date & Time] tab.

[Date Format]

Select the year, month, and day display format.

You can select from [YYYY/MM/DD] (year/month/day), [MM/DD/YYYY] (month/day/year), and [DD/MM/YYYY] (day/month/year).

The default setting is [MM/DD/YYYY].

[Date]

Configure the current date.

[Time]

Configure the current time.

- 4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

List of time zones

Time zone		DST
UTC-12:00	International Date Line West	
UTC-11:00	Coordinated Universal Time-11	
UTC-10:00	Hawaii	
UTC-09:00	Alaska	○
UTC-08:00	Baja California	○
UTC-08:00	Pacific Time (US and Canada)	○
UTC-07:00	Arizona	
UTC-07:00	Chihuahua, La Paz, Mazatlan	○
UTC-07:00	Mountain Time (US and Canada)	○
UTC-06:00	Guadalajara, Mexico City, Monterrey	○
UTC-06:00	Saskatchewan	
UTC-06:00	Central America	
UTC-06:00	Central Time (US and Canada)	○
UTC-05:00	Indiana (East)	
UTC-05:00	Bogota, Lima, Quito, Rio Branco	
UTC-05:00	Eastern Time (US and Canada)	○
UTC-04:30	Caracas	
UTC-04:00	Asuncion	○
UTC-04:00	Cuiaba	○
UTC-04:00	Santiago	
UTC-04:00	Georgetown, La Paz, Manaus, San Juan	
UTC-04:00	Atlantic Time (Canada)	○
UTC-03:30	Newfoundland	○
UTC-03:00	Cayenne, Fortaleza	
UTC-03:00	Greenland	○
UTC-03:00	Salvador	
UTC-03:00	Buenos Aires	○
UTC-03:00	Brasilia	○
UTC-03:00	Montevideo	○
UTC-02:00	Coordinated Universal Time-02	○
UTC-01:00	Azores	○
UTC-01:00	Cape Verde Is.	
UTC	Casablanca	
UTC	Dublin, Edinburgh, Lisbon, London	○
UTC	Monrovia, Reykjavik	
UTC	Coordinated Universal Time	
UTC+01:00	Amsterdam, Berlin, Bern, Rome, Stockholm, Vienna	○
UTC+01:00	Windhoek	○
UTC+01:00	Sarajevo, Skopje, Warsaw, Zagreb	○
UTC+01:00	Brussels, Copenhagen, Madrid, Paris	○

Time zone		DST
UTC+01:00	Belgrade, Bratislava, Budapest, Ljubljana, Prague	○
UTC+01:00	West Central Africa	
UTC+02:00	Athens, Bucharest	○
UTC+02:00	Amman	○
UTC+02:00	Istanbul	○
UTC+02:00	Jerusalem	○
UTC+02:00	Cairo	○
UTC+02:00	Kaliningrad (RTZ 1)	
UTC+02:00	Damascus	○
UTC+02:00	Tripoli	
UTC+02:00	Harare, Pretoria	
UTC+02:00	Beirut	○
UTC+02:00	Helsinki, Kiev, Riga, Skopje, Sofia, Tallinn, Vilnius	○
UTC+02:00	Eastern Europe	○
UTC+03:00	Kuwait, Riyadh	
UTC+03:00	Nairobi	
UTC+03:00	Baghdad	
UTC+03:00	Minsk	
UTC+03:00	Moscow, St. Petersburg, Volgograd (RTZ 2)	○
UTC+03:30	Tehran	○
UTC+04:00	Abu Dhabi, Muscat	
UTC+04:00	Izhevsk, Samara (RTZ 3)	
UTC+04:00	Yerevan	○
UTC+04:00	Tblisi	
UTC+04:00	Baku	○
UTC+04:00	Port Louis	○
UTC+04:30	Kabul	
UTC+05:00	Ashgabat, Tashkent	
UTC+05:00	Islamabad, Karachi	○
UTC+05:00	Yekaterinburg (RTZ 4)	
UTC+05:30	Sri Jayawardenepura	
UTC+05:30	Chennai, Kolkata, Mumbai, New Delhi	
UTC+05:45	Kathmandu	
UTC+06:00	Astana	
UTC+06:00	Dhaka	
UTC+06:00	Novosibirsk (RTZ 5)	○
UTC+06:30	Yangon (Rangoon)	
UTC+07:00	Krasnoyarsk (RTZ 6)	○
UTC+07:00	Bangkok, Hanoi, Jakarta	
UTC+08:00	Irkutsk (RTZ 7)	
UTC+08:00	Ulaanbaatar	○
UTC+08:00	Kuala Lumpur, Singapore	
UTC+08:00	Perth	○

Time zone		DST
UTC+08:00	Taipei	
UTC+08:00	Beijing, Chongqing, Hong Kong SAR, Urumqi	
UTC+09:00	Seoul	
UTC+09:00	Yakutsk (RTZ 8)	○
UTC+09:00	Osaka, Sapporo, Tokyo	
UTC+09:30	Adelaide	○
UTC+09:30	Darwin	
UTC+10:00	Vladivostok, Magadan (RTZ 9)	○
UTC+10:00	Canberra, Melbourne, Sydney	○
UTC+10:00	Guam, Port Moresby	
UTC+10:00	Brisbane	
UTC+10:00	Hobart	○
UTC+10:00	Magadan	
UTC+11:00	Solomon Islands, New Caledonia	
UTC+11:00	Chokurdakh (RTZ 10)	
UTC+12:00	Anadyr, Petropavlovsk-Kamchatsky (RTZ 11)	
UTC+12:00	Auckland, Wellington	○
UTC+12:00	Fiji	
UTC+12:00	Coordinated Universal Time+12	
UTC+13:00	Samoa	○
UTC+13:00	Nuku'alofa	

Function Settings

Configure initial settings for patient information, input signals, and other functions.

- 1 Select [Function Settings] in the [System Admin Settings] screen, and press the ENTER button.

The [Function Settings] screen appears.

- 2 Select a tab, and press the ENTER button.
- 3 Configure the settings in each tab as necessary.

For details on each tab, see the following.

[Patient Info.] Tab (page 64)

[Input Detection] Tab (page 64)

[Input Signal] Tab (page 64)

[Auto Live] Tab (page 65)

[General 1] Tab (page 65)

[General 2] Tab (page 66)

[Facility] Tab (page 66)

[CMS] Tab (page 67)

- 4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

[Patient Info.] Tab

Configure initial settings for patient information.

[Patient ID Length]

Select the number of digits to use for patient IDs.
You can select from 4 to 16 digits.
The default setting is [8].

[Add 0 to Patient ID]

Select whether to add zeros to patient IDs with less than the specified number of digits.
The default setting is [Use].

[Patient Name Field]

Select the number of display fields to use for the patient name.
If you select [2] or [3], the patient name field will be separated into first, middle, and last names.
The default setting is [1].

[Change Patient Info.]

Select whether to allow changes to patient information at any time.
The default setting is [Not Used].

[Input Detection] Tab

Configure the detection method for input signals.

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, the items in this tab will be grayed out.

[2D/3D]

Select whether to recognize the 2D/3D format of input signals automatically or whether to specify it manually.
The default setting is [Auto].

[Color Gamut]

Select whether to recognize the color gamut of input signals automatically or whether to specify it manually.
The default setting is [Auto].

[Input Signal] Tab

Configure initial settings for input signals.

[Region]

You can select from [NTSC] and [PAL].
The default setting is [NTSC].

[I/O Timing]

Select from [Normal] or [Pass-Through] for the input and output timing.
Select [Normal] as signal processing is usually performed under normal circumstances. There will be a slight delay in signal outputs when [Normal] is selected.
Select [Pass-Through] for simultaneous input and output timing without signal processing.
The default setting is [Normal].

[Input Detect Duration]

Select from 0 to 3 seconds for the duration for which input signals are detected.
The default setting is [3 sec].

Note

Image distortion may occur if a duration of 2 seconds or less is selected.

[3D Line by Line Signal]

When inputting 3D line-by-line signals, select [L First] (left image) or [R First] (right image) for the first line.
The default setting is [L First].

Tip

When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this item will be grayed out.

[Auto Live] Tab

Configure initial settings for the auto live function.

[Auto Live]

Select whether to use the auto live function.

When the auto live function is enabled and you capture a still image, the captured image is displayed for the duration specified with the [Duration] setting of [Auto Live Settings].

When using the auto live function, you can select the number of images to display during use ([Single] or [Multi]).

When [Multi] is selected, three or four captured images are displayed according to the display position setting. The default setting is [Not Used].

[Advanced] appears when [Single] or [Multi] is selected. For details on this setting, see “To configure auto live settings” (page 65).

Tip

If [Auto Live] is set to [Use] while [I/O Timing] is set to [Pass-Through], the auto live function will not function.

[Auto Live Still Image]

Displays the type of image that is displayed when a video is paused while using the auto live function.

This is [Frame] on this unit.

To configure auto live settings

Configure display duration and position settings for the captured images that are displayed when using the auto live function.

- 1 Select [Advanced], and press the ENTER button.

The [Auto Live Settings] screen appears.

- 2 Configure each setting.

[Duration]

Select the duration for which captured still images are displayed.

When [Single] is selected, select from 0.5 to 2 seconds for the duration. The default setting is [0.5 sec].

When [Multi] is selected, select from 1 to 5 seconds for the duration. The default setting is [3 sec].

[Position]

Select from [Top], [Bottom], [Left], and [Right] for the position in which captured still images are displayed.

The default setting is [Left].

The four most recent still images are displayed when [Top] or [Bottom] is selected, and the three most recent still images are displayed when [Left] or [Right] is selected.

- 3 When you finish configuration, select [Apply] and press the ENTER button.

The [Auto Live] tab appears again.

[General 1] Tab

Configure initial system settings.

[Show Status]

Select the check box for the status you want to display.

You can select from [Rec Status], [Play Status], and [Idle Status].

None of the check boxes are selected under default settings.

Tip

If [Show Rec Status] or [Show Idle Status] is set to [Use] while [I/O Timing] is set to [Pass-Through], the status information will not be displayed. In addition, status information will not be updated during the auto live function.

[Image List Display]

Select the category of thumbnails to display for videos in the [Image List] screen.

The default setting is [Title].

For details on [Advanced] settings, see “To configure advanced image list settings” (page 66).

[Auto Eject BD/DVD]

Select whether to automatically eject BD/DVD discs when they are full.

The default setting is [Not Used].

[Prohibit Remote Control]

Select whether to prohibit control of the unit via the remote control buttons.

The default setting is [Not Used].

[Menu Display (HDMI)]

Select whether to display menus on the video output from the HDMI connector.

The default setting is [Not Used].

Tips

- When the output resolution is SD, menus will not be displayed.
- When the output resolution is SD or 1920 × 1080, the idle status, recording status, playback status, and multiple auto live images will not be displayed.
- If [Menu Display (HDMI)] is set to [Use], input images from SDI output connectors A to D will not be displayed and black screens will be output instead.
- When [Ch1/Ch2 Rec Mode] is set to [Ch1 2K Rec] in the [Recording] tab, this item will be grayed out.

To configure advanced image list settings

Specify the selection state of recorded data when you perform copy or print operations in the [Image List] screen.

- 1 Select [Advanced] for [Image List Display], and press the ENTER button.

The [Advanced] settings screen appears.

- 2 Configure each setting, select [Apply] and press the ENTER button.

[Save to]

Select the selection state of recorded data when [Save to] is selected in the operation box.

You can select from [Deselect All], [Select All], [Video], and [Still].

The default setting is [Deselect All].

[Print]

Select the selection state of recorded data when [Print] is selected in the operation box.

You can select from [Deselect All] and [Still].

The default setting is [Deselect All].

[General 2] Tab

Configure initial system settings.

[Recovery]

Select whether to return to the action that was being performed when a recovery operation was executed.

The default setting is [Use].

Recovery of recorded data is the only action performed when [Not Used] is selected.

[Automatic Pause]

Select whether to automatically pause recording when there is no input signal.

The default setting is [Not Used].

If you select [Use], recording will pause if there is no input signal for 10 minutes.

[Mic Plug-In Power]

Select the voltage to use when supplying power to the microphone from the unit.

You can select from [2.5V], [3.2V], and [4.2V].

This item is available when [Mic Plug-In Power] is set to [Use] in the [User Settings] screen – [Select Input] tab.

The default setting is [2.5V].

[LCD Backlight Adjust]

Select the brightness of the front panel display.

You can select from [25%], [50%], [75%], and [100%].

The default setting is [75%].

[Prohibit USB/BD/DVD]

Select whether to prohibit transmission to external storage destinations other than servers.

The default setting is [Not Used].

The default state is transfer enabled.

[Facility] Tab

Configure initial settings for the facility.

[Facility Name]

Enter up to 32 characters for the facility name.

Use only alphanumeric characters, underscores (_), hyphens (-), and spaces.

[Display Logo]

Select whether to display a logo in the [MENU] screen.

The default setting is [Not Used].

If you want to set [Display Logo] to [Use], be sure to import the logo.

For details on importing, see “To import a logo” (page 66).

To import a logo

Import the logo file you want to display.

- 1 Copy the logo file onto a USB memory device, and insert the device into a USB port 1 on the front panel of the unit.
- 2 Select [Import], and press the ENTER button.
- 3 Press the ENTER button when the message that requests media to be inserted appears.

When import of the logo file is complete, a notification message will appear.

To prepare a logo file

To display a logo, prepare a logo file that can be imported to the unit and copy it onto a USB memory device beforehand.

1 Prepare the logo file.

The logo file must meet the following conditions.

File format: 24-bit bitmap

File name: MenuLogo.bmp

Size: 43 (height) × 400 (width) pixels

Background: RGB (102, 102, 102) (gray)

2 Copy the logo file onto a USB memory device.

Create a folder with the following name in the uppermost directory of the USB memory device, and copy the logo file.

Folder name: DisplayLogo

[CMS] Tab

Configure initial settings for the CMS mode. The CMS mode is used for transferring recorded data to the content management system.

[CMS Mode]

Select whether to use the CMS mode.

The default setting is [Not Used].

Tip

For details on the content management system, contact your local Sony representative.

Device Settings

Configure settings for external devices.

1 Select [Device Settings] in the [System Admin Settings] screen, and press the ENTER button.

The [Device] screen appears.

2 Select a tab, and press the ENTER button.

3 Configure the settings in each tab as necessary.

For details on each tab, see the following.

[Device 1] Tab (page 67)

[Device 2] Tab (page 69)

[Control] Tab (page 69)

4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

[Device 1] Tab

Configure settings for the USB devices connected to the unit.

For details on devices that can be connected to the unit, contact your local Sony representative.

[USB 3] to [USB 6]

Select the devices connected to USB ports 3 to 6 on the rear panel of the unit.

The available settings include devices that are compatible with the unit.

The default setting is [No Device].

If you select a printer, be sure to configure the [Printer Settings].

For details on this setting, see “To configure printer settings” (page 68).

Tips

- The unit can recognize the Sony UP-DR80MD printer automatically. To have the unit recognize UP-DR80MD printers automatically, set the USB port to which the printer will be connected to [No Device] in the device settings. Be aware that the printer will not be recognized automatically while recording or playback

is in progress or while the [Device Settings] screen is displayed.

- If the USB port to which a card reader is connected is changed while the unit is turned on, the card reader will not be recognized. In such cases, restart the unit.

If you selected [Bar Code Reader] or [Card Reader], be sure to configure the [Reader Settings].

For details on this setting, see “To configure reader settings” (page 69).

If you selected [Touch Panel / Mouse], be sure to configure the [Touch Panel / Mouse Settings].

For details on this setting, see “Configuring Touch Panel / Mouse Settings” (page 81).

To configure printer settings

When a UP-DR80MD is selected, configure the color balance, brightness, and other settings.

- 1 Select [Advanced] for the USB port to which the UP-DR80MD is connected, and press the ENTER button.

[Advanced] appears when one of [USB 3] to [USB 6] is set to [UP-DR80MD].

The [Printer Settings] screen appears

- 2 Display the [Color Balance] tab, and adjust the color balance for cyan, magenta, and yellow and red, green, and blue.

You can specify a value between -32 to +32 for each color.

- 3 Display the [Gray Balance] tab, and adjust the balance for achromatic colors.

You can specify a value between -50 to +50 for the X and Y values.

- 4 Display the [Bright] tab, and configure the [Sharpness], [Dark], [Gamma], and [Light] settings to adjust the brightness of printed images.

You can specify a value between -3 to +3 for [Sharpness], and a value between -32 to +32 for [Dark], [Gamma], and [Light].

- 5 Display the [Import] tab, and import a watermark file or logo file.

[Import Watermark]

Import a watermark file.

Copy the watermark file onto a USB memory device, and insert the device into a USB port on the unit.

For details on preparing a watermark file, see “To prepare a watermark file” (page 68).

When you select [Import Watermark] and press the ENTER button, a message requesting that you insert media appears. Select [OK] and press the ENTER button.

When import of the watermark file is complete, a notification message will appear.

[Import Logo]

Import a logo file.

Copy the logo file onto a USB memory device, and insert the device into a USB port on the unit.

For details on preparing a logo file, see “To prepare a logo file” (page 69).

When you select [Import Logo] and press the ENTER button, a message requesting that you insert media appears. Select [OK] and press the ENTER button.

When import of the logo file is complete, a notification message will appear.

- 6 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [Device 1] tab appears again.

To prepare a watermark file

To print watermarks, prepare a watermark file that can be imported to the unit and copy it onto a USB memory device beforehand.

- 1 Prepare the watermark file.

The watermark file must meet the following conditions.

• For A4 size paper

File format: 24-bit bitmap

File name: watermark_A4.bmp

Size: 2392 × 3400 pixels

Color:

Background: RGB (255, 255, 255) (white)

Watermark: RGB (0, 0, 0) (black)

• For letter size paper

File format: 24-bit bitmap

File name: watermark_85x11.bmp

Size: 2464 × 3192 pixels

Color:

Background: RGB (255, 255, 255) (white)

Watermark: RGB (0, 0, 0) (black)

- 2 Copy the watermark file onto a USB memory device.

Create a folder with the following name in the uppermost directory of the USB memory device, and copy the watermark file.

Folder name: Watermark

To prepare a logo file

To print logos, prepare a logo file on the unit and copy it onto a USB memory device beforehand.

- 1 Prepare the logo file.

The logo file must meet the following conditions.

File format: 24-bit bitmap

File name: printlogo_UP-DR80MD.bmp

Size: 94 (height) × 900 (width) pixels

- 2 Copy the logo file onto a USB memory device.

Create a folder with the following name in the uppermost directory of the USB memory device, and copy the logo file.

Folder name: PrintLogo

To configure reader settings

If you selected [Bar Code Reader] or [Card Reader], configure the settings for reading data. The procedure for configuring barcode reader and card reader settings is the same.

- 1 Select [Advanced] for the USB port to which the barcode reader or card reader is connected, and press the ENTER button.

[Advanced] appears when [Bar Code Reader] or [Card Reader] is selected.

The [Reader Settings] screen appears.

- 2 Display the [Patient ID/Name] tab, and configure each setting.

[Patient ID]

Specify the starting position and length for the patient ID.

[Patient Name]

Specify the starting position and length for the patient name.

- 3 Display the [Gender/DOB] tab, and configure each setting.

[Gender]

Specify the starting position, male, and female for the gender.

[Birth Date]

Specify the starting position for the date of birth and the date format.

You can select from [YYYY/MM/DD], [MM/DD/YYYY], [DD/MM/YYYY], [YYYYMMDD], [MMDDYYYY], and [DDMMYYYY], or leave the field empty for the date format.

- 4 Scan the card.

The data from the card appears in the [Read Data] area of the [Decode] tab.

Tip

If there is an empty field within the data, an asterisk (*) appears on the display.

- 5 Select [Decode], and press the ENTER button.

The decoding results are displayed in the [Decode] tab.

- 6 After confirming that the decoding results are correct, select [Apply] and press the ENTER button.

The configurations are saved, and the [Device 1] tab appears again.

[Device 2] Tab

Configure settings for the printer and keyboard connected to the unit.

[Printer Used]

Select the printer to be used.

This is fixed at [UP-DR80MD] on this unit.

[Keyboard Mode]

Select whether to use the soft keyboard (i.e., on-screen keyboard) or a hard keyboard in the [Recording Preparation] screen.

The default setting is [Soft Keyboard].

You can configure function key settings in [Advanced]. For details, see "Using Function Keys" (page 38).

[Control] Tab

Configure settings related to performing system controls for the unit from an external device.

[USB 7]

Specify the device connected to USB port 7.

The default setting is [No Device].

If you selected [System Control], you can verify the transmission speed under [System Control Setting].

For details on verifying the setting, see “To configure system control settings” (page 70).

[RS-232C]

Specify the device connected to the RS-232C connector. The default setting is [No Device].

If you selected [System Control], be sure to configure the transmission speed under [System Control Setting].

For details on this setting, see “To configure system control settings” (page 70).

[Remote Startup]

Select whether to link the power on/off operations of remote systems via USB port 7 and the RS-232C connector.

The default setting is [Do Not Start].

When [USB 7] or [RS-232C] is selected, the remote system's on/off operation will be linked to that of the unit via the selected connector.

To configure system control settings

Display or configure the transmission speed when performing system controls via USB port 7 or the RS-232C connector.

- 1 Select [Advanced] for [USB 7] or [RS-232C], and press the ENTER button.

[Advanced] appears when [USB 7] or [RS-232C] is set to [System Control].

The [System Control Setting] screen appears.

- 2 Select the transmission speed, select [Apply] and press the ENTER button.

For USB port 7, the transmission speed is fixed at 115200 bps.

For the RS-232C connector, you can select from 9600 bps, 19200 bps, 38400 bps, 57600 bps, and 115200 bps for the transmission speed.

The default setting is [9600 bps].

The configurations are saved, and the [Control] tab appears again.

[Contact Switch] Tab

Configure these settings when you want to control the unit via contact switch devices connected to the contact switch connectors on the rear panel of the unit.

- 1 Select the control items for [Contact Switch 1] to [Contact Switch 3].

Specify what will be controlled via the contact switch device for each contact switch.

The procedure for configuring [Contact Switch 1] to [Contact Switch 3] is the same.

The default setting is [Not Used].

[Not Used]

Disable use of the contact switch device.

[Record Video/Pause]

Control pause and resume functions while recording videos.

[Record Still]

Capture still images.

[Audio On/Off]

Enable or disable audio output.

[Set Index]

Add an index without saving a still image of the recording.

- 2 Select [Advanced] for the contact switch that will perform the control, and press the ENTER button.

The [Advanced] settings screen appears.

- 3 Specify the control timing, select [Apply] and press the ENTER button.

You can select from [Falling Edge], [Leading Edge], and [State].

The default setting is [Falling Edge].

[State] cannot be selected when the control item is set to [Record Still] or [Set Index].

Tip

The pause control of the [Record Video/Pause] function will not operate immediately after recording starts. Wait about 5 seconds before performing the control.

The [Contact Switch] tab appears again.

Password Settings

Configure password settings to prevent unauthorized changes to the [System Admin Settings].
When password protection is enabled, a password entry screen appears when [System Admin Settings] is selected in the [Settings] screen.

Note

The password is not configured under default conditions. When configuring the password for the first time, configure the new password in the [Change Password] screen without entering an old password. Be sure to write down the password in case you forget it. The [System Admin Settings] cannot be accessed without the password. Guard the password securely to maintain security.

To enable password protection

- 1 Select [Password Settings] in the [System Admin Settings] screen, and press the ENTER button.
The [Password Settings] screen appears.
- 2 Select the password lock box, and press the ENTER button.
- 3 Select [Use], and press the ENTER button.
The default setting is [Not Used].
- 4 Select [Settings], and press the ENTER button.
The configuration is saved.

To change the password

- 1 Select [Change] in the [Password Settings] screen, and press the ENTER button.
The [Change Password] screen appears.
- 2 Configure each setting, select [Apply] and press the ENTER button.

[Old Password]

Enter the old password.

[New Password]

Enter up to 255 characters for the new password. Any character that can be entered on the on-screen keyboard can be used in the password.

[Check New Password]

Reenter the new password.

The new password is saved.

Network Settings

Configure initial settings for the network and server.

- 1 Select [Network Settings] in the [System Admin Settings] screen, and press the ENTER button.

The [Network Settings] screen appears.

- 2 Select a tab, and press the ENTER button.

- 3 Configure the settings in each tab as necessary.

For details on each tab, see the following.

[Network] Tab (page 72)

[IP Address] Tab (page 72)

[DNS Server] Tab (page 72)

[File Server] Tab (page 72)

[Shared] Tab (page 73)

[NTP] Tab (page 73)

- 4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

[Network] Tab

Configure these settings when using a network.

[Network]

Select whether to use a network.

The default setting is [Not Used].

If you selected [Use], configure the settings in each tab as necessary.

[IP Address] Tab

[Obtain an IP address automatically]

Select [Use] to use the DHCP to obtain an IP address automatically, or select [Not Used] to assign a static IP address manually.

If you selected [Use], the IP address, subnet mask, and default gateway will be displayed as read-only.

If you selected [Not Used], be sure to configure the [IP Address], [Subnet Mask], and [Default Gateway] settings.

[IP Address]

Enter the IP address.

[Subnet Mask]

Enter the subnet mask.

[Default Gateway]

Enter the IP address for the default gateway.

If you are only using a local network and not connecting to other networks, do not configure this setting.

[DNS Server] Tab

[Obtain DNS server address automatically]

Select [Use] to obtain the IP address for the DNS server automatically, or select [Not Used] to specify the IP address manually.

If you selected [Not Used], be sure to configure the [Preferred DNS server] and [Alternate DNS server] settings.

[Preferred DNS server]

Enter the IP address for the preferred DNS server.

If you do not have nor need a preferred DNS server, do not configure this setting.

[Alternate DNS server]

Enter the IP address for the alternate DNS server.

If you do not have nor need an alternate DNS server, do not configure this setting.

[File Server] Tab

Configure these settings when specifying a server as the storage destination for recorded data.

[Server]

Select [Common Setting (Share)] to specify a server as the storage destination.

The default setting is [Not Used].

If you selected [Common Setting (Share)], be sure to configure the share name and other settings in the [Shared] tab.

For details on this setting, see “[Shared] Tab” (page 74).

[Server Name]

Enter the server name.

[Share Name]

Enter the share name.

Use only alphanumeric characters and symbols (excluding " = / \ [] : | < > + ; , ? *).

[Shared] Tab

Configure the necessary settings for using a server as the storage destination.

[User Name]

Enter a user name.

[Password]

Enter the password.

[NTP] Tab

Configure these settings to obtain the current time from the NTP server.

[NTP]

Select [Use] to obtain the current time from the NTP server.

The default setting is [Not Used].

If you selected [Use], be sure to configure the [Address] and [Port No.] settings.

[Address]

Enter the address for the NTP server.

[Port No.]

Enter the port number.

The default setting is [123].

Doctor List Registration

Register doctor names and initial settings for each doctor. The doctor information registered here will be used in the doctor lists and default settings that appear in other screens. You can register up to 100 doctors.

Registering Doctors

- 1 Select [Edit Doctor List] in the [System Admin Settings] screen, and press the ENTER button.

The [Edit Doctor List] screen appears.

- 2 Select [Add], and press the ENTER button.

The [Add Doctor] screen appears.

- 3 Select a tab, and press the ENTER button.

- 4 Configure the settings in each tab as necessary.

For details on each tab, see the following.

[General] Tab (page 73)

[File Server] Tab (page 74)

[Shared] Tab (page 74)

[Print] Tab (page 74)

[External Media] Tab (page 74)

- 5 When you finish configuration, select [Add] and press the ENTER button.

The doctor information is added.

[General] Tab

Configure general settings such as the doctor name and storage destination.

[No.]

Displays the doctor number.

A sequential number is assigned automatically and cannot be changed.

[Name]

Enter up to 64 characters for the doctor name.

[Video Quality]

Select the video recording quality.

You can select from [Best], [High], and [Standard].

[External Media]

Select the external storage media that will be used during simultaneous recording.

You can select from [Not Specified], [Auto-Detect], [BD/DVD], [USB 1], [USB 2], and [SERVER].

[File Server] Tab

Configure these settings when specifying a server as the storage destination.

[Server]

Select [System Settings] or [Share (CIFS)] for the server type.

The default setting is [Not Used].

If you selected [Share (CIFS)], be sure to configure the necessary settings in the [Shared] tab, in addition to the [Server Name].

For details on this setting, see “[Shared] Tab” (page 74).

[Server Name]

Enter the server name.

[Share Name]

Enter the share name.

Use only alphanumeric characters and symbols (excluding " = / \ [] : | < > + ; , ? *).

[Shared] Tab

Configure the necessary settings for when [Server] is set to [Share (CIFS)] in the [File Server] tab.

[User Name]

Enter a user name.

[Password]

Enter the password.

[Print] Tab

[Printer]

Displays the name of the connected printer.

[Layout]

Select the number of still images to print on each sheet.

To use the system settings, select [User Settings].

For the portrait orientation

You can select from 1 image (1 × 1), 2 images (1 × 2), 4 images (2 × 2), 6 images (2 × 3), 8 images (2 × 4), 9 images (3 × 3), 12 images (3 × 4), 15 images (3 × 5), or 18 images (3 × 6).

The default setting is [User Settings].

For the landscape orientation

You can select from 1 image (1 × 1), 2 images (2 × 1), 4 images (2 × 2), 6 images (3 × 2), 8 images (4 × 2), 9 images (3 × 3), 12 images (4 × 3), 15 images (5 × 3), or 18 images (6 × 3).

The default setting is [User Settings].

[External Media] Tab

Configure save settings for videos and still images.

[Save Video]

Select the check box of the format in which videos will be saved.

All of the check boxes are selected under default settings.

[Save Still]

Select the check box of the format in which still images will be saved.

All of the check boxes are selected under default settings.

Editing the Doctor List

- 1 Select [Edit Doctor List] in the [System Admin Settings] screen, and press the ENTER button.

The [Edit Doctor List] screen appears.

- 2 Edit the doctor information.

Sorting the Doctor List (page 75)

Editing Doctor Settings (page 75)

Deleting Doctors (page 75)

- 3 Select [Settings], and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

Sorting the Doctor List

- 1 Select the doctor you want to reposition in the [Edit Doctor List] screen, select [Up] or [Down], and then press the ENTER button.

The doctor entry moves up or down one position at a time.

Editing Doctor Settings

- 1 Select the doctor for which you want to change settings in the [Edit Doctor List] screen, select [Edit], and then press the ENTER button.

The [Edit Doctor] screen appears.

- 2 After editing the settings, select [Apply] and press the ENTER button.

The configuration is saved.

For details on each setting, see the following.

[General] Tab (page 73)

[File Server] Tab (page 74)

[Shared] Tab (page 74)

[Print] Tab (page 74)

[External Media] Tab (page 74)

Note

If you change the doctor name, you will not be able to use the previous doctor name to search for recorded data attributed to that name.

Deleting Doctors

- 1 Select the doctor you want to delete in the [Edit Doctor List] screen, select [Delete], and then press the ENTER button.
- 2 Select [OK] when the confirmation message appears, and press the ENTER button.

The doctor entry is deleted.

Editing Cases

Edit cases and categories for surgical procedures and examinations, and create the list of cases that will be used in the various screens.

- 1 Select [Edit Case List] in the [System Admin Settings] screen, and press the ENTER button.

The [Edit Case List] screen appears.

- 2 Edit the case information.

Registering Cases (page 76)

Deleting Cases (page 76)

Sorting the Case List (page 76)

Editing Case Names (page 76)

Registering Categories (page 76)

Selecting Categories (page 77)

Editing Categories (page 77)

- 3 Select [Apply] in the [Edit Case List] screen, and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

Registering Cases

Up to 100 case entries can be registered for a single category.

- 1 Select [Add] in the [Case] area of the [Edit Case List] screen, and press the ENTER button.

The [Add] screen appears.

- 2 Enter up to 20 characters for the case name, select [Add], and then press the ENTER button.

The case is added to the list in the [Edit Case List] screen.

Deleting Cases

- 1 Select the case you want to delete in the [Edit Case List] screen, select [Delete] in the [Case] area, and then press the ENTER button.

- 2 Select [OK] when the confirmation message appears, and press the ENTER button.

Sorting the Case List

- 1 Select the case you want to reposition in the [Edit Case List] screen, select [Up] or [Down] in the [Case] area, and then press the ENTER button.

The case entry moves up or down one position at a time.

Editing Case Names

- 1 Select the case you want to rename in the [Edit Case List] screen, select [Edit] in the [Case] area, and then press the ENTER button.

The [Edit] screen appears.

- 2 Rename the case, select [Apply], and then press the ENTER button.

The configurations are saved, and the [Edit Case List] screen appears again.

Registering Categories

- 1 Select [Edit] in the [Category] area of the [Edit Case List] screen, and press the ENTER button.

The [Edit Category List] screen appears.

- 2 Select [Add], and press the ENTER button.

The [Add] screen appears.

- 3 Enter up to 20 characters for the category name, select [Add], and then press the ENTER button.

A sequential number is assigned to [No.] automatically and cannot be changed.

- 4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [Edit Category List] screen appears again.

To delete a category

- 1 Select the category you want to delete in the [Edit Category List] screen, select [Delete], and then press the ENTER button.
- 2 Select [OK] when the confirmation message appears, and press the ENTER button.

Tip

If you delete a category, all procedures registered to that category will also be deleted.

Selecting Categories

- 1 Select the box in the [Category] area of the [Edit Case List] screen, and press the ENTER button.

The case category selection screen appears.
- 2 Select the case category, and press the ENTER button.

The configurations are saved, and the [Edit Case List] screen appears again.

Editing Categories

- 1 Select [Edit] in the [Category] area of the [Edit Case List] screen, and press the ENTER button.

The [Edit Category List] screen appears.
- 2 Edit the categories.

Sorting the Category List (page 77)
Editing Category Names (page 77)
- 3 Select [Apply] in the [Edit Category List] screen, and press the ENTER button.

The [Edit Case List] screen appears again.

Sorting the Category List

- 1 Select the category you want to reposition in the [Edit Category List] screen, select [Up] or [Down], and then press the ENTER button.

The category moves up or down one position at a time.

Editing Category Names

- 1 Select the category you want to rename in the [Edit Category List] screen, select [Edit], and then press the ENTER button.

The [Edit] screen appears.
- 2 Rename the category, select [Apply], and then press the ENTER button.

The configurations are saved, and the [Edit Category List] screen appears again.

Auto Delete Settings

Select whether to use the auto delete function, and configure the threshold value that will be used for the function.

- 1 Select [Auto Delete] in the [System Admin Settings] screen, and press the ENTER button.

The [Auto Delete Setting] screen appears.

- 2 Configure each setting, select [Apply] and press the ENTER button.

[Auto Delete]

Select whether to use the auto delete function.
The default setting is [Use].

[Auto Delete Value]

Select the threshold value that will be used for the auto delete function.

You can select from [10%], [20%], and [30%].
The default setting is [10%].

The configuration is saved.

DICOM Settings

Configure settings that are necessary for DICOM transmissions of still images.

- 1 Select [DICOM Settings] in the [System Admin Settings] screen, and press the ENTER button.

The [DICOM Settings] screen appears.

- 2 Select a tab, and press the ENTER button.

- 3 Configure the settings in each tab as necessary.

For details on each tab, see the following.

[MWM Server 1] Tab (page 78)

[MWM Server 2] Tab (page 79)

[C-Store] Tab (page 79)

[Local Station 1] Tab (page 80)

[Local Station 2] Tab (page 80)

- 4 When you finish configuration, select [Apply] and press the ENTER button.

The configurations are saved, and the [System Admin Settings] screen appears again.

[MWM Server 1] Tab

Configure MWM server settings.

[MWM]

Select whether to use the MWM server.
The default setting is [Not Used].

[Confirm] appears when [Use] is selected. Confirm communication with the MWM server. For details, see “To confirm communication with the MWM server” (page 79)

[Server Name]

Enter the server name.
This item is available when [MWM] is set to [Use].
Up to 253 characters can be entered.

[Port No.]

Enter the port number.
This item is available when [MWM] is set to [Use].
You can enter a value ranging from 1 to 65535 (up to 5 digits).
The default setting is [104].

[AE Title]

Enter the AE title.

This item is available when [MWM] is set to [Use].

Up to 16 characters can be entered.

[Retry]

Select whether to retry operations, such as accessing or searching MWL and confirming communication with the MWM server.

The default setting is [Use].

[Advanced] appears when [Use] is selected. For details on the settings, see “To configure the retry interval and number of retries” (page 79).

To confirm communication with the MWM server

Confirm whether communication with the MWM server is enabled.

- 1 After configuring each setting in the [MWM Server 1] tab, select [Confirm] for [MWM] and press the ENTER button.

Communication confirmation starts. An indicator appears during confirmation.

The results of confirmation appear when the process is complete.

- 2 Select [OK], and press the ENTER button.

Tip

If proper communication with the MWM server was not established, refer to the cause displayed on the screen and resolve the problem.

To configure the retry interval and number of retries

- 1 Select [Advanced] for [Retry], and press the ENTER button.

The [Advanced] settings screen for the retry appears.

- 2 Configure each setting, select [Apply] and press the ENTER button.

[Retry Interval (sec)]

Enter a value ranging from 1 to 120 seconds for the retry interval.

The default setting is 10 seconds.

[Number of Retries]

Enter a value ranging from 1 to 10 for the number of retries.

The default setting is 10.

[MWM Server 2] Tab

Configure MWM server settings.

[Patient Name Component]

Select the patient name component to be obtained.

The default setting is [1st Component].

[C-Store] Tab

Configure C-Store settings.

[C-Store]

Select whether to use C-Store.

The default setting is [Not Used].

[Confirm] appears when [Use] is selected. Confirm communication with the C-Store server. For details, see “To confirm communication with the C-Store server” (page 79).

[Server Name]

Enter the server name.

This item is available when [C-Store] is set to [Use].

Up to 253 characters can be entered.

[Port No.]

Enter the port number.

This item is available when [C-Store] is set to [Use].

You can enter a value ranging from 1 to 65535 (up to 5 digits).

The default setting is [104].

[AE Title]

Enter the AE title.

This item is available when [C-Store] is set to [Use].

Up to 16 characters can be entered.

[Retry]

Select whether to retry operations, such as confirming DICOM transmission and confirming communication with the C-Store server.

The default setting is [Use].

[Advanced] appears when [Use] is selected. For details on the settings, see “To configure the retry interval and number of retries” (page 79).

To confirm communication with the C-Store server

Confirm whether communication with the C-Store server is enabled.

- 1 After configuring each setting in the [C-Store] tab, select [Confirm] for [C-Store] and press the ENTER button.

Communication confirmation starts. An indicator appears during confirmation.
The results of confirmation appear when the process is complete.

- 2 Select [OK], and press the ENTER button.

Tip

If proper communication with the C-Store server was not established, refer to the cause displayed on the screen and resolve the problem.

[Local Station 1] Tab

[AE Title], [Station Name]

Up to 16 characters can be entered.

[Facility Name]

Up to 64 characters can be entered.

[Facility Address]

Up to 128 characters can be entered.

[Local Station 2] Tab

Configure local station settings.

[SOP Class UID]

Select [SC] or [VL] for the SOP Class UID.

The default setting is [SC].

[Character Set]

Select the character set.

You can select from [IR6], [IR13], [IR100], and [IR101].

The default setting is [IR6].

Overview

You can connect a commercially available touch panel monitor or mouse to the unit, and use them to operate the menu screens.

This chapter describes the settings for using a touch panel monitor or mouse, and the screens that appear during touch panel or mouse operations.

Note

If you set one of [USB 3] to [USB 6] to [Touch Panel / Mouse] in the [System Admin Settings] screen – [Device Settings] screen – [Device 1] tab in the touch panel / mouse mode while a compatible device is not connected, operation will become impossible. Should you inadvertently disable operation, you can restore the standard operation mode by performing the following.

- ① Connect a hard keyboard.
- ② While holding the o key (lower case), press the ⏻ (on/standby) switch, and release the o key when the unit turns off.
- ③ Press the ⏻ (on/standby) switch to exit standby mode.

Tips

- The minimum requirement for the touch panel monitor resolution is 1366 × 768.
- Be aware of the following when setting one of [USB 3] to [USB 6] to [Touch Panel / Mouse] in the [System Admin Settings] screen – [Device Settings] screen – [Device 1] tab.
 - You must restart the unit after applying the setting.
 - After applying the setting, the ▲, ▼, ◀, ▶, ENTER, and MENU buttons on the front panel of the unit and on the infrared remote control unit are disabled.
- For details on supported touch panel monitors, contact your local Sony representative.

Configuring Touch Panel / Mouse Settings

Configure the following settings to enable use of a touch panel monitor or mouse.

- 1 Select [Device Settings] in the [System Admin Settings] screen, and press the ENTER button.

For details on displaying the [System Admin Settings] screen, see “Displaying the [System Admin Settings] Screen” (page 61).

The [Device] screen appears.

- 2 Display the [Device 1] tab, and set the USB port to which the touch panel monitor or mouse will be connected ([USB 3] to [USB 6]) to [Touch Panel / Mouse].

- 3 Select [Advanced], and press the ENTER button.

The [Touch Panel / Mouse Settings] screen appears.

- 4 Set [Display Cursor] to [Use], select [Apply], and then press the ENTER button.

When [Use] is selected, a mouse cursor will appear on the screen while using touch panel mode or a mouse.

The default setting is [Not Used].

When you change this setting, a message indicating that you must restart the unit appears.

- 5 Connect the touch panel monitor or mouse, and follow the instructions on the screen and restart the unit.

Screen Displays in Touch Panel / Mouse Mode

The screens that appear in touch panel / mouse mode differ slightly from those that appear in standard mode. This section describes the main differences in the screens and buttons that appear.

[MENU] Screen



- 1 [Patient Info.]**
Displays the patient ID, patient name, gender, and date of birth.
- 2 [New Folder]**
Resets the patient information.
[Close Folder]
Use this to end recording, and display the [End Recording] screen.
- 3 [Recording Preparation] tab**
Enter patient information and configure the storage destination for recorded data before starting to record.
For details on each setting, see “Recording Preparation” (page 46).

Tip

When entering the date of birth in the touch panel monitor screen, the year is entered via the on-screen keyboard.

- 4 [Recording] tab**
Allows you to perform recording and playback controls for Ch2.

For details on each setting, see “Recording Preparation” (page 46).

Tips

- Ch1 playback cannot be performed during Ch2 recording or playback.
- Ch2 recording or playback cannot be performed during Ch1 playback.
- The items in the [Recording] tab will be grayed out during [Ch1 2K Rec] mode.

- 5 [Change Settings] tab**
Allows you to change the image quality and video input.
For details on each setting, see “Configuring Input Signal Settings” (page 34), “Configuring Image Quality Settings” (page 34), and “Configuring Recording Settings” (page 35).

- 6 [Link On] (Ch1/Ch2 link)**
Use this to enable Ch1/Ch2 simultaneous recording. This button is disabled when recording for Ch1 or Ch2 is in progress.
On (green): Recording and playback controls for Ch1 and Ch2 are linked.
Off (gray): Ch1 and Ch2 are not linked.

Tips

- When [CMS Mode] is set to [Use] in the [System Admin Settings] screen – [Function Settings] screen – [CMS] tab and [Ch1/Ch2 Rec Mode] is set to [Ch1/Ch2 Simul Rec] in the [User Settings] screen – [Recording] tab, this will be fixed at on (green).
- This button will be grayed out and disabled during [Ch1 2K Rec] mode.

7 [AUDIO] (Audio on/off)

Use this to record audio.

On (green): Audio is recorded.

Off (gray): Audio recording is stopped.

8 Status button

Displays the “[Status] Screen” (page 83).

9 [Recording List] button

Displays the [Recording List] screen.

10 HDD capacity / error indicator

Displays the internal HDD’s remaining capacity under normal conditions.

If an error occurs, the error that appears on the front panel display also appears here.

11 ⏻ (on/standby) button

Use this to set the unit to standby mode.

A confirmation message when you tap this.

12 Control buttons

Use these to control playback and recording.

The buttons function identically to the buttons on the front panel of the unit and on the infrared remote control unit.

For details, see “Names and Functions of Parts” (page 20).

Tips

- Ch1 playback cannot be performed during Ch2 recording or playback.
 - Ch2 recording or playback cannot be performed during Ch1 playback.
- 13 Ch1 image display and playback/recording status**
Displays the input image or playback image for Ch1 and various playback/recording status indicators.
- Playback/recording counter:** Displays the playback or recording time in “00:00:00” format.
- Still capture count:** Displays the number of still images you captured in “CAP:001” format.
- Rec mode:** Displays the current recording mode.

External media recording status: Displays the recording status for each external media in “■ <media name> -h -m” format.

Printer status: Displays the print status and the remaining number of print sheets.

Indicators: The indicators light according to the status of the corresponding media or printer.

■ (gray): Operating normally.

■ (green): Data is being written or printing is in progress.

■ (orange): An error has occurred.

14 [Input]

This is enabled when [2D/3D] is set to [Manual] in the [System Admin Settings] screen – [Function Settings] screen – [Input Detection] tab. When [Auto] is selected, the input signal is automatically recognized and displayed here.


This configures the Ch1 input signal.

[Ch1 Ext Storage]

Configures the external storage destination for Ch1 recorded data.

[Status] Screen

The [Status] screen displays the status of recording to external media, the conditions of connected devices, and other status information.

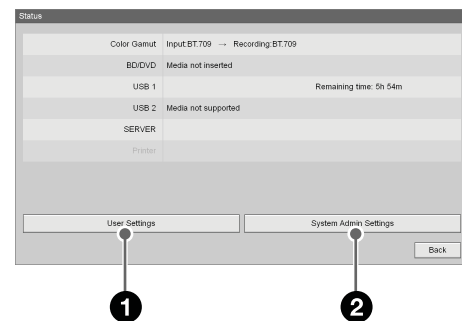

1 [User Settings]

Displays the [User Settings] screen.

For details, see “Displaying the [User Settings] Screen” (page 34).

2 [System Admin Settings]

Displays the [System Admin Settings] screen.


For details, see “Displaying the [System Admin Settings] Screen” (page 61).

Tip

[Color Gamut] information is not displayed during [Ch1 2K Rec] mode.

Video Playback Screen

This screen appears during video playback.
You can control playback using the buttons at the bottom.

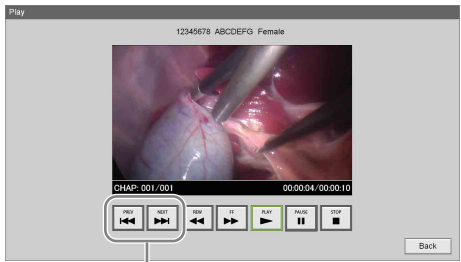


Chapter forward/back buttons

Tip

Ch1 video is played back during quick playback.

Still Image Playback Screen

This screen appears during still image playback.
You can play back the still images in sequence using the buttons at the bottom.
To play back 3D still images, tap the 3D data button.
To playback the left and right 2D still images, tap the [L] and [R] buttons.

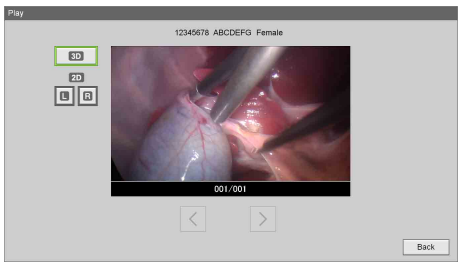

Error Messages

When contacting your local Sony representative, please indicate the error message and error ID.

Front panel display		Description	Solution
Error message	Error ID		
SYSTEM ERR.	0100 to 0199	Motherboard error	Restart the unit. (The unit will restart automatically for error ID 0600 to 0699 capture board errors.) If the error persists, contact your local Sony representative.
SYSTEM ERR.	0200 to 0299	Internal HDD error	
SYSTEM ERR.	0300 to 0399	Database error	
SYSTEM ERR.	0400 to 0499	Recovery failure	
SYSTEM ERR.	0500 to 0599	File data error	
SYSTEM ERR.	0600 to 0699	Capture board error	
SYSTEM ERR.	0700 to 0799	Optical drive error	
SYSTEM ERR.	0800 to 0899	Audio device error	
SYSTEM ERR.	0900 to 0999	Front panel error	
SYSTEM ERR.	1000 to 1099	Application error	
SYSTEM ERR.	1101	Activation error	Activation failed. Contact your local Sony representative.
SYSTEM ERR.	1102	Update error	Firmware update failed. Contact your local Sony representative.
SYSTEM ERR.	1200 to 1299	LCD controller error	An LCD controller error has occurred. The unit will restart automatically. If the error persists, contact your local Sony representative.
MAX REC DATA	—	Limit for recorded data reached	The limit for the number of recorded data entries has been reached. Delete any unnecessary recorded data.
MAX CHAPTER	—	Limit for chapters reached	The limit for the number of chapters has been reached. Stop the recording operation.
MAX STILL IMAGE	—	Limit for still images reached	Delete any unnecessary still images.
NO RECORD	—	No recorded data	Recorded data for quick playback does not exist.
HDD LACK	—	Insufficient space on internal HDD	The remaining space on the internal HDD is low. Delete any unnecessary recorded data.

Front panel display		Description	Solution
Error message	Error ID		
HDD FULL	—	Internal HDD full	There is no remaining space on the internal HDD. Delete any unnecessary recorded data.
KEY INHIBIT	—	Front panel buttons disabled	Button operations on the unit's front panel are disabled.
WARN. FAN	—	System/CPU fan warning	A fan rotation malfunction has occurred. Contact your local Sony representative.
WARN. TEMP.	—	System/CPU temperature warning	An internal temperature error has occurred. Check your operating environment.
WARN. HDD	—	Internal HDD warning or error	An internal HDD error has occurred. Contact your local Sony representative.
ERR.BATTERY	—	Battery dead	The battery is dead. Contact your local Sony representative.
WARN.HDD.2	—	Corrupted files warning	Corrupted files exist on the internal HDD. Perform recovery on the internal HDD. For details on how to perform HDD recovery, see <i>"Troubleshooting"</i> (page 87). If the error persists, contact your local Sony representative.
WARN.DB1	—	Database damaged	The database is damaged. Perform recovery on the database. For details on how to perform database recovery, see <i>"Troubleshooting"</i> (page 87). If the error persists, contact your local Sony representative.
WARN.DB2	—		
HDD RECOVERING ...	—	HDD recovery in progress	Internal HDD recovery is in progress. Wait until recovery is complete.
DB1 RECOVERING...	—	Database recovery in progress	Database recovery is in progress. Wait until recovery is complete.
DB2 RECOVERING...	—		
ERR.HDD	—	Internal HDD error	An internal HDD error has occurred. Contact your local Sony representative.
ERR.SSD	—	SSD error	Contact your local Sony representative.
WARN LCD	—	LCD backlight warning	The LCD usage time has exceeded 10,000 hours. Contact your local Sony representative.
PLAYBACK ERR.	—	Playback failed	Perform the playback operation again. If the error persists, contact your local Sony representative.

Troubleshooting

Be sure to check the following if a problem occurs. If the problem persists, contact your local Sony representative.

- The unit does not turn on when the ⏻ (on/standby) switch is pressed.
 - The main power switch on the rear panel of the unit is set to the ○ (off) position.
Set the main power switch to the I (on) position (*page 23, 24*).
 - Make sure that the power cord is properly connected to the power outlet.
- The buttons on the unit do not function.
Recording, playback, and other operations do not function properly.
 - Static electricity can have adverse effects on the functioning of the unit.
Turn off the ⏻ (on/standby) switch, disconnect the power cord, and then wait a moment before reconnecting the power cord and turning on the unit.
- The disc tray does not open when the ▲ EJECT button is pressed.
 - The disc tray may take a few seconds to open after you record data on the disc. This is because the unit is adding certain required information onto the disc. It is not a malfunction.
- The unit turns on, but images do not appear or are distorted.
 - Make sure that the cables are properly connected to the input/output connectors on the rear panel of the unit (*page 23*).
 - Make sure that the cables are not damaged.
- There is no audio, or the audio is noisy.
 - Make sure that the cables are properly connected to the input/output connectors on the rear panel of the unit (*page 23*).
 - Make sure that the cables are not damaged.
 - If the unit is in fast playback or fast reverse playback mode, press the ► PLAY button (*page 21*).
- The unit cannot record data.
 - Use media on which data can be recorded.
 - The remaining storage capacity on the media is insufficient.
Replace the media.
 - Use a BD/DVD disc that is compatible with the unit (*page 30*).
- The image freezes and audio is cut for brief (about one second) moments during playback.
 - This unit stops playback between titles, and resumes playback after about one second. It is not a malfunction.
- Data on the disc is damaged.
 - Contact your local Sony representative.
- “WARN.HDD.2” appears on the front panel display.
 - Corrupted files exist on the internal HDD. Press and hold the BACK button for at least 5 seconds to perform recovery on the internal HDD.
 - If the warning is displayed even after performing HDD recovery, the internal HDD is damaged. Contact your local Sony representative.

- “WARN.DB1” or “WARN.DB2” appears on the front panel display.
 - The database is damaged. Press and hold the BACK button for at least 5 seconds to perform database recovery.
 - If the warning is displayed even after performing database recovery, the database cannot be recovered. Contact your local Sony representative.

Note

- Do not turn off the unit while internal HDD or database recovery is in progress. The recovery process will start after the unit automatically restarts.
- Internal HDD or database recovery can take up to 1 hour each. If recovery does not complete within 1 hour, the unit will automatically shut down.

Licenses

DCMTK

Copyright (C) 1994-2011, OFFIS
e.V.

All rights reserved.

This software and supporting
documentation were developed by

OFFIS e.V.
R&D Division Health
Escherweg 2
26121 Oldenburg, Germany

Redistribution and use in source and
binary forms, with or without
modification, are permitted provided
that the following conditions
are met:

- Redistributions of source code must
retain the above copyright
notice, this list of conditions and the
following disclaimer.

- Redistributions in binary form must
reproduce the above copyright
notice, this list of conditions and the
following disclaimer in the
documentation and/or other
materials provided with the
distribution.

- Neither the name of OFFIS nor the
names of its contributors may be
used to endorse or promote products
derived from this software
without specific prior written
permission.

THIS SOFTWARE IS PROVIDED
BY THE COPYRIGHT HOLDERS
AND CONTRIBUTORS
"AS IS" AND ANY EXPRESS OR
IMPLIED WARRANTIES,
INCLUDING, BUT NOT
LIMITED TO, THE IMPLIED
WARRANTIES OF
MERCHANTABILITY AND
FITNESS FOR

A PARTICULAR PURPOSE ARE
DISCLAIMED. IN NO EVENT
SHALL THE COPYRIGHT
HOLDER OR CONTRIBUTORS
BE LIABLE FOR ANY DIRECT,
INDIRECT, INCIDENTAL,
SPECIAL, EXEMPLARY, OR
CONSEQUENTIAL DAMAGES
(INCLUDING, BUT NOT
LIMITED TO, PROCUREMENT
OF SUBSTITUTE GOODS OR
SERVICES; LOSS OF USE,
DATA, OR PROFITS; OR
BUSINESS INTERRUPTION)
HOWEVER CAUSED AND ON
ANY
THEORY OF LIABILITY,
WHETHER IN CONTRACT,
STRICT LIABILITY, OR TORT
(INCLUDING NEGLIGENCE OR
OTHERWISE) ARISING IN ANY
WAY OUT OF THE USE
OF THIS SOFTWARE, EVEN IF
ADVISED OF THE POSSIBILITY
OF SUCH DAMAGE.

Copyright (C) 1993/1994, OFFIS,
Oldenburg University and CERIU

This software and supporting
documentation were
developed by

Institut OFFIS
Bereich Kommunikationssysteme
Westerstr. 10-12
26121 Oldenburg, Germany

Fachbereich Informatik
Abteilung Prozessinformatik
Carl von Ossietzky Universitaet
Oldenburg
Ammerlaender Heerstr. 114-118
26111 Oldenburg, Germany

CERIU
Laboratoire SIM
Faculte de Medecine
2 Avenue du Pr. Leon Bernard
35043 Rennes Cedex, France

for CEN/TC251/WG4 as a
contribution to the Radiological
Society of North America (RSNA)
1993 Digital Imaging and
Communications in Medicine
(DICOM) Demonstration.

THIS SOFTWARE IS MADE
AVAILABLE, AS IS, AND
NEITHER OFFIS,
OLDENBURG UNIVERSITY NOR
CERIUM MAKE ANY
WARRANTY REGARDING
THE SOFTWARE, ITS
PERFORMANCE, ITS
MERCHANTABILITY OR
FITNESS FOR ANY
PARTICULAR USE, FREEDOM
FROM ANY COMPUTER
DISEASES OR ITS
CONFORMITY TO ANY
SPECIFICATION. THE
ENTIRE RISK AS TO QUALITY
AND PERFORMANCE OF THE
SOFTWARE
IS WITH THE USER.

Copyright of the software and
supporting documentation
is, unless otherwise stated, jointly
owned by OFFIS,
Oldenburg University and CERIU
and free access is hereby
granted as a license to use this
software, copy this
software and prepare derivative
works based upon this
software. However, any distribution
of this software
source code or supporting
documentation or derivative
works (source code and supporting
documentation) must
include the three paragraphs of this
copyright notice.

Copyright (C) 1993, RSNA and
Washington University

The software and supporting
documentation for the Radiological
Society of North America (RSNA)
1993 Digital Imaging and
Communications in Medicine
(DICOM) Demonstration were
developed
at the

Electronic Radiology
Laboratory
Mallinckrodt Institute of
Radiology
Washington University School
of Medicine

510 S. Kingshighway Blvd.
St. Louis, MO 63110
as part of the 1993 DICOM Central
Test Node project for, and
under contract with, the Radiological
Society of North America.

THIS SOFTWARE IS MADE
AVAILABLE, AS IS, AND
NEITHER RSNA NOR
WASHINGTON UNIVERSITY
MAKE ANY WARRANTY
ABOUT THE SOFTWARE, ITS
PERFORMANCE, ITS
MERCHANTABILITY OR
FITNESS FOR ANY
PARTICULAR
USE, FREEDOM FROM ANY
COMPUTER DISEASES OR ITS
CONFORMITY TO ANY
SPECIFICATION. THE ENTIRE
RISK AS TO QUALITY AND
PERFORMANCE OF
THE SOFTWARE IS WITH THE
USER.

Copyright of the software and
supporting documentation is
jointly owned by RSNA and
Washington University, and free
access
is hereby granted as a license to use
this software, copy this
software and prepare derivative
works based upon this software.
However, any distribution of this
software source code or
supporting documentation or
derivative works (source code and
supporting documentation) must
include the three paragraphs of
the copyright notice.

Copyright (c) 1998 Todd C. Miller
<Todd.Miller(at)courtesan.com>
All rights reserved.

Redistribution and use in source and
binary forms, with or without
modification, are permitted provided
that the following conditions
are met:

1. Redistributions of source code
must retain the above copyright
notice, this list of conditions and
the following disclaimer.

2. Redistributions in binary form
must reproduce the above copyright
notice, this list of conditions and
the following disclaimer in the
documentation and/or other
materials provided with the
distribution.
3. The name of the author may not be
used to endorse or promote products
derived from this software without
specific prior written permission.

THIS SOFTWARE IS PROVIDED
"AS IS" AND ANY EXPRESS OR
IMPLIED WARRANTIES,
INCLUDING, BUT NOT LIMITED
TO, THE IMPLIED
WARRANTIES OF
MERCHANTABILITY
AND FITNESS FOR A
PARTICULAR PURPOSE ARE
DISCLAIMED. IN NO EVENT
SHALL
THE AUTHOR BE LIABLE FOR
ANY DIRECT, INDIRECT,
INCIDENTAL, SPECIAL,
EXEMPLARY, OR
CONSEQUENTIAL DAMAGES
(INCLUDING, BUT NOT
LIMITED TO,
PROCUREMENT OF
SUBSTITUTE GOODS OR
SERVICES; LOSS OF USE, DATA,
OR PROFITS;
OR BUSINESS INTERRUPTION)
HOWEVER CAUSED AND ON
ANY THEORY OF LIABILITY,
WHETHER IN CONTRACT,
STRICT LIABILITY, OR TORT
(INCLUDING NEGLIGENCE OR
OTHERWISE) ARISING IN ANY
WAY OUT OF THE USE OF THIS
SOFTWARE, EVEN IF
ADVISED OF THE POSSIBILITY
OF SUCH DAMAGE.

Copyright 1988 Regents of the
University of California
Permission to use, copy, modify, and
distribute this software and
its documentation for any purpose
and without fee is hereby granted,
provided that the above copyright
notice appear in all copies. The
University of California makes no
representations about the

suitability of this software for any
purpose. It is provided "as
is" without express or implied
warranty.

The code for OFStandard::ftoa has
been derived
from an implementation which
carries the following copyright
notice:

Copyright (c) 1988 Regents of
the University of California.
All rights reserved.

Redistribution and use in source
and binary forms are permitted
provided that the above copyright
notice and this paragraph are
duplicated in all such forms and
that any documentation,
advertising materials, and other
materials related to such
distribution and use acknowledge
that the software was developed
by the University of California,
Berkeley. The name of the
University may not be used to
endorse or promote products
derived
from this software without
specific prior written permission.
THIS SOFTWARE IS
PROVIDED "AS IS" AND
WITHOUT ANY EXPRESS OR
IMPLIED WARRANTIES,
INCLUDING, WITHOUT
LIMITATION, THE IMPLIED
WARRANTIES OF
MERCHANTABILITY AND
FITNESS FOR A
PARTICULAR PURPOSE.

Copyright (C) 1999-2009
Contributors to log4cplus project.
All rights reserved.

Redistribution and use in source
and binary forms, with or without
modifica-
tion, are permitted provided that
the following conditions are met:

1. Redistributions of source code
must retain the above copyright
notice,

this list of conditions and the following disclaimer.

2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

THIS SOFTWARE IS PROVIDED ``AS IS" AND ANY EXPRESSED OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE APACHE SOFTWARE FOUNDATION OR ITS CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

Apache License
Version 2.0, January 2004
<http://www.apache.org/licenses/>

TERMS AND CONDITIONS FOR USE, REPRODUCTION, AND DISTRIBUTION

1. Definitions.

"License" shall mean the terms and conditions for use, reproduction, and distribution as defined by Sections 1 through 9 of this document.

"Licensor" shall mean the copyright owner or entity authorized by the copyright owner that is granting the License.

"Legal Entity" shall mean the union of the acting entity and all other entities that control, are controlled by, or are under common control with that entity. For the purposes of this definition, "control" means (i) the power, direct or indirect, to cause the direction or management of such entity, whether by contract or otherwise, or (ii) ownership of fifty percent (50%) or more of the outstanding shares, or (iii) beneficial ownership of such entity.

"You" (or "Your") shall mean an individual or Legal Entity exercising permissions granted by this License.

"Source" form shall mean the preferred form for making modifications, including but not limited to software source code, documentation source, and configuration files.

"Object" form shall mean any form resulting from mechanical transformation or translation of a Source form, including but not limited to compiled object code, generated documentation,

and conversions to other media types.

"Work" shall mean the work of authorship, whether in Source or Object form, made available under the License, as indicated by a copyright notice that is included in or attached to the work (an example is provided in the Appendix below).

"Derivative Works" shall mean any work, whether in Source or Object form, that is based on (or derived from) the Work and for which the editorial revisions, annotations, elaborations, or other modifications represent, as a whole, an original work of authorship. For the purposes of this License, Derivative Works shall not include works that remain separable from, or merely link (or bind by name) to the interfaces of, the Work and Derivative Works thereof.

"Contribution" shall mean any work of authorship, including the original version of the Work and any modifications or additions to that Work or Derivative Works thereof, that is intentionally submitted to Licensor for inclusion in the Work by the copyright owner or by an individual or Legal Entity authorized to submit on behalf of the copyright owner. For the purposes of this definition, "submitted" means any form of electronic, verbal, or written communication sent to the Licensor or its representatives, including but not limited to

communication on electronic mailing lists, source code control systems,

and issue tracking systems that are managed by, or on behalf of, the

Licensors for the purpose of discussing and improving the Work, but

excluding communication that is conspicuously marked or otherwise

designated in writing by the copyright owner as "Not a Contribution."

"Contributor" shall mean Licensor and any individual or Legal Entity

on behalf of whom a Contribution has been received by Licensor and

subsequently incorporated within the Work.

2. Grant of Copyright License. Subject to the terms and conditions of

this License, each Contributor hereby grants to You a perpetual, worldwide, non-exclusive, no-charge, royalty-free, irrevocable copyright license to reproduce, prepare Derivative Works of, publicly display, publicly perform, sublicense, and distribute the

Work and such Derivative Works in Source or Object form.

3. Grant of Patent License. Subject to the terms and conditions of

this License, each Contributor hereby grants to You a perpetual, worldwide, non-exclusive, no-charge, royalty-free, irrevocable (except as stated in this section) patent license to make, have made,

use, offer to sell, sell, import, and otherwise transfer the Work, where such license applies only to those patent claims licensable by such Contributor that are necessarily infringed by their

Contribution(s) alone or by combination of their Contribution(s)

with the Work to which such Contribution(s) was submitted. If You

institute patent litigation against any entity (including a cross-claim or counterclaim in a lawsuit) alleging that the Work or a Contribution incorporated within the Work constitutes direct

or contributory patent infringement, then any patent licenses

granted to You under this License for that Work shall terminate

as of the date such litigation is filed.

4. Redistribution. You may reproduce and distribute copies of the

Work or Derivative Works thereof in any medium, with or without modifications, and in Source or Object form, provided that You meet the following conditions:

(a) You must give any other recipients of the Work or Derivative Works a copy of this License; and

(b) You must cause any modified files to carry prominent notices stating that You changed the files; and

(c) You must retain, in the Source form of any Derivative Works that You distribute, all copyright, patent, trademark, and attribution notices from the Source form of the Work, excluding those notices that do not pertain to any part of the Derivative Works; and

(d) If the Work includes a "NOTICE" text file as part of its

distribution, then any Derivative Works that You distribute must

include a readable copy of the attribution notices contained within such NOTICE file, excluding those notices that do not

pertain to any part of the Derivative Works, in at least one of the following places: within a NOTICE text file distributed

as part of the Derivative Works; within the Source form or documentation, if provided along with the Derivative Works; or,

within a display generated by the Derivative Works, if and wherever such third-party notices normally appear. The contents

of the NOTICE file are for informational purposes only and do not modify the License. You may add Your own attribution

notices within Derivative Works that You distribute, alongside

or as an addendum to the NOTICE text from the Work, provided

that such additional attribution notices cannot be construed as modifying the License.

You may add Your own copyright statement to Your modifications and

may provide additional or different license terms and conditions

for use, reproduction, or distribution of Your modifications, or

for any such Derivative Works as a whole, provided Your use, reproduction, and distribution of the Work otherwise complies with

the conditions stated in this License.

5. Submission of Contributions.
Unless You explicitly state otherwise,

any Contribution intentionally submitted for inclusion in the Work

by You to the Licensor shall be under the terms and conditions of this License, without any additional terms or conditions.

Notwithstanding the above, nothing herein shall supersede or modify

the terms of any separate license agreement you may have executed

with Licensor regarding such Contributions.

6. Trademarks. This License does not grant permission to use the trade

names, trademarks, service marks, or product names of the Licensor,

except as required for reasonable and customary use in describing the

origin of the Work and reproducing the content of the NOTICE file.

7. Disclaimer of Warranty.
Unless required by applicable law or

agreed to in writing, Licensor provides the Work (and each Contributor provides its Contributions) on an "AS IS" BASIS,

WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied, including, without limitation, any warranties or conditions

of TITLE, NON-INFRINGEMENT, MERCHANTABILITY, or FITNESS FOR A

PARTICULAR PURPOSE. You are solely responsible for determining the

appropriateness of using or redistributing the Work and assume any

risks associated with Your exercise of permissions under this License.

8. Limitation of Liability. In no event and under no legal theory, whether in tort (including negligence), contract, or otherwise,

unless required by applicable law (such as deliberate and grossly

negligent acts) or agreed to in writing, shall any Contributor be liable to You for damages, including any direct, indirect, special,

incidental, or consequential damages of any character arising as a

result of this License or out of the use or inability to use the Work (including but not limited to damages for loss of goodwill, work stoppage, computer failure or malfunction, or any and all

other commercial damages or losses), even if such Contributor has been advised of the possibility of such damages.

9. Accepting Warranty or Additional Liability. While redistributing

the Work or Derivative Works thereof, You may choose to offer, and charge a fee for, acceptance of support, warranty, indemnity, or other liability obligations and/or rights consistent with this

License. However, in accepting such obligations, You may act only

on Your own behalf and on Your sole responsibility, not on behalf

of any other Contributor, and only if You agree to indemnify, defend, and hold each Contributor harmless for any liability

incurred by, or claims asserted against, such Contributor by reason

of your accepting any such warranty or additional liability.

END OF TERMS AND CONDITIONS

APPENDIX: How to apply the Apache License to your work.

To apply the Apache License to your work, attach the following boilerplate notice, with the fields enclosed by brackets "[]" replaced with your own identifying information. (Don't include

the brackets!) The text should be enclosed in the appropriate comment syntax for the file format. We also recommend that a

file or class name and description of purpose be included on the

same "printed page" as the copyright notice for easier identification within third-party archives.

Copyright 1999-2009
Contributors to log4cplus project

Licensed under the Apache License, Version 2.0 (the "License");

you may not use this file except in compliance with the License.

You may obtain a copy of the License at

<http://www.apache.org/licenses/LICENSE-2.0>

Unless required by applicable law or agreed to in writing, software

distributed under the License is distributed on an "AS IS" BASIS,

WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied.

See the License for the specific language governing permissions and

limitations under the License.

zlib

(C) 1995-2004 Jean-loup Gailly and Mark Adler

This software is provided 'as-is', without any express or implied warranty. In no event will the authors be held liable for any damages arising from the use of this software.

Permission is granted to anyone to use this software for any purpose, including commercial applications, and to alter it and redistribute it freely, subject to the following restrictions:

1. The origin of this software must not be misrepresented; you must not claim that you wrote the original software. If you use this software in a product, an acknowledgment in the product documentation would be appreciated but is not required.
2. Altered source versions must be plainly marked as such, and must not be misrepresented as being the original software.
3. This notice may not be removed or altered from any source distribution.

Jean-loup Gailly Mark Adler
jloup@gzip.org
madler@alumni.caltech.edu

If you use the zlib library in a product, we would appreciate *not* receiving lengthy legal documents to sign. The sources are provided for free but without warranty of any kind. The library has been entirely written by Jean-loup Gailly and Mark Adler; it does not include third-party code.

If you redistribute modified sources, we would appreciate that you include in the file ChangeLog history information documenting your changes. Please

read the FAQ for more information on the distribution of modified source versions.

PostgreSQL

PostgreSQL is released under the PostgreSQL License, a liberal Open Source license, similar to the BSD or MIT licenses.

PostgreSQL Database Management System
(formerly known as Postgres, then as Postgres95)

Portions Copyright (c) 1996-2013,
The PostgreSQL Global
Development Group

Portions Copyright (c) 1994, The
Regents of the University of
California

Permission to use, copy, modify, and distribute this software and its documentation for any purpose, without fee, and without a written agreement is hereby granted, provided that the above copyright notice and this paragraph and the following two paragraphs appear in all copies.

IN NO EVENT SHALL THE UNIVERSITY OF CALIFORNIA BE LIABLE TO ANY PARTY FOR DIRECT, INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, INCLUDING LOST PROFITS, ARISING OUT OF THE USE OF THIS SOFTWARE AND ITS DOCUMENTATION, EVEN IF THE UNIVERSITY OF CALIFORNIA HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

THE UNIVERSITY OF CALIFORNIA SPECIFICALLY DISCLAIMS ANY WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE

IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE SOFTWARE PROVIDED HEREUNDER IS ON AN "AS IS" BASIS, AND THE UNIVERSITY OF CALIFORNIA HAS NO OBLIGATIONS TO PROVIDE MAINTENANCE, SUPPORT, UPDATES, ENHANCEMENTS, OR MODIFICATIONS.

expat

Copyright (c) 1998, 1999, 2000 Thai Open Source Software Center Ltd and Clark Cooper

Copyright (c) 2001, 2002, 2003, 2004, 2005, 2006 Expat maintainers.

Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions:

The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT,

TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

freetype

The FreeType Project LICENSE

2006-Jan-27

Copyright 1996-2002, 2006 by
David Turner, Robert Wilhelm, and
Werner Lemberg

Introduction

=====

The FreeType Project is distributed in several archive packages; some of them may contain, in addition to the FreeType font engine, various tools and contributions which rely on, or relate to, the FreeType Project.

This license applies to all files found in such packages, and which do not fall under their own explicit license. The license affects thus the FreeType font engine, the test programs, documentation and makefiles, at the very least.

This license was inspired by the BSD, Artistic, and IJG (Independent JPEG Group) licenses, which all encourage inclusion

and use of free software in commercial and freeware products alike. As a consequence, its main points are that:

- o We don't promise that this software works. However, we will be interested in any kind of bug reports. ('as is' distribution)

- o You can use this software for whatever you want, in parts or

full form, without having to pay us. ('royalty-free' usage)

- o You may not pretend that you wrote this software. If you use it, or only parts of it, in a program, you must acknowledge somewhere in your documentation that you have used the

FreeType code. ('credits')

We specifically permit and encourage the inclusion of this software, with or without modifications, in commercial products.

We disclaim all warranties covering The FreeType Project and assume no liability related to The FreeType Project.

Finally, many people asked us for a preferred form for a credit/disclaimer to use in compliance with this license. We thus encourage you to use the following text:

""""

Portions of this software are copyright (c) <year> The FreeType Project (www.freetype.org). All rights reserved.

""""

Please replace <year> with the value from the FreeType version you actually use.

Legal Terms

=====

0. Definitions

Throughout this license, the terms 'package', 'FreeType Project', and 'FreeType archive' refer to the set of files originally distributed by the authors (David Turner, Robert Wilhelm, and Werner Lemberg) as the 'FreeType Project', be they named as alpha, beta or final release.

'You' refers to the licensee, or person using the project, where 'using' is a generic term including compiling the project's source code as well as linking it to form a 'program' or 'executable'.
This program is referred to as 'a program using the FreeType engine'.

This license applies to all files distributed in the original FreeType Project, including all source code, binaries and documentation, unless otherwise stated in the file in its original, unmodified form as distributed in the original archive.

If you are unsure whether or not a particular file is covered by this license, you must contact us to verify this.

The FreeType Project is copyright (C) 1996-2000 by David Turner, Robert Wilhelm, and Werner Lemberg. All rights reserved except as specified below.

1. No Warranty

THE FREETYPE PROJECT IS PROVIDED 'AS IS' WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. IN NO EVENT WILL ANY OF THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY DAMAGES CAUSED BY THE USE OR THE INABILITY TO USE, OF THE FREETYPE PROJECT.

2. Redistribution

This license grants a worldwide, royalty-free, perpetual and

irrevocable right and license to use, execute, perform, compile, display, copy, create derivative works of, distribute and sublicense the FreeType Project (in both source and object code forms) and derivative works thereof for any purpose; and to authorize others to exercise some or all of the rights granted herein, subject to the following conditions:

- o Redistribution of source code must retain this license file ('FTL.TXT') unaltered; any additions, deletions or changes to the original files must be clearly indicated in accompanying documentation. The copyright notices of the unaltered, original files must be preserved in all copies of source files.

- o Redistribution in binary form must provide a disclaimer that states that the software is based in part of the work of the FreeType Team, in the distribution documentation. We also encourage you to put an URL to the FreeType web page in your documentation, though this isn't mandatory.

These conditions apply to any software derived from or based on the FreeType Project, not just the unmodified files. If you use our work, you must acknowledge us. However, no fee need be paid to us.

3. Advertising -----

Neither the FreeType authors and contributors nor you shall use the name of the other for commercial, advertising, or promotional purposes without specific prior written permission.

We suggest, but do not require, that you use one or more of the

following phrases to refer to this software in your documentation or advertising materials: 'FreeType Project', 'FreeType Engine', 'FreeType library', or 'FreeType Distribution'.

As you have not signed this license, you are not required to accept it. However, as the FreeType Project is copyrighted material, only this license, or another one contracted with the authors, grants you the right to use, distribute, and modify it. Therefore, by using, distributing, or modifying the FreeType Project, you indicate that you understand and accept all the terms of this license.

4. Contacts -----

There are two mailing lists related to FreeType:

- o freetype@nongnu.org

Discusses general use and applications of FreeType, as well as future and wanted additions to the library and distribution.

If you are looking for support, start in this list if you haven't found anything to help you in the documentation.

- o freetype-devel@nongnu.org

Discusses bugs, as well as engine internals, design issues, specific licenses, porting, etc.

Our home page can be found at

<http://www.freetype.org>

--- end of FTL.TXT ---

memcached

Copyright (c) 2003, Danga Interactive, Inc.
All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.

Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

Neither the name of the Danga Interactive nor the names of its contributors may be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT OWNER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT

(INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

openldap

The OpenLDAP Public License
Version 2.8, 17 August 2003

Redistribution and use of this software and associated documentation ("Software"), with or without modification, are permitted provided that the following conditions are met:

1. Redistributions in source form must retain copyright statements and notices,
2. Redistributions in binary form must reproduce applicable copyright statements and notices, this list of conditions, and the following disclaimer in the documentation and/or other materials provided with the distribution, and
3. Redistributions must contain a verbatim copy of this document.

The OpenLDAP Foundation may revise this license from time to time. Each revision is distinguished by a version number. You may use this Software under terms of this license revision or under the terms of any subsequent revision of the license.

THIS SOFTWARE IS PROVIDED BY THE OPENLDAP FOUNDATION AND ITS CONTRIBUTORS "AS IS" AND ANY EXPRESSED OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY

AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE OPENLDAP FOUNDATION, ITS CONTRIBUTORS, OR THE AUTHOR(S) OR OWNER(S) OF THE SOFTWARE BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

The names of the authors and copyright holders must not be used in advertising or otherwise to promote the sale, use or other dealing in this Software without specific, written prior permission. Title to copyright in this Software shall at all times remain with copyright holders.

OpenLDAP is a registered trademark of the OpenLDAP Foundation.

Copyright 1999-2003 The OpenLDAP Foundation, Redwood City, California, USA. All Rights Reserved. Permission to copy and distribute verbatim copies of this document is granted.

OpenSSL

The OpenSSL toolkit stays under a dual license, i.e. both the conditions of the OpenSSL License and the original SSLeay license apply to the toolkit.

See below for the actual license texts. Actually both licenses are BSD-style Open Source licenses. In case of any license issues related to OpenSSL please contact openssl-core@openssl.org.

OpenSSL License

Copyright (c) 1998-2011 The OpenSSL Project. All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. All advertising materials mentioning features or use of this software must display the following acknowledgment:
"This product includes software developed by the OpenSSL Project for use in the OpenSSL Toolkit. (<http://www.openssl.org/>)"
4. The names "OpenSSL Toolkit" and "OpenSSL Project" must not be used to endorse or promote products derived from this software without prior written permission. For written permission, please contact openssl-core@openssl.org.

5. Products derived from this software may not be called "OpenSSL" nor may "OpenSSL" appear in their names without prior written permission of the OpenSSL Project.

6. Redistributions of any form whatsoever must retain the following acknowledgment:

"This product includes software developed by the OpenSSL Project for use in the OpenSSL Toolkit (<http://www.openssl.org/>)"

THIS SOFTWARE IS PROVIDED BY THE OpenSSL PROJECT "AS IS" AND ANY EXPRESSED OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE OpenSSL PROJECT OR ITS CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

Original SSLeay License

Copyright (C) 1995-1998 Eric Young (eay@cryptsoft.com)
All rights reserved.

This package is an SSL implementation written by Eric Young (eay@cryptsoft.com). The implementation was written so as to conform with Netscapes SSL.

This library is free for commercial and non-commercial use as long as the following conditions are adhered to. The following conditions apply to all code found in this distribution, be it the RC4, RSA, lhash, DES, etc., code; not just the SSL code. The SSL documentation included with this distribution is covered by the same copyright terms except that the holder is Tim Hudson (tjh@cryptsoft.com).

Copyright remains Eric Young's, and as such any Copyright notices in the code are not to be removed. If this package is used in a product, Eric Young should be given attribution as the author of the parts of the library used.

This can be in the form of a textual message at program startup or in documentation (online or textual) provided with the package.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. All advertising materials mentioning features or use of this software

must display the following acknowledgement:

"This product includes cryptographic software written by Eric Young (eay@cryptsoft.com)"
The word 'cryptographic' can be left out if the routines from the library being used are not cryptographic related :-).

4. If you include any Windows specific code (or a derivative thereof) from

the apps directory (application code) you must include an acknowledgement:

"This product includes software written by Tim Hudson (tjh@cryptsoft.com)"

THIS SOFTWARE IS PROVIDED BY ERIC YOUNG "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE AUTHOR OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

The licence and distribution terms for any publically available version or derivative of this code cannot be changed. i.e. this code cannot simply be copied and put under another distribution licence [including the GNU Public Licence.]

proj

Copyright (c) 2000, Frank Warmerdam

Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions:

The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT, TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE

SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

zlib

Copyright (C) 1995-2012 Jean-loup Gailly and Mark Adler

This software is provided 'as-is', without any express or implied warranty. In no event will the authors be held liable for any damages arising from the use of this software.

Permission is granted to anyone to use this software for any purpose, including commercial applications, and to alter it and redistribute it freely, subject to the following restrictions:

1. The origin of this software must not be misrepresented; you must not claim that you wrote the original software. If you use this software in a product, an acknowledgment in the product documentation would be appreciated but is not required.
2. Altered source versions must be plainly marked as such, and must not be misrepresented as being the original software.
3. This notice may not be removed or altered from any source distribution.

Jean-loup Gailly Mark Adler
jloup@gzip.org
madler@alumni.caltech.edu

gdal, libxm2, libxslt

Copyright (C) 1996 X Consortium

Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including

without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions:

The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL THE X CONSORTIUM BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT, TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

Except as contained in this notice, the name of the X Consortium shall not be used in advertising or otherwise to promote the sale, use or other dealings in this Software without prior written authorization from the X Consortium.

X Window System is a trademark of X Consortium, Inc.

Retrieved from "<http://directory.fsf.org/wiki/License:X11>"

Permission is granted to copy, distribute and/or modify this document under the terms of the GNU Free Documentation License, Version 1.3 or any later version published by the Free Software Foundation; with no Invariant

Sections, no Front-Cover Texts, and no Back-Cover Texts. A copy of the license is included in the page GNU Free Documentation License.

The copyright and license notices on this page only apply to the text on this page. Any software described in this text has its own copyright notice and license, which can usually be found in the distribution itself.

libedit, libevent, libmemcached, tiff

Libevent is available for use under the following license, commonly known as the 3-clause (or "modified") BSD license:

Copyright (c) 2000-2007 Niels Provos <provos@citi.umich.edu>
Copyright (c) 2007-2010 Niels Provos and Nick Mathewson

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. The name of the author may not be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE AUTHOR "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES

OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE AUTHOR BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

Portions of Libevent are based on works by others, also made available by them under the three-clause BSD license above. The copyright notices are available in the corresponding source files; the license is as above. Here's a list:

log.c:

Copyright (c) 2000 Dug Song
<dugsong@monkey.org>

Copyright (c) 1993 The Regents of the University of California.

strlcpy.c:

Copyright (c) 1998 Todd C. Miller
<Todd.Miller@courtesan.com>

win32select.c:

Copyright (c) 2003 Michael A. Davis <mike@datanerds.net>

evport.c:

Copyright (c) 2007 Sun Microsystems

ht-internal.h:

Copyright (c) 2002 Christopher Clark

minheap-internal.h:

Copyright (c) 2006 Maxim Yegorushkin
<maxim.yegorushkin@gmail.com>

The arc4module is available under the following, sometimes called the "OpenBSD" license:

Copyright (c) 1996, David Mazieres <dm@uun.org>
Copyright (c) 2008, Damien Miller <djm@openbsd.org>

Permission to use, copy, modify, and distribute this software for any purpose with or without fee is hereby granted, provided that the above

copyright notice and this permission notice appear in all copies.

THE SOFTWARE IS PROVIDED "AS IS" AND THE AUTHOR DISCLAIMS ALL WARRANTIES WITH REGARD TO THIS SOFTWARE INCLUDING ALL IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS. IN NO EVENT SHALL THE AUTHOR BE LIABLE FOR ANY SPECIAL, DIRECT, INDIRECT, OR CONSEQUENTIAL DAMAGES OR ANY DAMAGES WHATSOEVER RESULTING FROM LOSS OF USE, DATA OR PROFITS, WHETHER IN AN ACTION OF CONTRACT, NEGLIGENCE OR OTHER TORTIOUS ACTION, ARISING OUT OF OR IN CONNECTION WITH THE USE OR PERFORMANCE OF THIS SOFTWARE.

wxwidget

wxWidgets is currently licensed under the "wxWindows Licence" pending approval of the "wxWidgets Licence" which will be identical apart from the name.

The wxWindows Licence is essentially the L-GPL (Library General Public Licence), with an exception stating that derived works in binary form may be distributed on the user's own terms. This is a solution that satisfies those who wish to produce GPL'ed software using wxWidgets, and also those producing proprietary software.

wxWindows Library Licence,
Version 3.1

=====
Copyright (C) 1998-2005 Julian
Smart, Robert Roebing et al

Everyone is permitted to copy and
distribute verbatim copies
of this licence document, but
changing it is not allowed.

WXWINDOWS LIBRARY LICENCE TERMS AND CONDITIONS FOR COPYING, DISTRIBUTION AND MODIFICATION

This library is free software; you
can redistribute it and/or modify it
under the terms of the GNU Library
General Public Licence as published
by
the Free Software Foundation;
either version 2 of the Licence, or (at
your option) any later version.

This library is distributed in the
hope that it will be useful, but
WITHOUT ANY WARRANTY;
without even the implied warranty of
MERCHANTABILITY or
**FITNESS FOR A PARTICULAR
PURPOSE.** See the GNU Library
General Public Licence for more
details.

You should have received a copy of
the GNU Library General Public
Licence

along with this software, usually in
a file named COPYING.LIB. If not,
write to the Free Software
Foundation, Inc., 59 Temple Place,
Suite 330,
Boston, MA 02111-1307 USA.

EXCEPTION NOTICE

1. As a special exception, the
copyright holders of this library give
permission for additional uses of the
text contained in this release of
the library as licensed under the
wxWindows Library Licence,
applying
either version 3.1 of the Licence, or
(at your option) any later version of
the Licence as published by the
copyright holders of version
3.1 of the Licence document.

2. The exception is that you may
use, copy, link, modify and distribute
under your own terms, binary object
code versions of works based
on the Library.

3. If you copy code from files
distributed under the terms of the
GNU
General Public Licence or the GNU
Library General Public Licence into a
copy of this library, as this licence
permits, the exception does not
apply to the code that you add in this
way. To avoid misleading anyone as
to the status of such modified files,
you must delete this exception
notice from such code and/or adjust
the licensing conditions notice
accordingly.

4. If you write modifications of your
own for this library, it is your
choice whether to permit this
exception to apply to your
modifications.
If you do not wish that, you must
delete the exception notice from such
code and/or adjust the licensing
conditions notice accordingly.

ncurses

Copyright (c) 2006 Free Software
Foundation, Inc.

Permission is hereby granted, free of
charge, to any person obtaining a
copy of this software and associated
documentation files (the
"Software"), to deal in the Software
without restriction, including
without limitation the rights to use,
copy, modify, merge, publish,
distribute, distribute with
modifications, sublicense, and/or
sell copies of the Software, and to
permit persons to whom the
Software is furnished to do so,
subject to the following conditions:

The above copyright notice and this
permission notice shall be included
in all copies or substantial portions
of the Software.

THE SOFTWARE IS PROVIDED
"AS IS", WITHOUT WARRANTY
OF ANY KIND, EXPRESS OR
IMPLIED, INCLUDING BUT NOT
LIMITED TO THE WARRANTIES
OF MERCHANTABILITY,
FITNESS FOR A PARTICULAR
PURPOSE AND
NONINFRINGEMENT. IN NO
EVENT SHALL THE ABOVE
COPYRIGHT HOLDERS BE
LIABLE FOR ANY CLAIM,
DAMAGES OR OTHER
LIABILITY, WHETHER IN AN
ACTION OF CONTRACT, TORT
OR OTHERWISE, ARISING
FROM, OUT OF OR IN
CONNECTION WITH THE
SOFTWARE OR THE USE OR
OTHER DEALINGS IN THE
SOFTWARE.

Except as contained in this notice,
the name(s) of the above copyright
holders shall not be used in
advertising or otherwise to promote
the sale, use or other dealings in this
Software without prior written
authorization.

\$Id: AUTHORS,v 1.2 2006/10/28
21:44:52 tom Exp \$

These are the principal authors/
contributors of ncurses since 1.9.9e,
in decreasing order of their
contribution:

TD Thomas E. Dickey
JPF Juergen Pfeifer
ESR Eric S Raymond
AVL Alexander V Lukyanov
PB Philippe Blain
SV Sven Verdoolaege

krb5

Copyright (C) 1985-2012 by the
Massachusetts Institute of
Technology.

All rights reserved.

Export of software employing
encryption from the United States of
America may require a specific
license from the United States
Government. It is the
responsibility of any person or
organization
contemplating export to obtain
such a license before exporting.

WITHIN THAT CONSTRAINT,
permission to use, copy, modify, and
distribute
this software for any purpose and
without fee is hereby granted,
provided that the above copyright
notice appear in all copies and that
both that copyright notice and this
permission notice appear in
supporting documentation, and that
the name of M.I.T. not be used in
advertising or publicity pertaining to
distribution of the software
without specific, written prior
permission. Furthermore if you
modify
this software you must label your
software as modified software and
not
distribute it in such a fashion that it
might be confused with the

original MIT software. M.I.T. makes
no representations about the
suitability of this software for any
purpose. It is provided "as is"
without express or implied warranty.

Documentation components of this
software distribution are licensed
under a Creative Commons
Attribution-ShareAlike 3.0
Unported License.
([http://creativecommons.org/
licenses/by-sa/3.0/](http://creativecommons.org/licenses/by-sa/3.0/))

Individual source code files are
copyright MIT, Cygnus Support,
Novell,
OpenVision Technologies, Oracle,
Red Hat, Sun Microsystems,
FundsXpress, and others.

Project Athena, Athena, Athena
MUSE, Discuss, Hesiod, Kerberos,
Moirá,
and Zephyr are trademarks of the
Massachusetts Institute of
Technology
(MIT). No commercial use of these
trademarks may be made without
prior
written permission of MIT.

"Commercial use" means use of a
name in a product or other for-profit
manner. It does NOT prevent a
commercial firm from referring to
the
MIT trademarks in order to convey
information (although in doing so,
recognition of their trademark status
should be given).

curl

COPYRIGHT AND PERMISSION
NOTICE
Copyright (c) 1996 - 2011, Daniel
Stenberg, <daniel@haxx.se>.
All rights reserved.
Permission to use, copy, modify, and
distribute this software for any
purpose
with or without fee is hereby granted,
provided that the above copyright
notice and this permission notice
appear in all copies.

THE SOFTWARE IS PROVIDED
"AS IS", WITHOUT WARRANTY
OF ANY KIND, EXPRESS OR
IMPLIED, INCLUDING BUT NOT
LIMITED TO THE WARRANTIES
OF MERCHANTABILITY,
FITNESS FOR A PARTICULAR
PURPOSE AND
NONINFRINGEMENT OF THIRD
PARTY RIGHTS. IN
NO EVENT SHALL THE
AUTHORS OR COPYRIGHT
HOLDERS BE LIABLE FOR ANY
CLAIM,
DAMAGES OR OTHER
LIABILITY, WHETHER IN AN
ACTION OF CONTRACT, TORT
OR
OTHERWISE, ARISING FROM,
OUT OF OR IN CONNECTION
WITH THE SOFTWARE OR THE
USE
OR OTHER DEALINGS IN THE
SOFTWARE.

Except as contained in this notice,
the name of a copyright holder shall
not
be used in advertising or otherwise to
promote the sale, use or other
dealings
in this Software without prior written
authorization of the copyright
holder.

cyrus-sasl

Copyright (c) 1998-2003 Carnegie
Mellon University. All rights
reserved.

Redistribution and use in source and
binary forms, with or without
modification, are permitted provided
that the following conditions
are met:

1. Redistributions of source code
must retain the above copyright
notice, this list of conditions and
the following disclaimer.
2. Redistributions in binary form
must reproduce the above copyright

notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

3. The name "Carnegie Mellon University" must not be used to endorse or promote products derived from this software without prior written permission. For permission or any other legal details, please contact
Office of Technology Transfer
Carnegie Mellon University
5000 Forbes Avenue
Pittsburgh, PA 15213-3890
(412) 268-4387, fax: (412) 268-7395
tech-transfer@andrew.cmu.edu
4. Redistributions of any form whatsoever must retain the following acknowledgment:
"This product includes software developed by Computing Services at Carnegie Mellon University (<http://www.cmu.edu/computing/>)."

CARNEGIE MELLON UNIVERSITY DISCLAIMS ALL WARRANTIES WITH REGARD TO THIS SOFTWARE, INCLUDING ALL IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS, IN NO EVENT SHALL CARNEGIE MELLON UNIVERSITY BE LIABLE FOR ANY SPECIAL, INDIRECT OR CONSEQUENTIAL DAMAGES OR ANY DAMAGES WHATSOEVER RESULTING FROM LOSS OF USE, DATA OR PROFITS, WHETHER IN AN ACTION OF CONTRACT, NEGLIGENCE OR OTHER TORTIOUS ACTION, ARISING OUT OF OR IN CONNECTION WITH THE USE OR PERFORMANCE OF THIS SOFTWARE.

libjpeg

This software is based in part on the work of the Independent JPEG Group.

LibTIFF

Copyright (c) 1988-1997 Sam Leffler
Copyright (c) 1991-1997 Silicon Graphics, Inc.

Permission to use, copy, modify, distribute, and sell this software and its documentation for any purpose is hereby granted without fee, provided that (i) the above copyright notices and this permission notice appear in all copies of the software and related documentation, and (ii) the names of Sam Leffler and Silicon Graphics may not be used in any advertising or publicity relating to the software without the specific, prior written permission of Sam Leffler and Silicon Graphics.

THE SOFTWARE IS PROVIDED "AS-IS" AND WITHOUT WARRANTY OF ANY KIND, EXPRESS, IMPLIED OR OTHERWISE, INCLUDING WITHOUT LIMITATION, ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

IN NO EVENT SHALL SAM LEFFLER OR SILICON GRAPHICS BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES OF ANY KIND, OR ANY DAMAGES WHATSOEVER RESULTING FROM LOSS OF USE, DATA OR PROFITS, WHETHER OR NOT ADVISED OF THE POSSIBILITY OF DAMAGE, AND ON ANY THEORY OF

LIABILITY, ARISING OUT OF OR IN CONNECTION WITH THE USE OR PERFORMANCE OF THIS SOFTWARE.

GNU GPL/LGPL Licensed Software

This product includes the following software licensed under GNU General Public License (GPL) and GNU Lesser General Public License (LGPL) terms.

The attached GPL/LGPL conditions inform the user about their rights for procuring, altering, and redistributing the software source code.

For details on the source code, contact your local Sony representative. However, be aware that Sony will not answer questions concerning the data in the source code.

Package list

- geos
- gettext
- libgcrypt
- libgpg_error
- libiconv
- pthreads
- termcap
- unixODBC
- uuid
- wxwidget

GNU GENERAL PUBLIC LICENSE

Version 2, June 1991

Copyright (C) 1989, 1991 Free Software Foundation, Inc., 51 Franklin Street, Fifth Floor, Boston, MA 02110-1301 USA

Everyone is permitted to copy and distribute verbatim copies of this license document, but changing it is not allowed.

Preamble

The licenses for most software are designed to take away your freedom to share and change it. By contrast, the GNU General Public License is intended to guarantee your freedom to share and change free software--to make sure the software is free for all its users. This General Public License applies to most of the Free Software Foundation's software and to any other program whose authors commit to using it. (Some other Free Software Foundation software is covered by the GNU Lesser General Public License instead.) You can apply it to your programs, too.

When we speak of free software, we are referring to freedom, not price. Our General Public Licenses are designed to make sure that you have the freedom to distribute copies of free software (and charge for this service if you wish), that you receive source code or can get it if you want it, that you can change the software or use pieces of it in new free programs; and that you know you can do these things.

To protect your rights, we need to make restrictions that forbid anyone to deny you these rights or to ask you to surrender the rights.

These restrictions translate to certain responsibilities for you if you distribute copies of the software, or if you modify it.

For example, if you distribute copies of such a program, whether gratis or for a fee, you must give the recipients all the rights that you have. You must make sure that they, too, receive or can get the source code. And you must show them these terms so they know their rights.

We protect your rights with two steps: (1) copyright the software, and (2) offer you this license which gives you legal permission to copy, distribute and/or modify the software.

Also, for each author's protection and ours, we want to make certain that everyone understands that there is no warranty for this free software. If the software is modified by someone else and passed on, we want its recipients to know that what they have is not the original, so that any problems introduced by others will not reflect on the original authors' reputations.

Finally, any free program is threatened constantly by software patents. We wish to avoid the danger that redistributors of a free program will individually obtain patent licenses, in effect making the program proprietary. To prevent this, we have made it clear that any patent must be licensed for everyone's free use or not licensed at all.

The precise terms and conditions for copying, distribution and modification follow.

TERMS AND CONDITIONS FOR COPYING,

DISTRIBUTION AND MODIFICATION

0. This License applies to any program or other work which contains a notice placed by the copyright holder saying it may be distributed under the terms of this General Public License. The "Program", below, refers to any such program or work, and a "work based on the Program" means either the Program or any derivative work under copyright law: that is to say, a work containing the Program or a portion of it, either verbatim or with modifications and/or translated into another language. (Hereinafter, translation is included without limitation in the term "modification".) Each licensee is addressed as "you".

Activities other than copying, distribution and modification are not covered by this License; they are outside its scope. The act of running the Program is not restricted, and the output from the Program is covered only if its contents constitute a work based on the Program (independent of having been made by running the Program). Whether that is true depends on what the Program does.

1. You may copy and distribute verbatim copies of the Program's source code as you receive it, in any medium, provided that you conspicuously and appropriately publish on each copy an appropriate copyright notice and disclaimer of warranty; keep intact all the notices that refer to this License and to the absence of any warranty; and give any other recipients of the Program a copy of this License along with the Program.

You may charge a fee for the physical act of transferring a copy,

and you may at your option offer warranty protection in exchange for a fee.

2. You may modify your copy or copies of the Program or any portion of it, thus forming a work based on the Program, and copy and distribute such modifications or work under the terms of Section 1 above, provided that you also meet all of these conditions:
 - a) You must cause the modified files to carry prominent notices stating that you changed the files and the date of any change.
 - b) You must cause any work that you distribute or publish, that in whole or in part contains or is derived from the Program or any part thereof, to be licensed as a whole at no charge to all third parties under the terms of this License.
 - c) If the modified program normally reads commands interactively when run, you must cause it, when started running for such interactive use in the most ordinary way, to print or display an announcement including an appropriate copyright notice and a notice that there is no warranty (or else, saying that you provide a warranty) and that users may redistribute the program under these conditions, and telling the user how to view a copy of this License. (Exception: if the Program itself is interactive but does not normally print such an announcement, your work based on the Program is not required to print an announcement.)

These requirements apply to the modified work as a whole. If

identifiable sections of that work are not derived from the Program, and can be reasonably considered independent and separate works in themselves, then this License, and its terms, do not apply to those sections when you distribute them as separate works. But when you distribute the same sections as part of a whole which is a work based on the Program, the distribution of the whole must be on the terms of this License, whose permissions for other licensees extend to the entire whole, and thus to each and every part regardless of who wrote it.

Thus, it is not the intent of this section to claim rights or contest your rights to work written entirely by you; rather, the intent is to exercise the right to control the distribution of derivative or collective works based on the Program.

In addition, mere aggregation of another work not based on the Program with the Program (or with a work based on the Program) on a volume of a storage or distribution medium does not bring the other work under the scope of this License.

3. You may copy and distribute the Program (or a work based on it, under Section 2) in object code or executable form under the terms of Sections 1 and 2 above provided that you also do one of the following:
 - a) Accompany it with the complete corresponding machine-readable source code, which must be distributed under the terms of Sections 1 and 2 above on a medium customarily used for software interchange; or,
 - b) Accompany it with a written offer, valid for at least three years, to give any third party, for a charge no more than your cost of physically performing

source distribution, a complete machine-readable copy of the corresponding source code, to be distributed under the terms of Sections 1 and 2 above on a medium customarily used for software interchange; or,

- c) Accompany it with the information you received as to the offer to distribute corresponding source code. (This alternative is allowed only for noncommercial distribution and only if you received the program in object code or executable form with such an offer, in accord with Subsection b above.)

The source code for a work means the preferred form of the work for making modifications to it. For an executable work, complete source code means all the source code for all modules it contains, plus any associated interface definition files, plus the scripts used to control compilation and installation of the executable. However, as a special exception, the source code distributed need not include anything that is normally distributed (in either source or binary form) with the major components (compiler, kernel, and so on) of the operating system on which the executable runs, unless that component itself accompanies the executable.

If distribution of executable or object code is made by offering access to copy from a designated place, then offering equivalent access to copy the source code from the same place counts as distribution of the source code, even though third parties are not compelled to copy the source along with the object code.

4. You may not copy, modify, sublicense, or distribute the Program except as expressly provided under this License. Any attempt otherwise to copy,

modify, sublicense or distribute the Program is void, and will automatically terminate your rights under this License. However, parties who have received copies, or rights, from you under this License will not have their licenses terminated so long as such parties remain in full compliance.

5. You are not required to accept this License, since you have not signed it. However, nothing else grants you permission to modify or distribute the Program or its derivative works. These actions are prohibited by law if you do not accept this License. Therefore, by modifying or distributing the Program (or any work based on the Program), you indicate your acceptance of this License to do so, and all its terms and conditions for copying, distributing or modifying the Program or works based on it.

6. Each time you redistribute the Program (or any work based on the Program), the recipient automatically receives a license from the original licensor to copy, distribute or modify the Program subject to these terms and conditions. You may not impose any further restrictions on the recipients' exercise of the rights granted herein. You are not responsible for enforcing compliance by third parties to this License.

7. If, as a consequence of a court judgment or allegation of patent infringement or for any other reason (not limited to patent issues), conditions are imposed on you (whether by court order, agreement or otherwise) that contradict the conditions of this License, they do not excuse you from the conditions of this License. If you cannot distribute so as to satisfy simultaneously your obligations under this

License and any other pertinent obligations, then as a consequence you may not distribute the Program at all. For example, if a patent license would not permit royalty-free redistribution of the Program by all those who receive copies directly or indirectly through you, then the only way you could satisfy both it and this License would be to refrain entirely from distribution of the Program.

If any portion of this section is held invalid or unenforceable under any particular circumstance, the balance of the section is intended to apply and the section as a whole is intended to apply in other circumstances.

It is not the purpose of this section to induce you to infringe any patents or other property right claims or to contest validity of any such claims; this section has the sole purpose of protecting the integrity of the free software distribution system, which is implemented by public license practices. Many people have made generous contributions to the wide range of software distributed through that system in reliance on consistent application of that system; it is up to the author/donor to decide if he or she is willing to distribute software through any other system and a licensee cannot impose that choice.

This section is intended to make thoroughly clear what is believed to be a consequence of the rest of this License.

8. If the distribution and/or use of the Program is restricted in certain countries either by patents or by copyrighted interfaces, the original copyright holder who places the Program under this License may add an explicit geographical distribution limitation excluding those countries, so that distribution is permitted

only in or among countries not thus excluded. In such case, this License incorporates the limitation as if written in the body of this License.

9. The Free Software Foundation may publish revised and/or new versions of the General Public License from time to time. Such new versions will be similar in spirit to the present version, but may differ in detail to address new problems or concerns.

Each version is given a distinguishing version number. If the Program specifies a version number of this License which applies to it and "any later version", you have the option of following the terms and conditions either of that version or of any later version published by the Free Software Foundation. If the Program does not specify a version number of this License, you may choose any version ever published by the Free Software Foundation.

10. If you wish to incorporate parts of the Program into other free programs whose distribution conditions are different, write to the author to ask for permission. For software which is copyrighted by the Free Software Foundation, write to the Free Software Foundation; we sometimes make exceptions for this. Our decision will be guided by the two goals of preserving the free status of all derivatives of our free software and of promoting the sharing and reuse of software generally.

NO WARRANTY

11. BECAUSE THE PROGRAM IS LICENSED FREE OF CHARGE, THERE IS NO WARRANTY FOR THE PROGRAM, TO THE EXTENT PERMITTED BY APPLICABLE LAW. EXCEPT WHEN OTHERWISE STATED IN WRITING THE

COPYRIGHT HOLDERS AND/OR OTHER PARTIES PROVIDE THE PROGRAM "AS IS" WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE PROGRAM IS WITH YOU. SHOULD THE PROGRAM PROVE DEFECTIVE, YOU ASSUME THE COST OF ALL NECESSARY SERVICING, REPAIR OR CORRECTION.

12. IN NO EVENT UNLESS REQUIRED BY APPLICABLE LAW OR AGREED TO IN WRITING WILL ANY COPYRIGHT HOLDER, OR ANY OTHER PARTY WHO MAY MODIFY AND/OR REDISTRIBUTE THE PROGRAM AS PERMITTED ABOVE, BE LIABLE TO YOU FOR DAMAGES, INCLUDING ANY GENERAL, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OR INABILITY TO USE THE PROGRAM (INCLUDING BUT NOT LIMITED TO LOSS OF DATA OR DATA BEING RENDERED INACCURATE OR LOSSES SUSTAINED BY YOU OR THIRD PARTIES OR A FAILURE OF THE PROGRAM TO OPERATE WITH ANY OTHER PROGRAMS), EVEN IF SUCH HOLDER OR OTHER PARTY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.

END OF TERMS AND CONDITIONS

How to Apply These Terms to Your New Programs

If you develop a new program, and you want it to be of the greatest possible use to the public, the best way to achieve this is to make it free software which everyone can redistribute and change under these terms.

To do so, attach the following notices to the program. It is safest to attach them to the start of each source file to most effectively convey the exclusion of warranty; and each file should have at least the "copyright" line and a pointer to where the full notice is found.

```
<one line to give the program's
name and a brief idea of what it
does.>
Copyright (C) <year> <name of
author>
```

This program is free software; you can redistribute it and/or modify it under the terms of the GNU General Public License as published by the Free Software Foundation; either version 2 of the License, or (at your option) any later version.

This program is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the GNU General Public License for more details.

You should have received a copy of the GNU General Public License along with this program; if not, write to the Free Software Foundation, Inc., 51 Franklin Street, Fifth Floor, Boston, MA 02110-1301 USA.

Also add information on how to contact you by electronic and paper mail.

If the program is interactive, make it output a short notice like this when it starts in an interactive mode:

```
Gnomovision version 69,
Copyright (C) year name of
author Gnomovision comes with
ABSOLUTELY NO
WARRANTY; for details type
`show w'.
```

This is free software, and you are welcome to redistribute it under certain conditions; type `show c' for details.

The hypothetical commands `show w' and `show c' should show the appropriate parts of the General Public License. Of course, the commands you use may be called something other than `show w' and `show c'; they could even be mouse-clicks or menu items--whatever suits your program.

You should also get your employer (if you work as a programmer) or your school, if any, to sign a "copyright disclaimer" for the program, if necessary. Here is a sample; alter the names:

```
Yoyodyne, Inc., hereby
disclaims all copyright interest
in the program `Gnomovision'
(which makes passes at
compilers) written by James
Hacker.
```

```
<signature of Ty Coon>, 1 April
1989
Ty Coon, President of Vice
```

This General Public License does not permit incorporating your program

into proprietary programs. If your program is a subroutine library, you may consider it more useful to permit linking proprietary applications with the library. If this is what you want to do, use the GNU Lesser General Public License instead of this License.

GNU LESSER GENERAL PUBLIC LICENSE

Version 2.1, February 1999

Copyright (C) 1991, 1999 Free Software Foundation, Inc. 51 Franklin Street, Fifth Floor, Boston, MA 02110-1301 USA
Everyone is permitted to copy and distribute verbatim copies of this license document, but changing it is not allowed.

[This is the first released version of the Lesser GPL. It also counts as the successor of the GNU Library Public License, version 2, hence the version number 2.1.]

Preamble

The licenses for most software are designed to take away your freedom to share and change it. By contrast, the GNU General Public Licenses are intended to guarantee your freedom to share and change free software--to make sure the software is free for all its users.

This license, the Lesser General Public License, applies to some specially designated software packages--typically libraries--of the Free Software Foundation and other authors who decide to use it. You can use it too, but we suggest you first think carefully about whether this license or the ordinary General Public License is the better strategy to use in any particular case, based on the explanations below.

When we speak of free software, we are referring to freedom of use, not price. Our General Public Licenses are designed to make sure that you have the freedom to distribute copies

of free software (and charge for this service if you wish); that you receive source code or can get it if you want it; that you can change the software and use pieces of it in new free programs; and that you are informed that you can do these things.

To protect your rights, we need to make restrictions that forbid distributors to deny you these rights or to ask you to surrender these rights. These restrictions translate to certain responsibilities for you if you distribute copies of the library or if you modify it.

For example, if you distribute copies of the library, whether gratis or for a fee, you must give the recipients all the rights that we gave you. You must make sure that they, too, receive or can get the source code. If you link other code with the library, you must provide complete object files to the recipients, so that they can relink them with the library after making changes to the library and recompiling it. And you must show them these terms so they know their rights.

We protect your rights with a two-step method: (1) we copyright the library, and (2) we offer you this license, which gives you legal permission to copy, distribute and/or modify the library.

To protect each distributor, we want to make it very clear that there is no warranty for the free library. Also, if the library is modified by someone else and passed on, the recipients should know that what they have is not the original version, so that the original author's reputation will not be affected by problems that might be introduced by others.

Finally, software patents pose a constant threat to the existence of any free program. We wish to make sure that a company cannot effectively restrict the users of a free program by obtaining a restrictive license from a patent holder.

Therefore, we insist that any patent license obtained for a version of the library must be consistent with the full freedom of use specified in this license.

Most GNU software, including some libraries, is covered by the ordinary GNU General Public License. This license, the GNU Lesser General Public License, applies to certain designated libraries, and is quite different from the ordinary General Public License. We use this license for certain libraries in order to permit linking those libraries into non-free programs.

When a program is linked with a library, whether statically or using a shared library, the combination of the two is legally speaking a combined work, a derivative of the original library. The ordinary General Public License therefore permits such linking only if the entire combination fits its criteria of freedom. The Lesser General Public License permits more lax criteria for linking other code with the library.

We call this license the "Lesser" General Public License because it does Less to protect the user's freedom than the ordinary General Public License. It also provides other free software developers Less of an advantage over competing non-free programs. These disadvantages are the reason we use the ordinary General Public License for many libraries. However, the Lesser license provides advantages in certain special circumstances.

For example, on rare occasions, there may be a special need to encourage the widest possible use of a certain library, so that it becomes a de-facto standard. To achieve this, non-free programs must be allowed to use the library. A more frequent case is that a free library does the same job as widely used non-free libraries. In this case, there is little to gain by limiting the free library to free

software only, so we use the Lesser General Public License.

In other cases, permission to use a particular library in non-free programs enables a greater number of people to use a large body of free software. For example, permission to use the GNU C Library in non-free programs enables many more people to use the whole GNU operating system, as well as its variant, the GNU/Linux operating system.

Although the Lesser General Public License is Less protective of the users' freedom, it does ensure that the user of a program that is linked with the Library has the freedom and the wherewithal to run that program using a modified version of the Library.

The precise terms and conditions for copying, distribution and modification follow. Pay close attention to the difference between a "work based on the library" and a "work that uses the library". The former contains code derived from the library, whereas the latter must be combined with the library in order to run.

TERMS AND CONDITIONS FOR COPYING, DISTRIBUTION AND MODIFICATION

0. This License Agreement applies to any software library or other program which contains a notice placed by the copyright holder or other authorized party saying it may be distributed under the terms of this Lesser General Public License (also called "this License"). Each licensee is addressed as "you".

A "library" means a collection of software functions and/or data prepared so as to be conveniently linked with application programs

(which use some of those functions and data) to form executables.

The "Library", below, refers to any such software library or work which has been distributed under these terms. A "work based on the Library" means either the Library or any derivative work under copyright law: that is to say, a work containing the Library or a portion of it, either verbatim or with modifications and/or translated straightforwardly into another language. (Hereinafter, translation is included without limitation in the term "modification".)

"Source code" for a work means the preferred form of the work for making modifications to it. For a library, complete source code means all the source code for all modules it contains, plus any associated interface definition files, plus the scripts used to control compilation and installation of the library.

Activities other than copying, distribution and modification are not covered by this License; they are outside its scope. The act of running a program using the Library is not restricted, and output from such a program is covered only if its contents constitute a work based on the Library (independent of the use of the Library in a tool for writing it). Whether that is true depends on what the Library does and what the program that uses the Library does.

1. You may copy and distribute verbatim copies of the Library's complete source code as you receive it, in any medium, provided that you conspicuously and appropriately publish on each copy an appropriate copyright notice and disclaimer of warranty; keep intact all the notices that refer to this License and to the absence of any warranty; and distribute a copy of this License along with the Library.

You may charge a fee for the physical act of transferring a copy, and you may at your option offer warranty protection in exchange for a fee.

2. You may modify your copy or copies of the Library or any portion of it, thus forming a work based on the Library, and copy and distribute such modifications or work under the terms of Section 1 above, provided that you also meet all of these conditions:
 - a) The modified work must itself be a software library.
 - b) You must cause the files modified to carry prominent notices stating that you changed the files and the date of any change.
 - c) You must cause the whole of the work to be licensed at no charge to all third parties under the terms of this License.
 - d) If a facility in the modified Library refers to a function or a table of data to be supplied by an application program that uses the facility, other than as an argument passed when the facility is invoked, then you must make a good faith effort to ensure that, in the event an application does not supply such function or table, the facility still operates, and performs whatever part of its purpose remains meaningful.

(For example, a function in a library to compute square roots has a purpose that is entirely well-defined independent of the application. Therefore, Subsection 2d requires that any application-supplied function or table used by this function must be optional: if

the application does not supply it, the square root function must still compute square roots.)

These requirements apply to the modified work as a whole. If identifiable sections of that work are not derived from the Library, and can be reasonably considered independent and separate works in themselves, then this License, and its terms, do not apply to those sections when you distribute them as separate works. But when you distribute the same sections as part of a whole which is a work based on the Library, the distribution of the whole must be on the terms of this License, whose permissions for other licensees extend to the entire whole, and thus to each and every part regardless of who wrote it.

Thus, it is not the intent of this section to claim rights or contest your rights to work written entirely by you; rather, the intent is to exercise the right to control the distribution of derivative or collective works based on the Library.

In addition, mere aggregation of another work not based on the Library with the Library (or with a work based on the Library) on a volume of a storage or distribution medium does not bring the other work under the scope of this License.

3. You may opt to apply the terms of the ordinary GNU General Public License instead of this License to a given copy of the Library. To do this, you must alter all the notices that refer to this License, so that they refer to the ordinary GNU General Public License, version 2, instead of to this License. (If a newer version than version 2 of the ordinary GNU General

Public License has appeared, then you can specify that version instead if you wish.) Do not make any other change in these notices.

Once this change is made in a given copy, it is irreversible for that copy, so the ordinary GNU General Public License applies to all subsequent copies and derivative works made from that copy.

This option is useful when you wish to copy part of the code of the Library into a program that is not a library.

4. You may copy and distribute the Library (or a portion or derivative of it, under Section 2) in object code or executable form under the terms of Sections 1 and 2 above provided that you accompany it with the complete corresponding machine-readable source code, which must be distributed under the terms of Sections 1 and 2 above on a medium customarily used for software interchange.

If distribution of object code is made by offering access to copy from a designated place, then offering equivalent access to copy the source code from the same place satisfies the requirement to distribute the source code, even though third parties are not compelled to copy the source along with the object code.

5. A program that contains no derivative of any portion of the Library, but is designed to work with the Library by being compiled or linked with it, is called a "work that uses the Library". Such a work, in isolation, is not a derivative work of the Library, and therefore falls outside the scope of this License.

However, linking a "work that uses the Library" with the Library creates an executable that is a derivative of

the Library (because it contains portions of the Library), rather than a "work that uses the library". The executable is therefore covered by this License.

Section 6 states terms for distribution of such executables.

When a "work that uses the Library" uses material from a header file that is part of the Library, the object code for the work may be a derivative work of the Library even though the source code is not.

Whether this is true is especially significant if the work can be linked without the Library, or if the work is itself a library. The threshold for this to be true is not precisely defined by law.

If such an object file uses only numerical parameters, data structure layouts and accessors, and small macros and small inline functions (ten lines or less in length), then the use of the object file is unrestricted, regardless of whether it is legally a derivative work. (Executables containing this object code plus portions of the Library will still fall under Section 6.)

Otherwise, if the work is a derivative of the Library, you may distribute the object code for the work under the terms of Section 6.

Any executables containing that work also fall under Section 6, whether or not they are linked directly with the Library itself.

6. As an exception to the Sections above, you may also combine or link a "work that uses the Library" with the Library to produce a work containing portions of the Library, and distribute that work under terms of your choice, provided that the terms permit modification of the work for the customer's own use and reverse engineering for debugging such modifications.

You must give prominent notice with each copy of the work that the

Library is used in it and that the Library and its use are covered by this License. You must supply a copy of this License. If the work during execution displays copyright notices, you must include the copyright notice for the Library among them, as well as a reference directing the user to the copy of this License. Also, you must do one of these things:

- a) Accompany the work with the complete corresponding machine-readable source code for the Library including whatever changes were used in the work (which must be distributed under Sections 1 and 2 above); and, if the work is an executable linked with the Library, with the complete machine-readable "work that uses the Library", as object code and/or source code, so that the user can modify the Library and then relink to produce a modified executable containing the modified Library. (It is understood that the user who changes the contents of definitions files in the Library will not necessarily be able to recompile the application to use the modified definitions.)
- b) Use a suitable shared library mechanism for linking with the Library. A suitable mechanism is one that (1) uses at run time a copy of the library already present on the user's computer system, rather than copying library functions into the executable, and (2) will operate properly with a modified version of the library, if the user installs one, as long as the modified version is interface-compatible with

the version that the work was made with.

- c) Accompany the work with a written offer, valid for at least three years, to give the same user the materials specified in Subsection 6a, above, for a charge no more than the cost of performing this distribution.
- d) If distribution of the work is made by offering access to copy from a designated place, offer equivalent access to copy the above specified materials from the same place.
- e) Verify that the user has already received a copy of these materials or that you have already sent this user a copy.

For an executable, the required form of the "work that uses the Library" must include any data and utility programs needed for reproducing the executable from it. However, as a special exception, the materials to be distributed need not include anything that is normally distributed (in either source or binary form) with the major components (compiler, kernel, and so on) of the operating system on which the executable runs, unless that component itself accompanies the executable.

It may happen that this requirement contradicts the license restrictions of other proprietary libraries that do not normally accompany the operating system. Such a contradiction means you cannot use both them and the Library together in an executable that you distribute.

7. You may place library facilities that are a work based on the Library side-by-side in a single library together with other library facilities not covered by this License, and distribute such a combined library, provided that the separate distribution of

the work based on the Library and of the other library facilities is otherwise permitted, and provided that you do these two things:

- a) Accompany the combined library with a copy of the same work based on the Library, uncombined with any other library facilities. This must be distributed under the terms of the Sections above.
 - b) Give prominent notice with the combined library of the fact that part of it is a work based on the Library, and explaining where to find the accompanying uncombined form of the same work.
8. You may not copy, modify, sublicense, link with, or distribute the Library except as expressly provided under this License. Any attempt otherwise to copy, modify, sublicense, link with, or distribute the Library is void, and will automatically terminate your rights under this License. However, parties who have received copies, or rights, from you under this License will not have their licenses terminated so long as such parties remain in full compliance.
9. You are not required to accept this License, since you have not signed it. However, nothing else grants you permission to modify or distribute the Library or its derivative works. These actions are prohibited by law if you do not accept this License. Therefore, by modifying or distributing the Library (or any work based on the Library), you indicate your acceptance of this License to do so, and all its terms and conditions for copying, distributing or modifying the Library or works based on it.

10. Each time you redistribute the Library (or any work based on the Library), the recipient automatically receives a license from the original licensor to copy, distribute, link with or modify the Library subject to these terms and conditions. You may not impose any further restrictions on the recipients' exercise of the rights granted herein. You are not responsible for enforcing compliance by third parties with this License.

11. If, as a consequence of a court judgment or allegation of patent infringement or for any other reason (not limited to patent issues), conditions are imposed on you (whether by court order, agreement or otherwise) that contradict the conditions of this License, they do not excuse you from the conditions of this License. If you cannot distribute so as to satisfy simultaneously your obligations under this License and any other pertinent obligations, then as a consequence you may not distribute the Library at all. For example, if a patent license would not permit royalty-free redistribution of the Library by all those who receive copies directly or indirectly through you, then the only way you could satisfy both it and this License would be to refrain entirely from distribution of the Library.

If any portion of this section is held invalid or unenforceable under any particular circumstance, the balance of the section is intended to apply, and the section as a whole is intended to apply in other circumstances.

It is not the purpose of this section to induce you to infringe any patents or other property right claims or to contest validity of any such claims; this section has the sole purpose of protecting the integrity of the free

software distribution system which is implemented by public license practices. Many people have made generous contributions to the wide range of software distributed through that system in reliance on consistent application of that system; it is up to the author/donor to decide if he or she is willing to distribute software through any other system and a licensee cannot impose that choice.

This section is intended to make thoroughly clear what is believed to be a consequence of the rest of this License.

12. If the distribution and/or use of the Library is restricted in certain countries either by patents or by copyrighted interfaces, the original copyright holder who places the Library under this License may add an explicit geographical distribution limitation excluding those countries, so that distribution is permitted only in or among countries not thus excluded. In such case, this License incorporates the limitation as if written in the body of this License.

13. The Free Software Foundation may publish revised and/or new versions of the Lesser General Public License from time to time. Such new versions will be similar in spirit to the present version, but may differ in detail to address new problems or concerns.

Each version is given a distinguishing version number. If the Library specifies a version number of this License which applies to it and "any later version", you have the option of following the terms and conditions either of that version or of any later version published by the Free Software Foundation. If the Library does not specify a license version number, you may choose any version ever published by the Free Software Foundation.

14. If you wish to incorporate parts of the Library into other free programs whose distribution conditions are incompatible with these, write to the author to ask for permission. For software which is copyrighted by the Free Software Foundation, write to the Free Software Foundation; we sometimes make exceptions for this. Our decision will be guided by the two goals of preserving the free status of all derivatives of our free software and of promoting the sharing and reuse of software generally.

NO WARRANTY

15. BECAUSE THE LIBRARY IS LICENSED FREE OF CHARGE, THERE IS NO WARRANTY FOR THE LIBRARY, TO THE EXTENT PERMITTED BY APPLICABLE LAW. EXCEPT WHEN OTHERWISE STATED IN WRITING THE COPYRIGHT HOLDERS AND/OR OTHER PARTIES PROVIDE THE LIBRARY "AS IS" WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE LIBRARY IS WITH YOU. SHOULD THE LIBRARY PROVE DEFECTIVE, YOU ASSUME THE COST OF ALL NECESSARY SERVICING, REPAIR OR CORRECTION.
16. IN NO EVENT UNLESS REQUIRED BY APPLICABLE LAW OR AGREED TO IN WRITING WILL ANY COPYRIGHT HOLDER, OR ANY OTHER

PARTY WHO MAY MODIFY AND/OR REDISTRIBUTE THE LIBRARY AS PERMITTED ABOVE, BE LIABLE TO YOU FOR DAMAGES, INCLUDING ANY GENERAL, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OR INABILITY TO USE THE LIBRARY (INCLUDING BUT NOT LIMITED TO LOSS OF DATA OR DATA BEING RENDERED INACCURATE OR LOSSES SUSTAINED BY YOU OR THIRD PARTIES OR A FAILURE OF THE LIBRARY TO OPERATE WITH ANY OTHER SOFTWARE), EVEN IF SUCH HOLDER OR OTHER PARTY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.

END OF TERMS AND CONDITIONS

How to Apply These Terms to Your New Libraries

If you develop a new library, and you want it to be of the greatest possible use to the public, we recommend making it free software that everyone can redistribute and change. You can do so by permitting redistribution under these terms (or, alternatively, under the terms of the ordinary General Public License).

To apply these terms, attach the following notices to the library. It is safest to attach them to the start of each source file to most effectively convey the exclusion of warranty; and each file should have at least the "copyright" line and a pointer to where the full notice is found.

<one line to give the library's name and a brief idea of what it does.>

Copyright (C) <year> <name of author>

This library is free software; you can redistribute it and/or modify it under the terms of the GNU Lesser General Public License as published by the Free Software Foundation; either version 2.1 of the License, or (at your option) any later version.

This library is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the GNU Lesser General

Public License for more details.

You should have received a copy of the GNU Lesser General Public License along with this library; if not, write to the Free Software Foundation, Inc., 51 Franklin Street, Fifth Floor, Boston, MA 02110-1301 USA

Also add information on how to contact you by electronic and paper mail.

You should also get your employer (if you work as a programmer) or your school, if any, to sign a "copyright disclaimer" for the library, if necessary. Here is a sample; alter the names:

Yoyodyne, Inc., hereby disclaims all copyright interest in the library `Frob' (a library for tweaking knobs) written by James Random Hacker.

<signature of Ty Coon>, 1
April 1990
Ty Coon, President of Vice

That's all there is to it!

Specifications

Recording Features

Recording Video Format

MPEG-4 AVC/H.264

Recording Audio Format

LPCM, AAC LC

Recording File Format

XAVC S, MP4

Recording Media

Internal HDD (4TB)

External USB Storage

Network (CIFS)

DVD-R

BD-R/BD-R DL

BD-RE/BD-RE DL

Input Resolution

4096 × 2160

3840 × 2160

1920 × 1080

Recording Resolution

3840 × 2160

1920 × 1080

Recording Bit Rate (4K)

150 Mbps (Best)

100 Mbps (High)

60 Mbps (Standard)

Recording Bit Rate (2K)

24 Mbps (Best)

18 Mbps (High)

12 Mbps (Standard)

Connectors

Input Connectors

3G-SDI (BNC type) (4)

AUDIO (Stereo mini jack) (1)

MIC (Stereo mini jack) (1)

AC Inlet (3-pin) (1)

Output Connectors

3G-SDI (BNC type) (4)

HDMI (Type A) (1)

AUDIO (Stereo mini jack) (1)

Other Interfaces

USB 3.0 (Type A) (2)

USB 2.0 (Type A) (4)

USB 2.0 (Type B) (1)

Network (RJ-45, 1000 Base-T/100 Base) (1)

REMOTE RS-232C (D-sub 9-pin) (1)

REMOTE contact switch (stereo mini jack) (4)

Equipotential

General

Power Requirements

100 V to 240 V AC, 50/60 Hz

Input Current

1.25 A to 0.52 A

Operating Temperature

5 °C to 40 °C (41 °F to 104 °F)

Operating Humidity

20% to 80% (Maximum wet-bulb temperature: 30 °C (86 °F)) (no condensation allowed)

Operating Pressure

700 hPa to 1,060 hPa

Storage and Transport Temperature

−20 °C to +60 °C (−4 °F to +140 °F)

Storage and Transport Humidity

20% to 90% (Maximum wet-bulb temperature: 30 °C (86 °F)) (no condensation allowed)

Storage and Transport Pressure

700 hPa to 1,060 hPa

Mass

Approx. 6.5 kg (Approx. 14 lb. 5.3 oz.)

Dimensions

305.0 × 329.0 × 115.5 mm (12 1/8 × 13 × 4 5/8 in.)
(W × D × H) (including longest protrusions)

Supplied Items

Before Using this Unit (1)

CD-ROM (Instructions for Use, PROTOCOL
MANUAL) (1)

Warranty booklet (1)

Service Contact List (1)

Infrared remote control unit (RM-M010) (1)

Information for Customers in Europe (1)

Separately-Sold Accessories


Foot Switch (FS-24)

Caution

The FS-24 has an Ingress Protection rating of IPX3. Therefore, do not operate it environments exposed to splashing liquids (e.g., surgical operating rooms). For safety, use a device with a rating of IPX6 or higher when operating in such environments.

REMOTE contact switch connectors

Connector specifications (stereo mini jack)

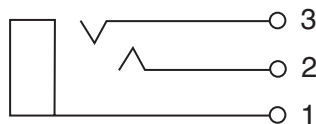

1 GND

2 TRIGGER SIGNAL (TTL)

If the status changes for at least 100 msec, the action configured in the REMOTE contact switch settings will start.

3 BUSY STATUS (TTL)

HIGH level will be initiated according to the REMOTE contact switch settings.

Medical Specifications

Protection against electric shock:

Class I

Protection against harmful ingress of water:

Ordinary

Degree of safety in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide:

Not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide

Mode of operation:

Continuous

Design and specifications are subject to change without notice.

Notes

- Always make a test recording, and verify that it was recorded successfully.
SONY WILL NOT BE LIABLE FOR DAMAGES OF ANY KIND INCLUDING, BUT NOT LIMITED TO, COMPENSATION OR REIMBURSEMENT ON ACCOUNT OF FAILURE OF THIS UNIT OR ITS RECORDING MEDIA, EXTERNAL STORAGE SYSTEMS OR ANY OTHER MEDIA OR STORAGE SYSTEMS TO RECORD CONTENT OF ANY TYPE.
- Always verify that the unit is operating properly before use. SONY WILL NOT BE LIABLE FOR DAMAGES OF ANY KIND INCLUDING, BUT NOT LIMITED TO, COMPENSATION OR REIMBURSEMENT ON ACCOUNT OF THE LOSS OF PRESENT OR PROSPECTIVE PROFITS DUE TO FAILURE OF THIS UNIT, EITHER DURING THE WARRANTY PERIOD OR AFTER EXPIRATION OF THE WARRANTY, OR FOR ANY OTHER REASON WHATSOEVER.
- SONY WILL NOT BE LIABLE FOR CLAIMS OF ANY KIND MADE BY USERS OF THIS UNIT OR MADE BY THIRD PARTIES.
- SONY WILL NOT BE LIABLE FOR THE LOSS, REPAIR, OR REPRODUCTION OF ANY DATA RECORDED ON THE INTERNAL STORAGE SYSTEM, RECORDING MEDIA, EXTERNAL STORAGE SYSTEMS OR ANY OTHER MEDIA OR STORAGE SYSTEMS.
- SONY WILL NOT BE LIABLE FOR THE TERMINATION OR DISCONTINUATION OF ANY SERVICES RELATED TO THIS UNIT THAT MAY RESULT DUE TO CIRCUMSTANCES OF ANY KIND.

Index

Numerics

4K Convert Modes 40

A

Auto delete function 50, 78

[Auto Live] tab 65

Auto live 65

B

BD/DVD indicator 21

C

Capturing still images 44, 51

Ch1/Ch2 simultaneous recording 52

Close folder 45

[CMS] tab 67

Compatible discs 30

[Contact Switch] tab 70

Contact switches 70

[Control] tab 69

Copying 57

[C-Store] tab 79

D

Deleting data 59

[Device 1] tab 67

[Device 2] tab 69

Device settings 67

DICOM settings 78

[DNS Server] tab 72

E

Editing cases 76

Editing categories 77

Editing the case list 76

Editing the doctor list 75

[End Recording] screen 45

Error messages 85

External media 47

F

[Facility] tab 66

[File Server] tab 72, 74

File servers 72, 74

Front panel 20

Function keys 38

Function settings 64

G

[General] tab 73

[General 1] tab 65

[General 2] tab 66

H

Hard keyboards 47

HDD indicator 21

I

Image list 54, 55

Image quality settings 34

Image search 54

[Input Detection] tab 64

[Input Signal] tab 64

Input signal settings 34

Input signals 34

[IP Address] tab 72

K

Keyboard mode 47

L

Language and time settings 62

[Local Station 1] tab 80

[Local Station 2] tab 80

Logos 66, 69

M

Mouse 81

[MWM Server 1] tab 78

[MWM Server 2] tab 79

N

[Network] tab 72

Network settings 72

[NTP] tab 73

O

On/standby switch 32

On-Screen Keyboard 29

P

Password settings 71

Patient ID 64

[Patient Info.] tab 64

Patient Information 58

Playback 53

[Print] tab 74

Print settings 37

Printer settings 68

Printing 56

Protecting data 59

Q

Quick playback 44, 53

R

Reader settings 69

Rear panel 23

Recording 43, 50

Recording list 54, 55

Recording preparation 46

Registering doctors 73

Registering patients 48

Remote control unit 25

Removing protection 59

S

Save settings 36

SERVER indicator 21

[Shared] tab 73, 74

Status displays 26

Still image capture 51

Still image directory 51

System administrator settings 61

System control settings 70

T

Test recording 49

Thumbnail display 54

Touch panel and monitor 81

Troubleshooting 87

U

USB indicators 21

User settings 34

V

Video data directory 50


W

Watermark files 68

Watermarks 37



EU: Sony Europe B.V.
Da Vincilaan 7-D1, 1930 Zaventem, Belgium
UK: Sony Europe B.V.
The Heights, Brooklands, Weybridge,
Surrey KT13 0XW, United Kingdom
CH: Sony Europe B.V., Hoofddorp,
Schlieren/Switzerland Branch
Wiesenstrasse 5, 8952 Schlieren, Switzerland

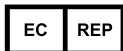

Sony Belgium, bijkantoor van
Sony Europe B.V.
Da Vincilaan 7-D1, 1930 Zaventem,
Belgium



Sony Europe B.V.
The Heights, Brooklands, Weybridge,
Surrey KT13 0XW, United Kingdom



Sony Corporation
1-7-1 Konan Minato-ku Tokyo,
108-0075 Japan