

Baaske network isolator

MED MI 1005 MB

10/100/1000 MBit/s



2 MOPP
MEANS OF
PROTECTION

5
YEARS
WARRANTY

5kV AC
10kV DC



Wall-mounted medical network isolator


The MED MI 1005 MB passive network isolator is specifically designed for wall mounting in medical environments. It provides galvanic network isolation for electronic medical devices and systems starting directly at the wall outlet. Once it is mounted to the wall, the isolator is well-protected by the housing. Accidental impacts can be absorbed, protecting the isolator and network cable connector from damage. The MED MI 1005 MB network isolator delivers dielectric strength up to 5 kV, up to 10 kV DC for 10 seconds under high-stress conditions. Coupled with Harting PushPull technology, an absolutely stable connection. No installation of drivers or additional software necessary because the MED MI 1005 MB network isolator behaves fully transparently on the data network.

Use in medical technology: In accordance with the IEC EN 60601-1 standard for medical electrical devices and systems, the MI1005 protects the network interface by two independent protective measures (2 MOPP) to the patient and the network connection is safely galvanically isolated.

Benefits

- › Saleable complete set
- › Simply plug in, no installation required
- › Protects medical devices and patients from excessive leakage currents in the data network
- › No keystone module required
- › No installation of drivers or additional software required
- › does not require its own power supply
- › Long-term availability for min. 10 years
- › Durable: 1 fault in 21 000 000 hrs.

Datas and facts.

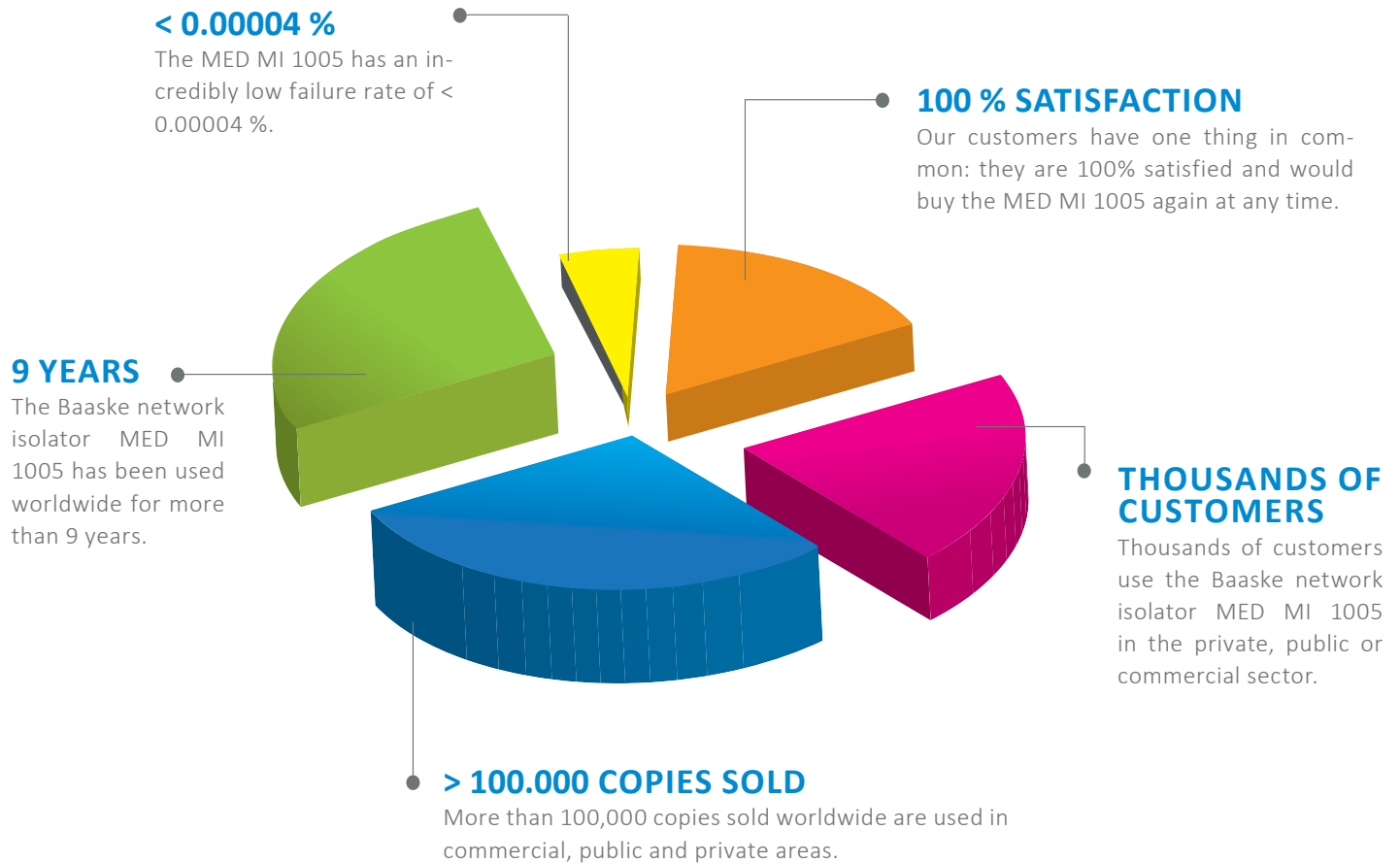


Standards and certificates

Description	
UL Recognized Component	✓
IEC 60601-1	:2005 + A1:2012 Ed. 3.1 (guttet section 17 - EMV)
IEC 60601-1-2	:2014 Ed. 4
EN 60601-1	:2006 + A1:2013 Ed. 3.1 (guttet section 17 - EMV)
EN 60601-1-2	:2015 Ed. 4
Low-voltage directive	2014/35/EU
EMC directive	14/30/EU
RoHS	✓
SVHC / REACH	✓
Lead-free	✓

We will gladly send you the corresponding EU declarations of conformity on request.

Applications

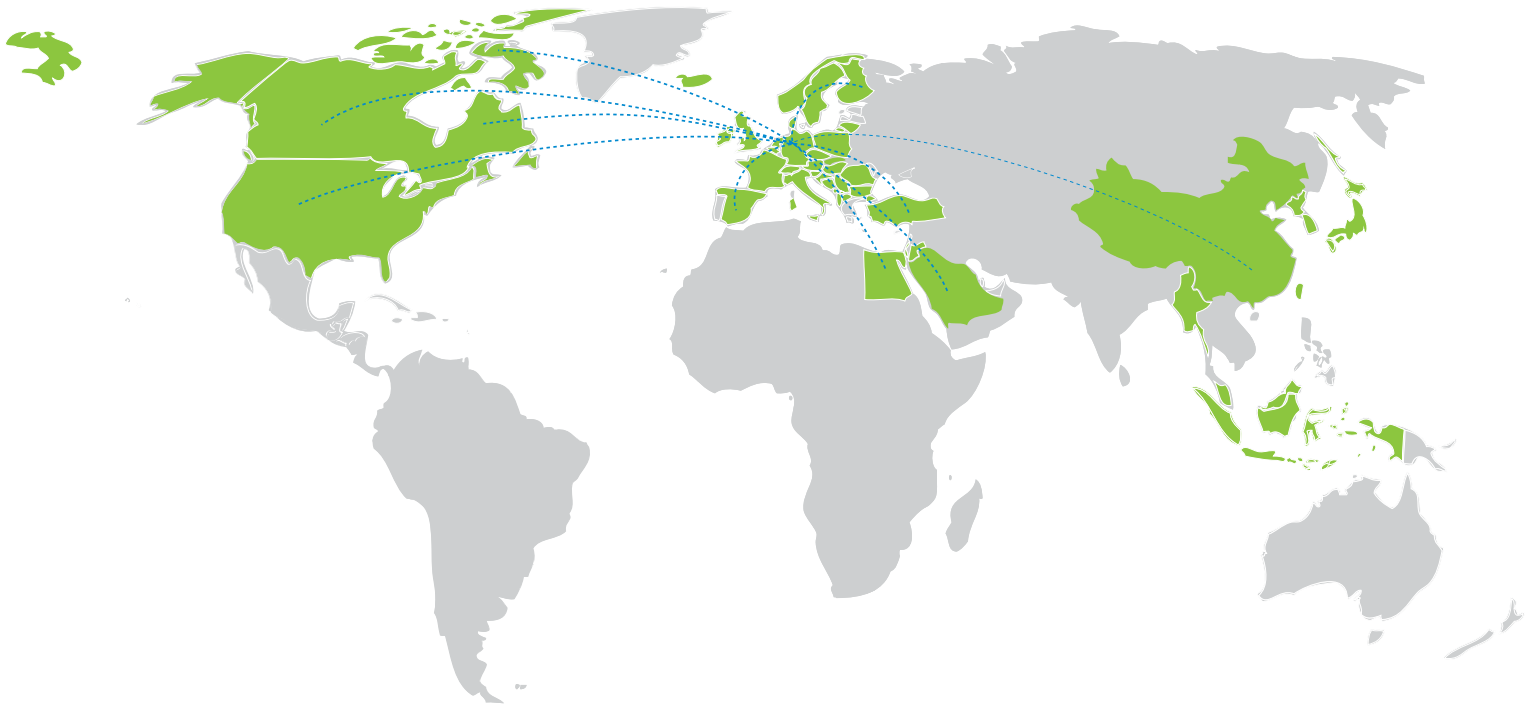


The Baaske [network isolator MED MI 1005](#) is used in more than 30 countries in the private, public or commercial sector. It reliably protects electrical devices against excessive voltages, ground loops and signal noise from the network cable. The usual areas of application are:


Application areas



In medicine: The MED MI 1005 network isolators protect patients, medical devices and users from excessive leakage currents that can flow over the data network in accordance with standards.

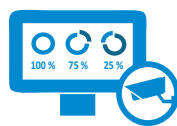

Measuring and monitoring equipment: Protects sensitive measuring and monitoring equipment that is connected to a control station via Ethernet interfaces from interference voltages and potential differences.

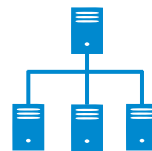


Server systems / building services engineering: To avoid failures or potential equalization currents, network isolators protect server systems that are connected over long distances using copper cabling.




Audio applications: Network isolators reliably protect against overvoltages, low-frequency alternating voltages (mains hum) and eliminate noise in the sound.

Types of construction and technical details

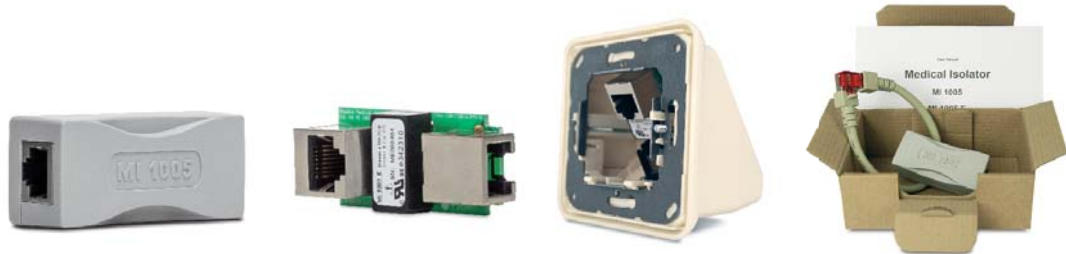


Technical Data	Flexible	Build In Isolator	Wall mounting	Complete set
Name	Netzwerk Isolator MED MI 1005	MI 1005 E Medical Build In Isolator	Netzwerk Isolator MED MI 1005 MB	Network Isolator MED MI 1005 Retail
Color case	light grey	without housing	creamy	light grey
Dimensions (DWH)	65x29x23 mm	51x23,4x18.2 mm	81x81x57 mm	65x29x23 mm
Weight approx	50 g	20 g	106 g	50 g
Operating mode	continuous operation	continuous operation	continuous operation	continuous operation
Dielectric strength	5000 V 50/60 Hz for 1 Min.			
Signal and shielding	10000 V 50/60 Hz for 10 s.			
Data throughput	10/100/1000 MBit/s			
Ports	2x RJ45 IN/OUT			
DC Stress Test (10 Sec)	10 kV			
Network Specifications	IEEE 803.2 ab 10/100/1000-BaseT; Twisted-Pair; auto-conf (behaves completely transparent in the network)			
Material	UL94V-0 / RoHS conform / REACH			
Produkt Rating (UL)	passive, electronic components, insulation DI (250 V AC / 300 V DC) tested at 5 kV AC			
Operating time between failures	0,21 x 10 ^ 8 (1 failures in 21 000 000 h)			
Return loss	min. 10 dB			
Insertion loss	max. 1,1 dB			
Standards	IEC 60601-1 (3rd Edition); IEC 60601-1-2; UL 60601 (UL Listed E342310); IEEE 802.3			
Temperature	-10 °C up to +70 °C (operation/ storage/ transport)			
Relative air humidity	10 % up to 90 % (non-condensing (operation/ storage/ transport)			
Air pressure	700 hPa up to 1060 hPa (operation/ storage/ transport)			
Warranty	5 years			
Scope of delivery	Netzwerk Isolator MED MI 1005	MI 1005 E Medical Build In Isolator	1x mounting plate; 1x cover; 1x frame; 1x MI 1005 E Medical Built In Isolator; 1x retaining plate; 1x Harting PushPull; 2x screws; 2x nut; 1x assembly instructions; 1x manual	Netzwerk Isolator MED MI 1005 , Patch Cable 25 cm, Manual, Package and serial number
Item No.	2005674	2006633	2007704	2006484

Circuit

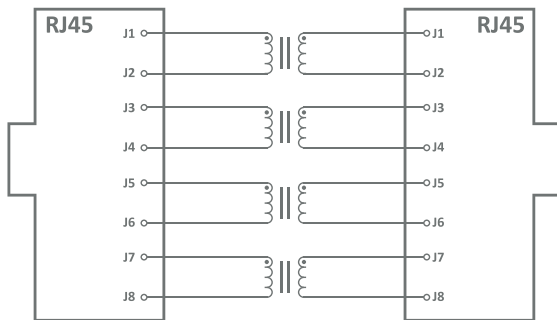

Frequency range MED MI 1005

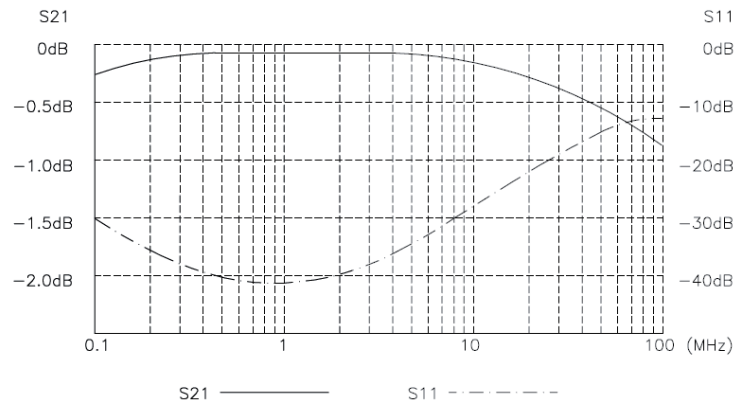

Harting PushPull Accessories



Mounting housing Flat

Black plastic angled, with integrated seal + retaining clip; mounting holes for M2, 5, IP 65/67

Item No. 2007020



Mounting housing Compact

Black plastic angled, with integrated seal + retaining clip; mounting holes for M2, 5, IP 65/67

Item No. 2006634



Mounting housing EasyInstall

Black plastic angled, with retaining clip; mounting holes for M3, IP 65/67

Item No. 2006635



CAT6 patch cable 0.25 m grey PIMF, RoHs

Category 6 S / FTP patch cable > twisted foil-shielded with braided shield (UTP) > halogen-free jacket > shielded, molded RJ45 connector > colored, molded kink protection grommet > length on sleeve > protection of the snap-in nose sleeve > EIA / TIA B.2-1, EN > 50173 ISO / IEC 11801 Class E

0.25 m	Item No. 2005670
0.50 m	Item No. 2006294
2.00 m	Item No. 2006296
5.00 m	Item No. 2006298